

Contiplex® C Set

NRFit®

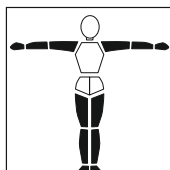Peripheral
Spinal
Epidural


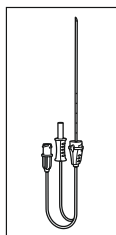
en	Continuous peripheral nerve block catheter set	4
de	Katheterset für kontinuierliche periphere Nervenblockaden	6
bg	Комплект катетри за непрекъснати периферни нервни блокове	8
cs	Set katétru na kontinuální blokádu periferních nervů	10
da	Katetersæt til kontinuerlig perifer nerveblokade	12
el	Σετ καθετήρα συνεχή περιφερικού αποκλεισμού νεύρων	14
es	Set de catéter para bloqueo nervioso periférico continuo	16
et	Pideva perifeerse närviblokaadi kateetrikomplekt	18
fi	Jatkuvatoiminen ääreishermostojen puudutuskatetrisarja	20
fr	Ensemble de cathéter pour bloc nerveux périphérique continu	22
hu	Katéterkészlet folyamatos perifériás idegblokádhoz	24
it	Set catetere per blocco continuo dei nervi periferici	26
lt	Nuolatinės periferinio nervo blokados kateterio rinkinys	28
lv	Ilgstošas perifēro nervu bloķēšanas katetru komplekts	30
nl	Katheterset voor ononderbroken perifere zenuwblockade	32
no	Katetersett til kontinuerlige, perifere nerveblokader	34
pl	Zestaw cewnika do ciągłej blokady nerwów obwodowych	36
pt	Conjunto de cateter de bloqueio de nervo periférico contínuo	38
ro	Set cateter pentru bloc nervos periferic continuu	40
ru	Набор постоянных катетеров для блокады периферийных нервов	42
sk	Súprava katétra na kontinuálnu blokádu periférneho nervu	44
sl	Komplet s katetrom za kontinuirano blokado perifernega živčevja	46
sv	Kateterset för kontinuerlig perifer nervblockad	48
th	ชุดสายสวนสำหรับการระงับความรู้สึกประสาทส่วนปลายแบบต่อเนื่อง	50
tr	Sürekli periferik sinir bloğu kateter seti	52
zh	连续外周神经阻滞导管套件	54

B | BRAUN

en Line drawings of the individual components
de Strichzeichnungen der Komponenten
bg Скици на отделните съставни части
cs Perokresby jednotlivých částí
da Stregtegninger af de enkelte dele
el Γραμμικά σχέδια των ξεχωριστών στοιχείων
es Esquemas de los componentes individuales
et Üksikute komponentide skeemid
fi Osien merkit

fr Croquis des différents constituants
hu Az egyes alkotórészek vonalas rajzai
it Schemi dei singoli componenti
lt Kiekviena detalė pažymėta brūkšniu kodu
lv Atsevišķo materiāla sastāvdaļu svītrkodsi
nl Tekeningen van de verschillende onderdelen
no Strektegninger av de enkelte bestanddeler
pl Rysunki pojedynczych części składowych
pt Traços da marca de cada componente






ro Desene liniare ale componentelor
ru Рисунки отдельных компонентов
sk Perokresbay jednotlivých častí
sl Shematski prikazi posameznih sestavnih delov
sv Bilder på de enskilda komponenterna
th ภาพส่วนประกอบต่าง ๆ
tr Komponentlerin ayrı çizimleri
zh 单独组件的图示

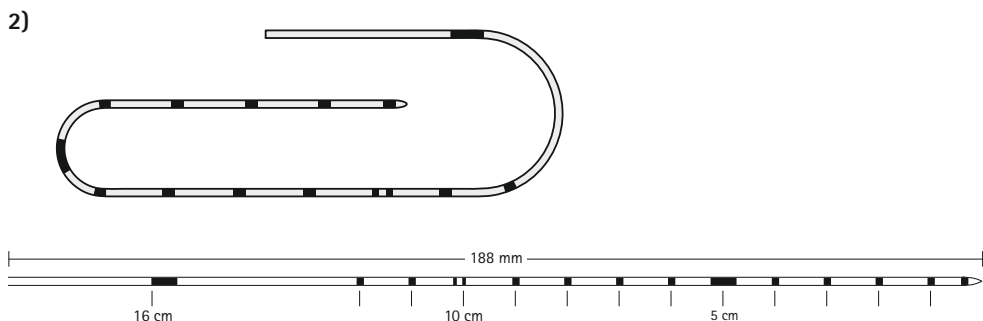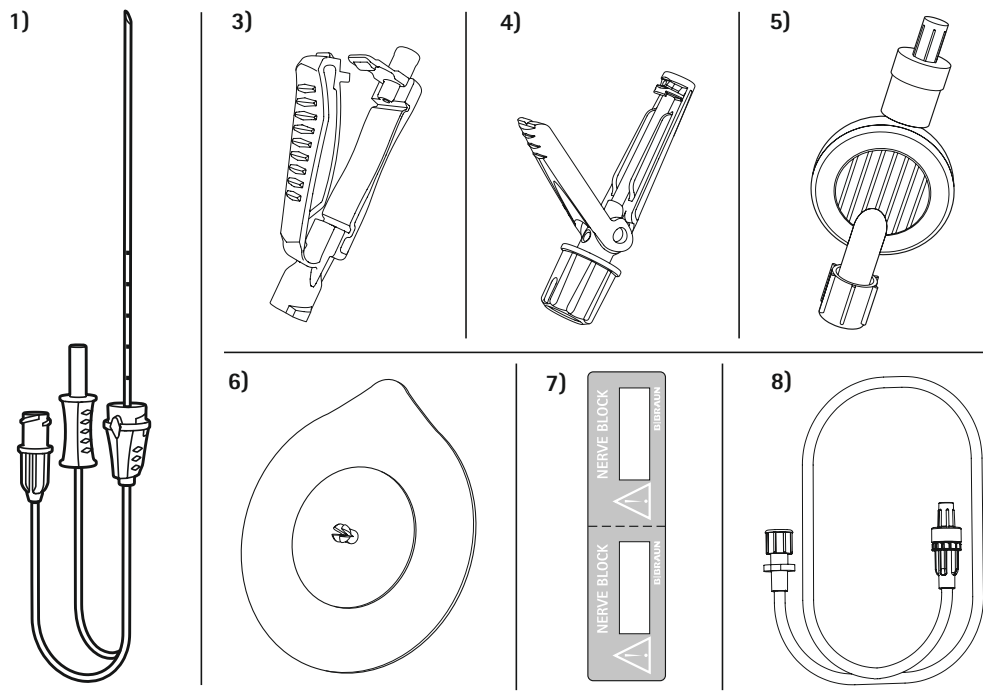en Directions
de Anwendungshinweise
bg Указания за употреба
cs Pokyny pro použití
da Betjeningsvejledning
el Υποδείξεις εφαρμογής
es Aplicaciones
et Kasutusjuhendid
fi Käyttöohje

fr Conseils d'utilisation
hu Használati útmutatás
it Modalità d'uso
lt Įspėjimai
lv Norādījumi par lietošanu
nl Gebruiksaanwijzingen
no Bruk
pl Przepisy stosowania
pt Instruções de uso

ro Instrucțiuni de folosire
ru Указания по применению
sk Aplikačné pokyny
sl Napotki za uporabo
sv Användningsanvisningar
th วิธีการใช้
tr Kullanım talimatları
zh 说明

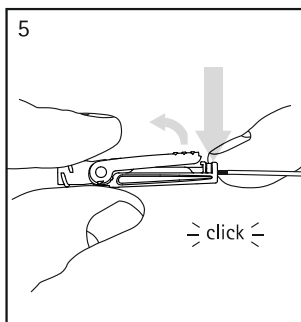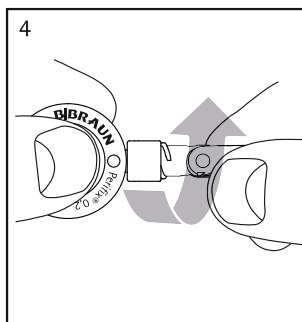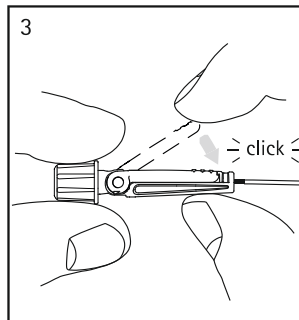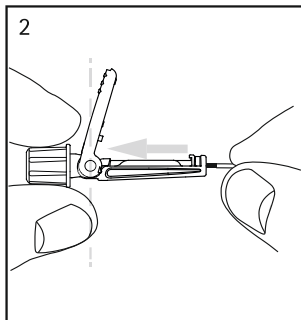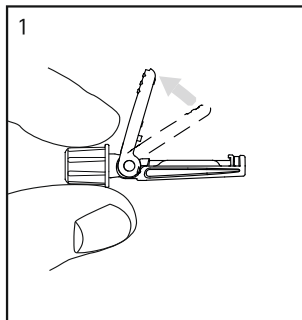en Instructions for Use

Contiplex® C Set NRFit®

- Contiplex® C needle NRFit® with integrated wire (Fig. 1)
- Contiplex® ONE Catheter (open tip) (Fig. 2)
- C-Grip (positioning device) (Fig. 3)
- Perifix® Catheter Connector NRFit® (Fig. 4)
- Perifix® Filter (Bacteria filter) NRFit® (Fig. 5)
- Perifix® PinPad (Fig. 6)
- Catheter label (Fig. 7)
- Extension tube NRFit® (Fig. 8)

For further technical characteristics such as needle and catheter size, see product packaging.

NRFit®-connector

The female NRFit® connector of the Contiplex® C devices is a standardized connector. It can be connected to syringes with male NRFit® connector. The syringes can be used for flushing the needle, aspirating and infusing anesthetics.

As risk mitigation measure for continuous plexus and peripheral nerve blocks, the catheter label included in the set mentions the respective Regional Anesthesia application.

Electrode connector

The Contiplex® product range contains a standardized electrode connector. The electrode connector can be connected to a nerve stimulator which is intended to be used for peripheral nerve block procedures. In order to achieve the optimum connection between the stimulation needle and the nerve stimulator, thereby ensuring patient safety, only use B. Braun products with other B. Braun products. B. Braun Melsungen AG does not assume any liability where its products are combined with equipment from other manufacturers.

Assessment of Magnetic Resonance (MR) compatibility

Contiplex® C needle NRFit®: The Device has not been evaluated for safety in the MR environment, since it is not expected to enter the MR Environment in a typical application. Performing an MR exam on a person who has this medical device may result in injury or device malfunction.

Perifix® Catheter Connector NRFit® / Perifix® Filter NRFit® / Perifix® PinPad: MR safe.

Contiplex® ONE Catheter: MR conditional (Please check Operating Instructions section for additional conditions).

NOTE: Components of the set, once connected, take the highest level of MR compatibility assessment (e.g., the Perifix® Catheter Connector NRFit® - Contiplex® ONE Catheter system is considered to be MR conditional).

Sterility

Sterilized by EO.

Do not use if package is damaged.

Visually inspect the sterile barrier system for breaches and integrity prior to use.

Intended purpose

Continuous peripheral nerve block catheter set.

Indication

Plexus and peripheral nerve blocks using electrical nerve stimulation as block anesthesia for surgical procedures (upper or lower extremities), as well as pain treatment, analgesia for early mobilization and physiotherapy or acute and continuous analgesia (e.g. cancer diseases). The Contiplex® C system is indicated for peripheral nerve blocks until a puncture depth of 6 cm.

Patient population

All patients in need for anesthesia or pain therapy by plexus or peripheral nerve blocks. If, despite adequate explanation, the patient refuses plexus anesthesia, their wishes should be respected. Likewise, mentally handicapped patients and those with psychiatric problems need careful preoperative assessment.

Intended user

The product should only be used by healthcare workers that have been adequately trained in this technique.

Contraindications

- Known pre-existing neurological diseases
- Infection / inflammation at the site of insertion/puncturing
- Excessive preoperative anxiety (not responsive to interventions) and inability to tolerate positioning
- Pre-existing coagulopathies, either endogenous (e.g. thrombocytopenia) or iatrogenic (e.g. warfarin therapy) due to increased risk of hematoma
- Anatomic anomalies that make needle placement difficult via ultrasound, nerve stimulation and anatomical landmarks
- Hepatic disease that may interfere with clearance
- Failure to obtain patient's consent
- Known hypersensitivity to one or more materials used
- Contraindications for axillary blocks: lymphangitis and lymphadenectomies
- Stimulation is not to be carried out during defibrillation or cardioversion
- For information on other contraindications against providing continuous regional anesthesia, please refer to the specialist literature

Residual Risks/Side effects

- Nervous lesions
- Pneumothorax
- Pain and hematoma at the site of insertion of the needle
- Incomplete block
- Local anesthetics related toxicity: cardiac, neurologic, allergic
- Incorrect epidural or intrathecal injection of local anesthetic in the case of peripheral blocks in the area of the spinal column
- Infection (sepsis)

Precautions

- Visually inspect the sterile barrier system for breaches and integrity prior to use.
- Use strictly aseptic techniques.
- Do not use if package is damaged.
- When particularly large volumes of local anesthetic are being injected, aspiration tests should be carried out intermittently to ensure that, in the meantime, an incorrect intravascular placement has not accidentally occurred.
- Monitor all connections of the catheter set-up for tightness and disconnection during whole anesthesia and analgesia procedure (e.g., pain treatment).
- Only use with Stimuplex® HNS11 or HNS12 nerve stimulators.

Warnings

- Single use only. Do not re-use. Do not resterilize.
- Re-use of single-use devices creates a potential risk of patient or user. It may lead to contamination and/or impairment of functional capability. Contamination and/or limited functionality of the device may lead to injury, illness or death of the patient.
- If paresthesia is caused inadvertently as a result of direct contact with the nerve, the stimulation needle must not, under no circumstance, be inserted any further.
- An unusually high injection pressure may indicate that an injection has been administered intraneurally. Discontinue the injection, otherwise irreversible nerve damage may occur.
- Products with NRFit® connector are not compatible with Luer devices. They are only suitable for products with NRFit® connector.
- Do not wet Perifix® Filter NRFit® 0.2 µm with disinfectants containing alcohol as this may result in cracks in the housing.
- Do not use syringes smaller than 10 ml as they may cause the filter housing to be damaged as a result of excessive pressure from the syringe.
- Never pull the catheter through the needle as it may otherwise shear off!
- The catheter must be inserted into the connector as far as possible. Injection is impossible otherwise and the catheter can become disconnected.
- Do not use any adhesive spray which contains organic solvents to fix the catheter.
- In order to achieve the optimum connection between the stimulation needle and the nerve stimulator, thereby ensuring patient safety, only use B. Braun products with other B. Braun products. B. Braun Melsungen AG does not assume any liability where its products are combined with equipment from other manufacturers.
- Use the properly stored product by the date stated on the packaging.
- Contiplex® ONE Catheter: MR conditional (Please check Operating Instructions section for additional conditions).

Operating instructions

- Disinfect block area thoroughly. Repeat this at least two times. Allow time to act between each application and cover with a sterile fenestrated drape.
- Position the patient and inject a skin wheal of local anesthetic at the puncture site and if necessary, infiltrate the puncture track.
- Connect the nerve stimulator's electrode cable to the stimulation needle and the skin electrode.
- Set the initial nerve stimulator parameters considering e.g. the respective block type and patient condition (e.g. diabetes).
- Prime needle with solution.
- Disconnect the C-Grip from the lock connection of the needle hub and forward the C-Grip toward the end of the needle at 2 cm mark. Fix the C-Grip at 2 cm by pressing on the C-Grip clip and perform skin penetration. To advance Contiplex® C, release pressure from the clip, push back until 2 cm from skin level, fixate C-Grip again and slowly advance towards the target structure. The distance between positioning clip (C-Grip) and skin level must not exceed 2 cm to ensure optimal stability.
- Hold the stimulation needle by the hub and slowly advance towards the target structure.
- If clear, typical muscle contractions are felt in the target nerve's area of distribution the current is reduced until the muscle contractions are hardly visible.
- Reduce the current and optimize needle position until muscle contractions occur at lower current levels. The tip of the needle has reached the target position when noticeable contractions occur at a current of approximately 0.2–0.5 mA (the exact current depend on e.g. block type and patient condition (e.g. diabetes)).
- The stimulation needle is now situated at appropriate distance from the nerve.
- The stimulation needle should not be placed any nearer the nerve to avoid nerve lesions.
- Aspirate for possible intravascular placement.
- After a negative aspiration control result via the stimulation needle to exclude any incorrect intravascular placement, a bolus of local anesthetic is injected via the stimulation needle. The position can be also controlled under ultrasonographic visualization.
- Once the local anesthetic has been injected, the needle cannot be repositioned as it is usually no longer possible to get the nerves stimulated and the patient no longer experiences pain in the event of the needle's inadvertent intraneural placement.
- To place the catheter, use the C-Grip to fixate the catheter at up to 2 cm above skin level and withdraw the needle from the catheter lumen by pulling the needle at the needle hub. The needle glides out smoothly due to the coating on the needle shaft.

To advance the catheter over the needle tip, use one hand to hold the stimulation needle on the hub in position and forward the catheter using the C-Grip until desired endpoint. The graduation on the steel needle indicates up until 5 cm how much catheter is advanced over the needle tip.

☉ Illustrations on page 3.

- Open the catheter connector lid as far as possible. [Fig. 1]
- Insert the catheter into the funnel shaped opening on the catheter connector and advance it until a clear stop can be felt. [Fig. 2]
- **NOTE:** The catheter must be inserted into the connector as far as possible in order to avoid any leakage and/or disconnection.
- Close the catheter connector's lid by pressing firmly on the top. The lid closes with an audible click. [Fig. 3]
- After removal of the protective cap attach the connector to the Perifix® Filter with a twist action. [Fig. 4]
- Secure the filter to the patient using the B. Braun Perifix® PinPad.
- If necessary: Re-opening the catheter connector: Hold the catheter connector firmly with two fingers and press down hard on the top cover side at the end of the catheter connector. The cover pops up with an audible clicking sound. [Fig. 5]
- Close the injection system with a sterile cap and fixate the catheter on the patient skin (e.g. with the Perifix® Catheter Fixation).

Contiplex® ONE Catheter – A patient with this device can be safely scanned in an MR system meeting the following conditions:

- Static magnetic field of 1.5 Tesla and 3 Tesla, **with**
 - Maximum spatial field gradient of 2,500 G/cm (25 T/m)
 - Maximum force product of 50,000,000 G²/cm (50 T²/m)
 - Maximum whole body averaged (WBA) specific absorption rate (SAR) of 2W/kg (Normal Operating Mode) without time limit in continuous scanning.

Duration of use

The stimulation needle must only be used for detecting the target nerve and, if necessary, for injecting an initial dose of local anesthetic. It must then be removed. The catheter has a possible indwelling period of up to 30 days.

Disposal

Disposal according to local guidelines and clinical protocols.

Storage and handling conditions

Please refer to information on packaging.

Keep away from sunlight.

Keep dry.

Notice to the user

If, during the use of this product or as a result of its use, a serious incident has occurred, please report it to the manufacturer and/or its authorized representative and to your national authority.

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B.Braun (Thailand) Ltd.,