
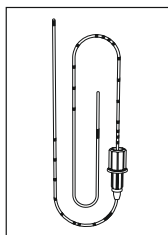
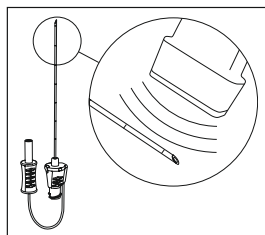
Contiplex® S Ultra 360®

NRFit®


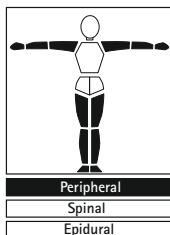
en	Continuous peripheral nerve block catheter set with ultrasound-visible stimulation needles	4
de	Katheterset für kontinuierliche periphere Nervenblockaden mit im Ultraschall sichtbaren Stimulationskanülen	6
bg	Комплект катетри за продължително блокиране на периферните нерви с видими под ултразвук игли за стимулация	8
cs	Souprava katétru na kontinuální blokádu periferních nervů se stimulačními jehlami viditelnými pod ultrazvukem	10
da	Katetersæt til kontinuerlig perifer nerveblokade med ultralydssynlige stimuleringsnåle	12
el	Σετ καθετήρα για το συνεχή αποκλεισμό νευρών με υπερηχογραφικά ορατές βελόνες διέγερσης	14
es	Set de catéter para bloqueo nervioso periférico continuo con agujas de estimulación por ultrasonidos visible	16
et	Püsiva perifeerse närviblokaadi kateetrikomplekt koos ultrahelikontrastsete stimulsiooninõeltega	18
fi	Jatkuvatoiminen ääreishermostojen puudutuskatetrisarja ultraäänessä näkyvillä stimulointineuloilla	20
fr	Cathéter de bloc nerveux périphérique continu équipé d'aiguilles de stimulation visibles aux ultrasons	22
hr	Kateter set za kontinuiranu blokadu perifernog živca s iglama za stimulaciju vidljivim putem ultrazvuka	24
hu	Katéterkészlet folyamatos perifériás idegblokádhöz ultrahangon látható stimulációs tűkkel	26
it	Catetere per il blocco periferico continuo del nervo applicato con aghi stimolanti visibili agli ultrasuoni	28
ko	초음파 가시 자극 바늘이 있는 연속 말초 신경 차단 카테터 세트	30
lt	Nuolatinio periferinių nervų blokavimo kateterio rinkinys su ultragarsu matomomis stimuliuojimo adatomis	32
lv	Ilgstošas periferō nervu bloķēšanas katetru komplekts ar ultraskaņā saskatāmām stimulēšanas adatām	34
nl	Kateterbestek voor continue zenuwblokade met ultrasoon-waarneembare stimulatienaalden	36
no	Katetersett med nåler til stimulering som er synlige på ultralyd, til kontinuerlige, perifer nerveblokader	38
pl	Zestaw cewnika z igłami do wykonywania ciągłych blokad nerwów obwodowych pod kontrolą USG	40
pt	BR Catéter de bloqueio do nervo periférico contínuo definido com agulhas de estimulação de ultrassom-visível	42
pt	PT Conjunto de cateter para bloqueio de nervo periférico contínuo com agulhas de estimulação visíveis por ultrassom	44
ro	Set cateter pentru bloc nervos periferic continuu cu ace de neurostimulare cu vizibilitate ecografică	46
ru	Набор катетеров с экзогенными стимуляционными иглами для продленной блокады периферического нерва	48
sk	Súprava katétra na kontinuálnu blokádu periférneho nervu so stimulačnými ihlami viditeľnými pod ultrazvukom	50
sl	Katetrski sistem za kontinuirano blokado perifernih živcev z ultrazvočno vidnimi stimulacijskimi iglami	52
sv	Komplet kateterna för kontinualno blockering periferinier nerava sa stimulationsnålarna i iglarna vid ultraljud	54
th	ชุดสายสวนสำหรับการระงับความรู้สึกประสาทส่วนปลายแบบต่อเนื่อง พร้อมเข็มกระตุ้นชนิดมองเห็นด้วยอัลตราซาวด์ได้	58
tr	Ultrasonda görülebilir stimülasyon iğneleri sürekli periferik sinir bloğu kateter seti	60
vi	Bộ kim gây tê thần kinh ngoại biên liên tục với kim kích thích có khả năng hiển thị tốt dưới siêu âm	62
zh	带超声可视刺激针的连续外周神经阻断导管套装	64

B | BRAUN

en Line drawings of the individual components
de Strichzeichnungen der Komponenten
bg Скици на отделните съставни части
cs Perokresby jednotlivých částí
da Stregtegninger af de enkelte dele
el Γραμμικά σχέδια των ξεχωριστών στοιχείων
es Esquemas de los componentes individuales
et Üksikute komponentide skeemid
fi Osien merkit
fr Croquis des différents constituants
hr Prikaz pojedinih komponenti






hu Az egyes alkotórészek vonalas rajzai
it Schemi dei singoli componenti
ko 개별 구성품 라인 그림
lt Kiekviena detalė pažymėta brūkšniu kodu
lv Atsevišķo materiāla sastāvdaļu svītrkodi
nl Tekeningen van de verschillende onderdelen
no Strektegninger av de enkelte bestanddeler
pl Rysunki pojedynczych części składowych
pt BR Traços da marca de cada componente
pt PT Traços da marca de cada componente
ro Desene liniare ale componentelor

ru Рисунки отдельных компонентов
sk Perokresby jednotlivých častí
sl Shematski prikazi posameznih sestavnih delov
sr Skice pojedinačnih komponenti
sv Bilder på de enskilda komponenterna
th ภาพส่วนประกอบต่าง ๆ
tr Komponentlerin ayrı çizimleri
vi Bản vẽ nét của từng thành phần riêng lẻ
zh 单独组件的图示

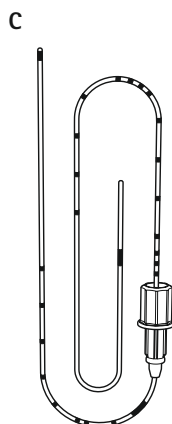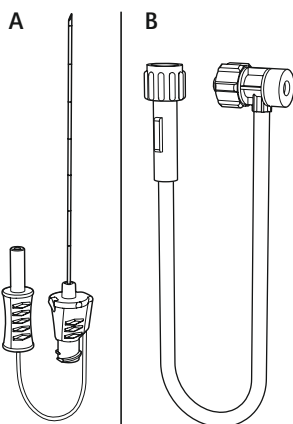Kanülenlänge Needle length	Abstand Spitze -> Stop Block Distance Catheter tip -> Stop Block
Contiplex S	Measure
50 mm	106 mm
100 mm	156 mm
150 mm	206 mm






Catheter-Ø	Measure
0,84 mm	29,5 ± 0,5 mm

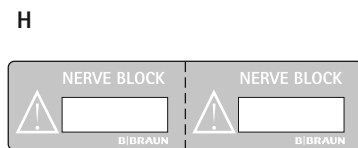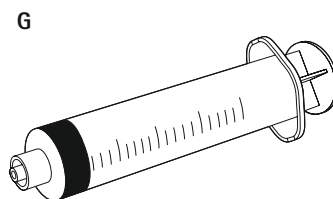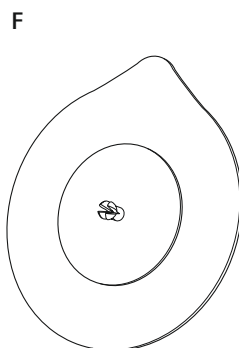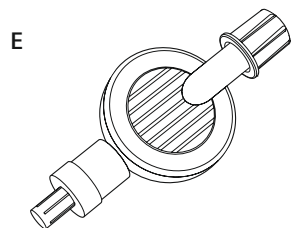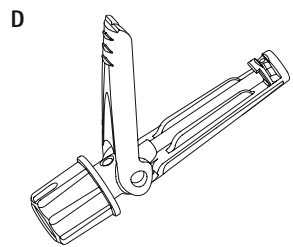en Directions
de Anwendungshinweise
bg Указания за употреба
cs Pokyny pro použití
da Betjeningsvejledning
el Υποδείξεις εφαρμογής
es Aplicaciones
et Kasutusjuhendid
fi Käyttöohje
fr Conseils d'utilisation
hr Upute

hu Használati útmutatás
it Modalità d'uso
ko 사용 방법
lt Įspėjimai
lv Norādījumi par lietošanu
nl Gebruiksaanwijzingen
no Bruk
pl Przepisy stosowania
pt BR Instruções de uso
pt PT Instruções de uso
ro Instrucțiuni de folosire

ru Указания по применению
sk Aplikčné pokyny
sl Napotki za uporabo
sr Uputstva
sv Användningsanvisningar
th วิธีการใช้
tr Kullanım talimatları
vi Chỉ dẫn
zh 说明

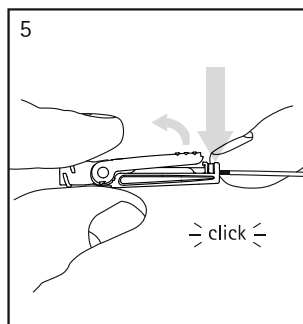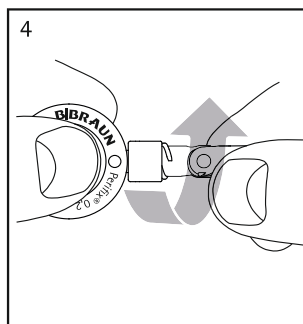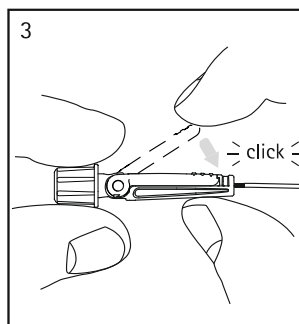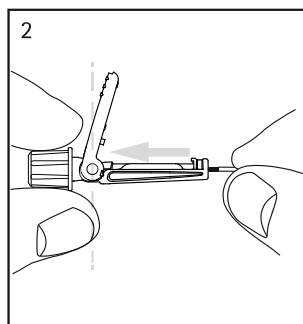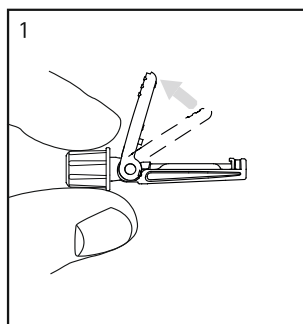en Instructions for Use

Contiplex® S Ultra 360® NRFit® Set

- Contiplex® S Ultra 360® needle NRFit® with integrated wire (Fig. A)
- Sideport assembly with tubing NRFit® (Fig. B)
- Contiplex® ONE Catheter 20G with threading assist guide (Fig. C)
- Perifix® Catheter Connector NRFit® (Fig. D)
- Perifix® Filter (bacteria filter) NRFit® (Fig. E)
- Perifix® PinPad (Fig. F)
- Omnifix® syringe Type NRFit® 5 ml (Fig. G)
- Catheter label (Fig. H)

For further technical characteristics such as needle and catheter size, see product packaging.

NRFit®-connector

The female NRFit® connector of the Contiplex® devices is a standardized connector. It can be connected to syringes with male NRFit® connector. The syringes can be used for flushing the needle, aspirating and infusing anesthetics.

As risk mitigation measure for continuous plexus and peripheral nerve blocks, the catheter label included in the set mentions the respective Regional Anesthesia application.

Electrode connector

The Contiplex® product range contains a standardized electrode connector. The electrode connector can be connected to a nerve stimulator which is intended to be used for peripheral nerve block procedures. In order to achieve the optimum connection between the stimulation needle and the nerve stimulator, thereby ensuring patient safety, only use B. Braun products with other B. Braun products. B. Braun Melsungen AG does not assume any liability where its products are combined with equipment from other manufacturers.

Assessment of Magnetic Resonance (MR) compatibility

Contiplex® S Ultra 360® needle NRFit®/Omnifix® syringe type NRFit®: The Device has not been evaluated for safety in the MR environment, since it is not expected to enter the MR Environment in a typical application. Performing an MR exam on a person who has this medical device may result in injury or device malfunction.

Perifix® Catheter Connector NRFit®/ Perifix® Filter NRFit®/ Perifix® PinPad: MR safe.

Contiplex® ONE Catheter: MR conditional (Please check Operating Instructions section for additional conditions).

NOTE: Components of the set, once connected, take the highest level of MR compatibility assessment (e.g., the Perifix® Catheter Connector NRFit®- Contiplex® ONE Catheter system is considered to be MR conditional). The device contains Cobalt (Co) as indicated on the packaging.

Following components of this set contain cobalt (Co): Contiplex® S Ultra 360® NRFit® needle. Based on the biological safety evaluation and the risk benefit analysis, no risk for the patient/user is expected and no precautionary measures in this regard are needed.

Sterility

Sterilized by EO.

Do not use if package is damaged.

Visually inspect the sterile barrier system for breaches and integrity prior to use.

Intended purpose

Continuous peripheral nerve block catheter set with ultrasound-visible stimulation needles.

Indication

Plexus and peripheral nerve blocks using electrical nerve stimulation and/or ultrasound guidance as block anesthesia for surgical procedures (upper or lower extremities), as well as pain treatment, analgesia for early mobilization and physiotherapy or acute and continuous analgesia (e.g. cancer diseases).

Patient population

All patients in need for anesthesia or pain therapy by plexus or peripheral nerve blocks. If, despite adequate explanation, the patient refuses plexus anesthesia, their wishes should be respected. Likewise, mentally handicapped patients and those with psychiatric problems need careful preoperative assessment

Intended user

The product should only be used by healthcare workers that have been adequately trained in this technique.

Contraindications

- Known pre-existing neurological diseases
- Infection / inflammation at the site of insertion
- Excessive preoperative anxiety (not responsive to interventions) and inability to tolerate positioning
- Pre-existing coagulopathies, either endogenous (e.g. thrombocytopenia) or iatrogenic (e.g. warfarin therapy) due to increased risk of hematoma
- Anatomic anomalies that make needle placement difficult via ultrasound, nerve stimulation and anatomical landmarks
- Hepatic disease that may interfere with clearance
- Failure to obtain patient's consent
- Known hypersensitivity to one or more materials used
- Contraindications for axillary blocks: lymphangitis and lymphadenectomies
- Stimulation is not to be carried out during defibrillation or cardioversion
- For information on other contraindications against providing continuous regional anesthesia, please refer to the specialist literature

Residual Risks/Side effects

- Nervous lesions
- Pneumothorax
- Pain and hematoma at the site of insertion of the needle
- Incomplete block
- Local anesthetics related toxicity: cardiac, neurologic, allergic
- Incorrect epidural or intrathecal injection of local anesthetic in the case of peripheral blocks in the area of the spinal column
- Infection (sepsis)

Precautions

- Visually inspect the sterile barrier system for breaches and integrity prior to use.
- Use strictly aseptic techniques.
- Do not use if package is damaged.
- When particularly large volumes of local anesthetic are being injected, aspiration tests should be carried out intermittently to ensure that, in the meantime, an incorrect intravascular placement has not accidentally occurred.
- Monitor all connections of the catheter set-up for tightness and disconnection during whole anesthesia and analgesia procedure (e.g., pain treatment).
- Only use with Stimuplex® HNS11 or HNS12 nerve stimulators.

Warnings

- Single use only. Do not re-use. Do not re-sterilize.
- Re-use of single-use devices creates a potential risk of patient or user. It may lead to contamination and/or impairment of functional capability. Contamination and/or limited functionality of the device may lead to injury, illness or death of the patient.
- If paresthesia is caused inadvertently as a result of direct contact with the nerve, the stimulation needle must not, under no circumstance, be inserted any further.
- An unusually high injection pressure may indicate that an injection has been administered intraneurally. Discontinue the injection, otherwise irreversible nerve damage may occur.
- Do not wet Perifix® Filter NRFit® 0.2 µm with disinfectants containing alcohol as this may result in cracks in the housing.
- Do not use syringes smaller than 10 ml as they may cause the filter housing to be damaged as a result of excessive pressure from the syringe.
- Never pull the catheter through the needle as it may otherwise shear off!
- The catheter must be inserted into the connector as far as possible. Injection is impossible otherwise and the catheter can become disconnected.
- Do not use any adhesive spray which contains organic solvents to fix the catheter.
- In order to achieve the optimum connection between the stimulation needle and the nerve stimulator, thereby ensuring patient safety, only use B. Braun products with other B. Braun products. B. Braun

Melsungen AG does not assume any liability where its products are combined with equipment from other manufacturers.

- Use the properly stored product by the date stated on the packaging.
- As this medical device uses an alternative small-bore connector design different from those specified in the ISO 80369 series, there is a possibility that a misconnection can occur between this medical device and a medical device using a different alternative small-bore connector, which can result in a hazardous situation causing harm to the patient. Special measures need be taken by the user to mitigate these reasonable foreseeable risks.
- Products with NRFit® connector are not compatible with luer devices. They are only suitable for products with NRFit® connector.
- Contiplex® ONE Catheter: MR conditional (Please check Operating Instructions section for additional conditions).

Operating instructions

- Disinfect block area thoroughly. Repeat this at least two times. Allow time to act between each application and cover with a sterile fenestrated drape.
- Position the patient and inject a skin wheal of local anesthetic at the puncture site. If necessary, infiltrate the puncture track.
- Connect the nerve stimulator's electrode cable to the stimulation needle and the skin electrode.
- Set the initial nerve stimulator parameters considering e.g. the respective block type and patient condition (e.g. diabetes).
- Attach the sideport assembly to the needle hub.
- Prime tubing and needle with solution.
- Hold the stimulation needle by the hub and slowly advance towards the target structure.
- If clear, typical muscle contractions are felt in the target nerve's area of distribution the current is reduced until the muscle contractions are hardly visible.
- Reduce the current and optimize needle position until muscle contractions occur at lower current levels. The tip of the needle has reached the target position when noticeable contractions occur at a current of approximately 0.2–0.5 mA (the exact current depends on e.g. block type and patient condition (e.g. diabetes)).
- The stimulation needle is now situated an appropriate distance from the nerve.
- The stimulation needle should not be placed any nearer the nerve to avoid nerve lesions.
- If ultrasound is used distinct reflexes of the tip geometry and the specific pattern of the distal shaft area of the needle should be visible and indicate the needle position in relation to the nerve. Please note that the visibility may depend on a number of factors such as proper setting of the ultrasound unit, patient related conditions of

the tissue, artefacts and proper alignment of the US probe and the needle.

- Aspirate for possible intravascular placement.
- Following negative aspiration, a test dose of local anesthetic can be injected through the extension tubing. Muscle contractions should cease within 5 to 10 seconds. Prior to indwelling catheter placement, the desired dose of anesthetic agent can be injected through the tubing of the sideport assembly.
- Once the local anesthetic has been injected, the needle cannot be repositioned as it is usually no longer possible to get the nerves stimulated and the patient no longer experiences pain in the event of the needle's inadvertent intraneural placement.
- Attach the plastic threading assist guide to the sideport assembly. Insert the catheter via the threading assist guide through the sideport assembly and the needle.
- The stop block marking on the catheter indicates that the catheter has reached the needle tip (set up: stimulation needle – sideport assembly – threading assist guide).
- Insert catheter to the desired depth. Hold the catheter in place and slowly remove the stimulation needle and sideport assembly.
- Once the catheter has been positioned, connect the remaining catheter set components.

☞ Illustrations on page 3.

- Open the catheter connector lid as far as possible. [Fig. 1]
- Insert the catheter into the funnel shaped opening on the catheter connector and advance it until a clear stop can be felt. [Fig. 2]

NOTE: The catheter must be inserted into the connector as far as possible in order to avoid any leakage and/or disconnection.

- Close the catheter connector's lid by pressing firmly on the top. The lid closes with an audible click. [Fig. 3]
- After removal of the protective cap attach the connector to the Perifix® Filter with a twist action. [Fig. 4]
- Secure the filter to the patient using the B. Braun Perifix® PinPad.
- If necessary: Re-opening the catheter connector: Hold the catheter connector firmly with two fingers and press down hard on the top cover side at the end of the catheter connector. The cover pops up with an audible clicking sound. [Fig. 5]
- Close the injection system with a sterile cap and fixate the catheter on the patient skin (e.g. with the Perifix® Catheter Fixation).

Contiplex® ONE Catheter – A patient with this device can be safely scanned in an MR system meeting the following conditions:

- Static magnetic field of 1.5 Tesla and 3 Tesla, **with**
 - Maximum spatial field gradient of 2,500 G/cm (25 T/m)

- Maximum force product of 50,000,000 G²/cm (50 T²/m)
- Maximum whole body averaged (WBA) specific absorption rate (SAR) of 2W/kg (Normal Operating Mode) without time limit in continuous scanning.

Duration of use

The stimulation needle must only be used for detecting the target nerve and, if necessary, for injecting an initial dose of local anesthetic. It must then be removed. The catheter has a possible indwelling period of up to 30 days.

Disposal

Disposal according to local guidelines and clinical protocols.

Storage and handling conditions

Please refer to information on packaging.

Keep away from sunlight.

Keep dry.

Notice to the user

If, during the use of this product or as a result of its use, a serious incident has occurred, please report it to the manufacturer and/or its authorized representative and to your national authority.

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B.Braun (Thailand) Ltd.,