
Perican® / Perican® Paed

Epidural Needle for Single Shot Technique and
Epidural Catheter Insertion

Epican® Paed

Epidural Needle for Pediatric Caudal Anesthesia (Single Shot)

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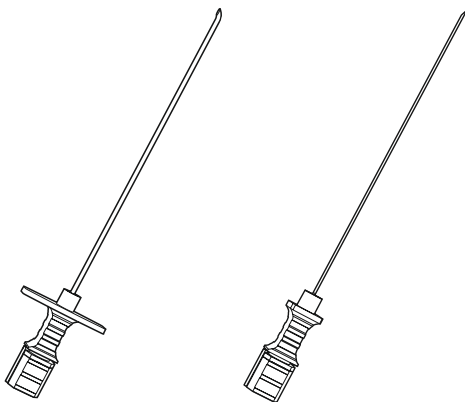en Instructions for use

Perican® and Epican® epidural needles (Luer versions) can be used with:

- Products with Luer connectors (e.g. syringes)
- Regional Anesthesia catheters and threading assist guide (only valid for Perican®)

Epidural needles can be used with devices which are compatible or in conformance with Luer connectors (e.g. ISO 80369-7).

Assessment of Magnetic Resonance (MR) compatibility
Perican®/Perican® Paed/Epican® Paed epidural needles: The device has not been evaluated for safety in the MR environment, since it is not expected to enter the MR environment in a typical application. Performing an MR exam on a person who has this medical device may result in injury or device malfunction.

The device contains Cobalt (Co) as indicated on the packaging.

Based on the biological safety evaluation and the risk benefit analysis, no risk for the patient/user is expected and no precautionary measures in this regard are needed.

Sterility

Sterilized using ethylene oxide.

Do not use if package is damaged.

Visually inspect the sterile barrier system for breaches and integrity prior to use.

Intended purpose

Perican® / Perican® Paed:

Epidural needle for Single Shot Technique and Epidural Catheter Insertion

Epican® Paed:

Epidural needle for Pediatric Caudal Anesthesia (Single Shot)

Indication

Used for epidural anesthesia performed in adults (lumbar, thoracic, cervical) and in pediatrics (caudal, lumbar, thoracic) during surgery and epidural analgesia for regional pain treatment. If necessary, anesthesia can be prolonged during surgery. Analgesia is possible at any time after surgery (pain treatment).

Patient population

- Adults
- Pregnant women
- Pediatric patients

Intended User

The products should only be used by healthcare workers that have been adequately trained in this technique.

Contraindications

Absolute contraindications include patient refusal, severe uncorrected hypovolemia, increased intracranial pressure, infection at the site of injection and known hypersensitivity to local anesthetics.

Relative contraindications include coagulation disorders, anticoagulant therapy, anatomical abnormalities of the vertebral column, chronic low back pain, neurological diseases (e.g. nerve lesions, polyneuropathy, multiple sclerosis), cardiovascular

diseases (inability to increase cardiac output) and sepsis.

Additional contraindications for cervical epidural anesthesia are sympathetic blockade in patients with congestive heart failure or severe aortic or mitral stenosis or congestive lung disease.

Do not perform cervical epidural anesthesia for pediatric anesthesia!

Please consult medical textbooks for standard contraindications to performing epidural anesthesia procedures.

Residual Risks/Side effects

Known risks of epidural anesthesia techniques include, in particular, hypotension, dyspnea/apnea, cardiac arrest, bradycardia, hypothermia, urinary retention, toxic effects of local anesthetic/adjunct (please refer to the Summary of Product Characteristics (SPC) of the drug), inadequate block, temporary or permanent neurological complications and local or systemic infections (postdural puncture headache (PDPH), paresthesia, meningitis/cauda equine syndrome, transverse myelitis or anterior spinal artery syndrome with permanent paralysis, epidural hematoma/epidural space infection, abscess, skin and subcutaneous tissue infection), inadvertent puncture of blood vessels.

Additional risks for caudal anesthesia (technique of choice for epidural catheter placement in babies and children) are back pain, intraosseous injections with subsequent pelvic damage and venous air embolism,

and in very rare cases the development of epidermoid tumors.

Additional risks for cervical epidural anesthesia include spinal cord damage, T1 – T4 sympathectomy, phrenic nerve paresis or paralysis, decrease in heart rate by blocking sympathetic heart innervations and modulation of the vagal nerve and blood flow reduction in the vertebrobasilar artery system.

In case of catheter placement:

When applying the **Catheter technique** please refer to the corresponding instructions for use of Perifix® provided with the catheter.

Precautions

- Trained staff and sufficient emergency equipment, in particular artificial ventilation and hemodynamic support, must be nearby.
- During and after cervical epidural anesthesia, cardiovascular and respiratory monitoring is important to exclude cardio-pulmonary insufficiency.
- When applying the catheter technique, please refer to the Instructions for use of the Perifix® set provided with the catheter.
- Use strictly aseptic techniques.
- Do not use if package is damaged.
- Do not use if product is incomplete or damaged.
- Visually inspect the sterile barrier system for breaches and integrity prior to use.

Warning

- Single use only. Do not re-use. Do not re-sterilize.
- Re-use of single-use devices creates a potential risk of patient or

user. It may lead to contamination and/or impairment of functional capability. Contamination and/or limited functionality of the device may lead to injury, illness or death of the patient.

- If used with an epidural catheter (catheter technique), never pull the catheter back through the needle as it may shear off!
- The tip of the needle becomes blunt or can even split after bone contact. The needle might bend and eventually break when applying too strong force during advancement. Needles with damaged tip increase the risk of intrathecal and intravascular placement.
- If any resistance is felt during advancement of the epidural needle carefully correct the orientation of the needle but never apply strong forces in order to overcome the obstacles.
- In order to achieve the optimum connection between catheter and epidural needle, thereby ensuring patient safety, only use B. Braun products with other B. Braun products. B. Braun Melsungen AG does not assume any liability where its products are combined with equipment from other manufacturers.

Operating Instructions

Directions for using Perican® / Perican® Paed / Epican® Paed epidural needles

1. Detailed information of the patient is mandatory. Optimal positioning of the patient during puncture and injection. Lumbar

epidural anesthesia can be carried out with the patient being in a lateral decubitus (preferred) or sitting position. For both positions it is important to minimize lumbar lordosis and to identify the midline. Check emergency equipment (intubation kit, drugs) and ensure continuous monitoring of the patient.

2. Identification of the landmarks and the approach area.
3. Skin disinfection, local anesthesia of the injection channel.
4. Introduction of the epidural needle (Tuohy needle = Perican®, Crawford needle = Epican®), including the stylet, through the "median" or "lateral" approach. After passing the interspinous ligament, the needle must be advanced carefully, millimeter by millimeter, in the direction of the ligamentum flavum.
5. The stylet is removed and a low friction syringe (Loss of Resistance; LOR - syringe) is attached for the identification of the epidural space with the "Loss of Resistance technique" (LOR). Two different LOR techniques are used, either with air or with sodium saline 0,9 %.
6. Aspiration test.
7. Test dose: After negative aspiration test, a test dose of the local anesthetic drug is recommended.
8. Waiting period: If there is no anesthetic effect after 5–10 minutes, the main dose of the local anesthetic is injected as described in the manufacturer's instructions (**Single Shot technique**) or

- the catheter is inserted into the epidural space (**Catheter technique**). For continuous catheter technique please see corresponding instructions for use (Perifix®).
9. After complete injection of the drug withdraw the needle and cover the puncture site with sterile plaster (**Single Shot technique**). When applying the **Catheter technique** please refer to the corresponding instructions for use of Perifix® provided with the catheter.
10. The level of local anesthesia is to be carefully monitored by e.g. cold devices. After reaching complete block of the motoric and/or sensoric and/or sympathetic nerves in the desired dermatome, surgery can start.

Duration of use

Epidural needles are used for a few minutes for puncture or for the injection of drugs into the epidural space or epidural catheter placement.

The needles are withdrawn from the patient immediately after completion of the procedure.

Disposal

Disposal according to local guidelines and clinical protocols.

Storage and handling conditions

Keep away from sunlight

Keep dry

Clinical benefits & Performances

The summary of safety and clinical performance is available in the European database on medical devices (EUDAMED) at ec.europa.eu/tools/eudamed using the following Basic

UDI-DI: 40392390000023912R for Perican® / Perican® Paed, 40392390000027953M for Epican® Paed and can be requested to the Manufacturer.

Notice to the user

All serious incidents related to this product must be reported to B. Braun and the competent authority in the country where the product is being used.

If further Instructions for Use are required, they can be requested from the manufacturer or can be obtained from the B. Braun homepage: <https://eifu.bbraun.com/>.

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AU For sterilization type see
primary packaging
B. Braun Australia Pty. Ltd.
Level 5, 7 – 9 Irvine Place
Bella Vista NSW 2153
Australia

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CA Distributed by:
B.Braun Medical Inc.
Bethlehem
PA 18018-3524

HN AGENCIA MATAMOROS S.
de R.L

ID Distributed by:
PT. B. Braun Medical
Indonesia
Jakarta, Indonesia

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MY Authorised representative:
B. Braun Medical Industries
Sdn. Bhd.
Bayan Lepas Free Industrial
Zone, 11900 Penang, Malaysia.
www.bbraun.com

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B. Braun Medical AG