
Pencan / Spinocan®

NRFit®

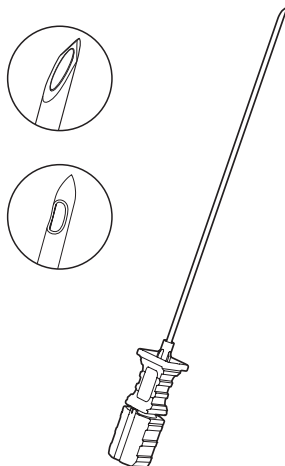B | BRAUN

LLDorder 7334 - grunewald 231433

en	Spinal needle for intrathecal injections and diagnostic lumbar puncture	4
de	Spinalkanüle zur intrathekalen Injektion und diagnostischen Lumbalpunktion	6
bg	Спинална игла за интратекални инжекции и диагностична лумбална пункция	9
cs	Spinální jehla pro intratekální injekce a diagnostickou lumbální punkci	11
da	Spinalkanyle til intratekale injektioner og diagnostisk lumbalpunktur	13
el	Βελόνα σφυονωτιαίας παρακέντησης για ενδορραχιαία ένεση και διαγνωστική σφυονωτιαία παρακέντηση	16
es	Aguja espinal para inyecciones intratecales y punción lumbar diagnóstica	19
et	Spinaalnõel intratekaalseteks injektsioonideks ja diagnostilisteks lumbaalpunktsioonideks	21
fi	Spinaalineula intratekaalisia injektioita ja diagnostisia lannepistoja varten	24
fr	Aiguille spinale pour injections intrathécales et ponction lombaire diagnostique	26
hr	Spinalna igla za intratekalne injekcije i dijagnostičku lumbalnu punkciju	29
hu	Spinális tű intratekális injekciókhoz és diagnosztikai célú lumbálpunkcióhoz	31
id	Jarum spinal untuk injeksi intratekal dan pungsi lumbal diagnostik	34
it	Ago spinale per iniezioni intratecali e puntura lombare diagnostica	36
ko	척수강내 주사 및 진단상 요추 천자를 위한 척수침	39
lt	Spinalinė adata injekcijoms į intratekalinę ertmę ir diagnostinėms punkcijoms į juosmeninę sritį	41
lv	Mugurkaula adata intratekālām injekcijām un diagnostiskai lumbālai punkcijai	43

nl	Spinale naald voor intrathecale injecties en diagnostische lumbale puncties	46
no	Spinalnål for intratekale injeksjoner og diagnostisk lumbalpunksjon	48
pl	Igła do wstrzyknięć dooponowych i diagnostycznych nakłuć lędźwiowych	51
pt	Agulha espinal para injeções intratecais e punção lombar de diagnóstico	53
ro	Ac spinal pentru injectii intratecale și punctie lombară în scop de diagnostic	56
ru	Спинальная игла для интратекального введения и диагностической люмбальной пункции	58
sk	Spinálna ihla na intratekálne injekcie a diagnostickú lumbálnu punkciu	61
sl	Spinalna igla za intratekalno injiciranje in diagnostično lumbalno punkcijo	64
sr	Spinalna igla za intratekalne injekcije i dijagnostičku lumbalnu punkciju	66
sv	Spinalnål för intratekala injektioner och diagnostisk lumbalpunktion	68
th	เข็มเจาะไขสันหลังสำหรับการฉีดยาทางช่องไขสันหลังและการเจาะน้ำไขสันหลังเพื่อตรวจวินิจฉัย	71
tr	Intratekal enjeksiyonlar ve tanısal lomber ponksiyon için spinal iğne	73
vi	Kim gây tê tủy sống dùng để tiêm vào nội tủy và chọc dò chẩn đoán vùng thắt lưng	76
zh	TW 用于鞘内注射和诊断性腰椎穿刺的脊椎穿刺针	78

en Instructions for use

Device description

The spinal needle can be attached to e. g. a spinal introducer needle, syringes (e. g. Omnifix®), a three-way-stopcock or other medical devices with a matching connector which are compatible or in conformance with NREFit® connectors (e.g., ISO 80369-6).

Sterility

Sterilized by EO

Do not use if package is damaged
Visually inspect the sterile barrier system for breaches and integrity prior to use

Intended purpose

Spinal needle for intrathecal injections and diagnostic lumbar puncture

Indication

- Lumbar puncture (diagnosis)
- Spinal anesthesia (injection of local anesthetics into the subarachnoid space for operative procedures and pain management)
- Injection of drugs into the subarachnoid space according to their approved Summary of product characteristics

Patient population

Spinal anesthesia and/or diagnostic lumbar puncture are used for adults, pregnant women and pediatric patients.

The patient groups for injection of drugs into the subarachnoid space

have to follow the Summary of product characteristics of the drug.

Intended User

The product should only be used by physicians or other specially qualified Health Care Professionals (e. g. CRNA) that have been adequately trained in the technique of Regional Anesthesia applied.

Contraindications

The recognized contraindications for spinal anesthesia, lumbar puncture and in case of delivery of drugs into the intrathecal space must be observed. The respective Summary of product characteristics and contraindication of the drugs injected have to be considered.

Contraindications are:

- Blood coagulation disorders / endogenous or iatrogenic coagulopathy
- Skin infection at and near the puncture / injection site
- Sepsis
- History of hypersensitivity to the local anesthetic or the drug have to be considered as well as any of the materials employed
- Patient refusal
- Severe untreated or decompensated hypovolemia
- Increased intracranial pressure
- Shock
- Anatomical deformities of the patient's back
- Amniocentesis
- Acute cerebral or spinal cord disease

Residual Risks/Side effects

The known risks of lumbar puncture for diagnosis or injections include:

- Post dural puncture headache
- Low cerebrospinal fluid pressure syndrome
- Neurological disorders such as a consequence of epidural hematoma or abscess formation
- Anterior spinal artery syndrome
- Cauda equina syndrome
- Herniation
- Pain
- Infections
- Mispositioning
- Epidermoid tumor
- Retroperitoneal abscess

In addition, general risks of spinal anesthesia (such as acute toxicity to the local anesthetic, hypotension and bradycardia, hypothermia, nausea and vomiting, total spinal anesthesia, cardiac arrest, apnoe, urinary retention, hearing loss) or the respective therapy applied according to the Summary of product characteristics of the injected drug may occur.

Leakage of cytotoxic or biological fluids can be an additional hazard to patients, visitors or staff.

Spinal anesthesia, spinal injections and lumbar puncture should only be performed with the appropriate equipment and trained staff to manage such complications.

Precautions

- Visually inspect the sterile barrier system for breaches and integrity prior to use.

- The local anesthetic should not be injected without the evidence of cerebrospinal fluid.
- Mind the prescribing information of the applied drug.
- Products with NRFit® connector are not compatible with Luer devices. They are only suitable for products with NRFit® connector!

Warning

- Single use only. Do not re-use. Do not resterilize.
- Re-use of single-use devices creates a potential risk of patient or user. It may lead to contamination and/or impairment of functional capability. Contamination and/or limited functionality of the device may lead to injury, illness or death of the patient.
- Avoid applying excessive force during advancement as the spinal needle might bend and eventually break.
- The tip of the spinal needle may become blunt or deform after bone contact. A needle with a damaged tip increases the risk of post dural puncture headache due to extended perforation of the dura mater.
- Do not use if package is damaged.
- In order to achieve the optimum connection between spinal needle and other medical devices or accessories, thereby ensuring patient safety, only use B. Braun products

with other B. Braun products. B. Braun Melsungen AG does not assume any liability where its products are combined with equipment from other manufacturers.

- MR unsafe.

Operating Instructions

When performing spinal anesthesia/lumbar puncture, aseptic conditions are required in accordance with national guidelines and/or hospital protocols.

1. Orient the patient according to his/her condition and/or the planned puncture/injection site. Choose between lateral decubitus, sitting or prone position. Locate the midline and follow it during needle placement.
2. Select a suitable spinal needle. It is recommended to use an introducer for spinal needles that are 24G or less in diameter. Check emergency equipment (intubation kit, drugs) and ensure continuous monitoring of the patient.
3. Disinfect puncture site.
4. Inject 1-1.5 cc of local anesthetic into the skin at the selected puncture site of the spinal needle to anesthetize the tissue and the supraspinous and intraspinal ligaments.
5. The introducer needle (for needles 24G or smaller) is inserted at the puncture site until reaching the interspinous ligament. The introducer

er is fixed with thumb and index finger of the left hand (for right-hand user).

6. The spinal needle including the stylet is advanced through the introducer needle. The typical "dura click" occurs after the ligamentum flavum has been passed and the subarachnoid space has been reached.
7. Remove the stylet and check for CSF (cerebrospinal fluid) backflow to confirm that needle has punctured the dura mater and is in the subarachnoid space.
8. After confirmation either withdraw spinal fluid in case needed for diagnostic purpose or connect the syringe with the drug to the spinal needle and inject the drug according to the Summary of product characteristics.
9. Pull back spinal needle and cover the insertion site with a dressing.
10. In case of spinal anesthesia validate success of injection (e.g. by checking if patient feels warm legs and loses ability to lift legs and/or by prick test).

Duration of use

Spinal needles and Spinal Introducer needles are intended for transient use for less than 60 minutes. In clinical use the spinal needles are applied only for a few minutes.

Disposal

Disposal acc. to local guidelines and/or clinical protocols.

Storage and handling conditions

Please refer to the information on the packaging.

Keep away from sunlight.

Keep dry.

Clinical benefits & Performances

The summary of safety and clinical performance is available via EUDAMED at ec.europa.eu/tools/eudamed and can be requested to the manufacturer.

Notice to the user

If, during the use of this product or as a result of its use, a serious incident has occurred, please report it to the manufacturer and/or its authorized representative and to your national authority.

If further Instructions for Use are required, they can be requested from the manufacturer or can be obtained from the B. Braun homepage: <https://eifu.bbBraun.com/>.

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B. Braun Aesculap Japan Co. Ltd.