

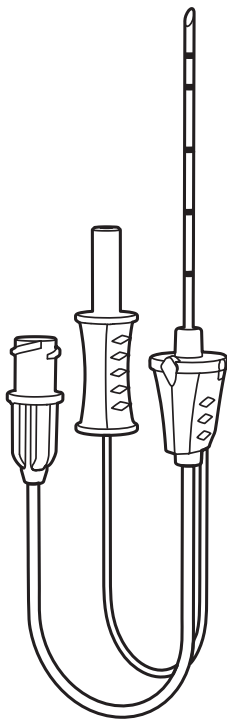
Stimuplex® Ultra 360®

NRFit®

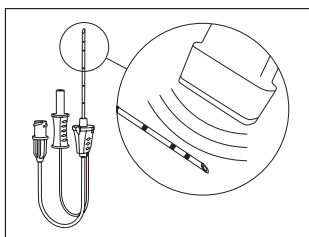B | BRAUN

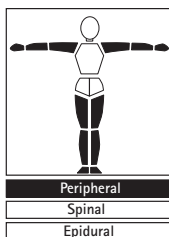en	Ultrasound-visible stimulation needle for peripheral nerve blocks	4
de	Ultraschallsichtbare Stimulationskanüle für periphere Nervenblockaden	6
bg	Канюла за периферна нервна блокада, подходяща за ултразвуково стимулиране	8
cs	Stimulační jehla pro blokádu periferních nervů viditelná pod ultrazvukem	10
da	Ultralydsynlig stimuleringsnål til perifere nerveblokker	12
el	Βελόνα διέγερσης για περιφερειακούς αποκλεισμούς νεύρων, κατάλληλη για χρήση με υπέρηχο	14
es	Cánula de neuroestimulación visible mediante ecografía para bloqueos nerviosos periféricos	16
et	Ultraheliks sobiv perifeersetele närviblokkide stimulatsiooninõel	18
fi	Ultraäänikäyttöön sopiva stimulointineula ääreishermostojen johtopuudutuksiin	20
fr	Aiguille de stimulation visible sous échographie pour blocs nerveux périphériques	22
hr	Stimulacijska igla za izvođenje ultrazvukom vođenih blokada perifernih živaca	24
hu	Ultrahangot visszaverő stimulációs tű perifériás idegblokkoláshoz	26
it	Ago di stimolazione ecogena per blocchi dei nervi periferici	28
lt	Ultragarsu kontroliuojama stimuliavimo adata, skirta periferinių nervų blokadoms	30
lv	Ar ultraskaņu redzama stimulācijas adata perifērā nerva blokādei	32
nl	Stimulatienaald voor perifere zenuwblokkades, zichtbaar bij echografie	34
no	Ultralydegnet stimuleringsnål for perifere nerveblokker	36

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Device description

NRFit®-connector

The female NRRFit® connector of the Stimuplex® Ultra 360® device is a standardized connector. It can be connected to syringes with male NRRFit® connector. The syringes can be used for flushing the needle, aspirating and infusing anesthetics.

Electrode connector

The Stimuplex® Ultra 360® contains a standardized electrode connector. The electrode connector can be connected to a nerve stimulator which is intended to be used for peripheral nerve block procedures. In order to achieve the optimum connection between the stimulation needle and the nerve stimulator, thereby ensuring patient safety, only use B. Braun products with other B. Braun products. B. Braun Melsungen AG does not assume any liability where its products are combined with equipment from other manufacturers.

Assessment of Magnetic Resonance (MR) compatibility

Stimuplex® Ultra 360® Needle NRRFit®: The Device has not been evaluated for safety in the MR environment, since it is not expected to enter the MR Environment in a typical application. Performing an MR exam on a person who has this medical device may result in injury or device malfunction.

For further technical characteristics such as needle size, see product packaging.

Sterility

Sterilized by EO.

Do not use if package is damaged.

Visually inspect the sterile barrier system for breaches and integrity prior to use.

Intended purpose

Ultrasound-visible stimulation needle for peripheral nerve blocks.

Indication

Plexus and peripheral nerve blocks using electrical nerve stimulation and/or ultrasound respectively for precise detection of definite nerve areas. Nerve blocks are used as anesthesia for surgical procedures as well as for pain therapy involving the upper and lower extremities with various block techniques.

Patient population

All patients in need for anesthesia or pain therapy by plexus or peripheral nerve blocks. If, despite adequate explanation, the patient refuses plexus anesthesia, their wishes should

be respected. Likewise, mentally handicapped patients and those with psychiatric problems need careful preoperative assessment.

Intended User

The product should only be used by healthcare workers that have been adequately trained in this technique.

Contraindications

- Known pre-existing neurological diseases
- Infection / inflammation at the desired site of insertion / puncturing
- Excessive preoperative anxiety (not responsive to interventions) and inability to tolerate positioning
- Pre-existing coagulopathies, either endogenous (e.g. thrombocytopenia) or iatrogenic (e.g. warfarin therapy) due to increased risk of hematoma
- Anatomic anomalies that make needle placement difficult via ultrasound, nerve stimulation and anatomical landmarks
- Hepatic disease that may interfere with clearance
- Failure to obtain patient's consent
- Known hypersensitivity to one or more materials used
- Contraindications for axillary blocks: lymphangitis and lymphadenectomies
- Stimulation is not to be carried out during defibrillation or cardioversion
- For information on other contraindications against providing regional anesthesia, please refer to the specialist literature

Residual Risks/Side effects

- Nervous lesions
- Pneumothorax
- Pain and hematoma at the site of insertion of the needle
- Incomplete block
- Local anesthetics related toxicity: cardiac, neurologic, allergic
- Incorrect epidural or intrathecal injection of local anesthetic in the case of peripheral blocks in the area of the spinal column
- Infection (sepsis)

Precautions

- Once the local anesthetic has been injected, the needle cannot be repositioned as it is usually no longer possible to get the nerves stimulated and the patient no longer experiences pain in the event of the needle inadvertent intraneural placement.

WARNING

- Single use only. Do not re-use. Do not resterilize.

- Re-use of single-use devices creates a potential risk of patient or user. It may lead to contamination and/or impairment of functional capability. Contamination and/or limited functionality of the device may lead to injury, illness or death of the patient.
- If paresthesia is caused inadvertently as a result of direct contact with the nerve, the stimulation needle must not, under no circumstance, be inserted any further.
- An unusually high injection pressure may indicate that an injection has been administered intraneurally. Discontinue the injection, otherwise irreversible nerve damage may occur.
- In order to achieve the optimum connection between the stimulation needle and the nerve stimulator, thereby ensuring patient safety, only use B. Braun products with other B. Braun products. B. Braun Melsungen AG does not assume any liability where its products are combined with equipment from other manufacturers.
- Visually inspect the sterile barrier system for breaches and integrity prior to use.
- Do not use if package is damaged.
- Use the properly stored product by the date stated on the packaging.
- Products with NRFit® connector are not compatible with luer devices. They are only suitable for products with NRFit® connector!

Operating Instructions

- Use strictly under aseptic conditions.
- Disinfect block area thoroughly. Repeat this at least two times. Allow time to act between each application and cover with a sterile fenestrated drape.
- Position the patient and inject a skin wheal of local anesthetic at the puncture site and if necessary, infiltrate the puncture track.
- Connect the nerve stimulator's electrode cable to the stimulation needle and the skin electrode.
- Set the initial nerve stimulator parameters considering e.g. the respective block type and patient condition (e.g. diabetes).
- Prime tubing and needle with solution.
- Hold the stimulation needle by the hub and slowly advance towards the target structure.
- If clear, typical muscle contractions are felt in the target nerve's area of distribution the current is reduced until the muscle contractions are hardly visible.
- Reduce the current and optimize needle position until muscle contractions occur at lower current levels. The tip of the needle has reached the target position when noticeable contractions occur at a current of approximately 0.2 – 0.5 mA (the exact current depend

on e.g. block type and patient condition (e.g. diabetes).

- The stimulation needle is now situated an appropriate distance from the nerve.
- The stimulation needle should not be placed any nearer the nerve to avoid nerve lesions.
- If ultrasound is used distinct reflexes of the tip geometry and the specific pattern of the distal shaft area of the needle should be visible and indicate the needle position in relation to the nerve. Please note that the visibility may depend on a number of factors such as proper setting of the ultrasound unit, patient related conditions of the tissue, artefacts and proper alignment of the US probe and the needle.
- Aspirate for possible intravascular placement.
- Following negative aspiration, a test dose of local anesthetic can be injected through the extension tubing.
- Once the local anesthetic has been injected, the needle cannot be repositioned as it is usually no longer possible to get the nerves stimulated and the patient no longer experiences pain in the event of the needle's inadvertent intraneural placement.

Duration of use

Sterile device for single use.

The Stimuplex® Ultra 360® nerve block needle NRFit® must only be used for detecting the target nerve and for injecting the local anesthetic. It must then be removed.

Disposal

Disposal according to local guidelines and clinical protocols.

Storage and handling conditions

Please refer to information on packaging. Keep away from sunlight. Keep dry.

Notice to the user

If, during the use of this product or as a result of its use, a serious incident has occurred, please report it to the manufacturer and/or its authorized representative and to your national authority.

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B. Braun Australia Pty Ltd
B. Braun Aesculap Japan Co., Ltd.