


# Instructions Manual

## Beauty instrument

F825.N.1/2

THERMOTHERAPY

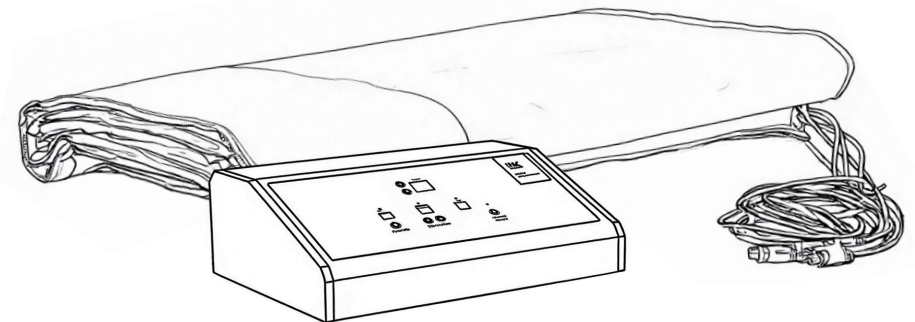

**WK**  
Imported by  
**WEEKKO BARCELONA SL**

Pol. Ind. La Vall d'an  
Camí de Garreta, 92  
08600 Berga  
Barcelona - Spain  
T. +34 93 824 90 20  
E. info@weekko.com

VERSION 1.0

Thank you sincerely for choosing our product. This manual is compiled to help customer service staff use this instrument better. We hope our company's products and services will bring you a distinguished experience.

This user manual contains all general information about the instrument. It is the starting point for new users to begin using the instrument, as it provides a general introduction to the system, explains how it works, and details the instrument's installation and various operation methods.

### Usage Precautions

1. Read the user manual carefully before turning on the machine. Operators must be familiar with the user manual and operation methods, and the instrument must be managed by a dedicated person.
2. Do not use a power strip to avoid fire hazards.
3. Always turn off the power switch when replacing the fuse to prevent fire and electric shock.
4. Must use a power supply of 220V±22V, 50HZ±1HZ and an effectively grounded socket.
5. Using a damaged or modified power cord may cause fire or electric shock.
6. If the instrument is loose, falling, or damaged, please contact the company in a timely manner.
7. Do not use the instrument in unstable environments such as tilting, vibration, or impact.
8. No water or debris should enter the interior of the instrument. In case of foreign objects entering, please contact the company promptly.
9. The power plug must be unplugged before moving the instrument.
10. When moving the instrument, try to avoid vibration.
11. The instrument should be placed on a stable surface with at least 5cm of space around it to ensure air circulation and facilitate heat dissipation of the instrument.
12. Do not cover the instrument with anything to prevent blockage of the air vents.
13. Avoid using the instrument in the following environments: direct sunlight; high humidity environments; near water; vibrating or tilting places; strong magnetic field environments; places with electromagnetic waves or impulse voltage; places where chemicals are stored; places with corrosive gases.


### Product Packaging List

| Name                  | Image Quantity |
|-----------------------|----------------|
| Main Unit             | 1              |
| Fat Reduction Blanket | 1              |
| Power Cord            | 1              |

### Product Specifications

- Input Voltage: AC220V-240V
- Product Frequency: 50/60HZ
- Product Power: 500W
- Packaging Size: 45x33x17CM
- Gross Weight: 5.25KG
- Net Weight: 4.45KG
- 
- Packaging Size: 85x45x12CM
- Gross Weight: 5.75KG
- Net Weight: 4.95KG

## Operation Instructions

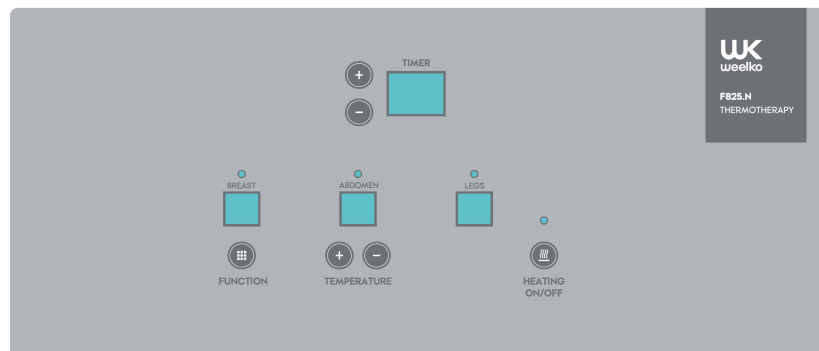

1. Connect the device to a 220V power supply. A beep will sound, indicating that the device is powered on and in standby mode.
2. Press the ON/OFF button to activate the heating function.
3. Working Time Adjustment: The digital tube displays the default setting of 20 minutes upon startup, and the duration can be adjusted between 1–60 minutes.
4. Temperature & Body Part Selection:
  - The default temperature is 25°C upon startup, with a maximum adjustable range of up to 65°C.
  - Press the FUNCTION button to select one of four body part combinations: BREAST, ABDOMEN, LEGS, or all three together (BREAST+ABDOMEN+LEGS).
  - Press the TEMPERATURE button to set the preset temperature for the selected body part(s). After 3 seconds, the active display frame will show the real-time temperature.

14. This instrument shall not be placed or used near radios, televisions, copiers, fax machines, etc.
15. Do not use it together with other instruments (such as microwaves) to avoid mutual interference between instruments and cause misoperation.
16. If an abnormal situation occurs to the instrument or the patient during operation, the operator shall immediately move the operating handle away from the patient, press the stop switch, and stop using the instrument.
17. The instrument shall be correctly used by skilled operators. Personnel other than operators are not allowed to use the instrument.
18. After operation, place the control handle back on the control handle bracket.
19. Clean the probe of the operating handle after use for next use.

## Precautions for Troubleshooting, Storage, and Inspection

1. Do not use the instrument if it malfunctions. In case of smoke, burning smell, or inability to operate, continued use may cause fire or electric shock. Please quickly press the stop button, turn off the power, unplug the power plug, and contact the company or distributor immediately.
2. Personnel other than technical staff are absolutely not allowed to repair or modify the instrument. The interior of the instrument must never be opened.
3. The instrument and its components must be inspected regularly.
4. Do not use cleaning tools that may damage the instrument's protective coating, such as brushes and brooms.
5. If the instrument is not used for more than 2 weeks, it must be carefully inspected before reuse.

### Contraindicated Groups (Do Not Use If You Belong to the Following Groups)

- People with cardiac pacemakers, defibrillators, or metal implants (such as steel plates, bone nails) in the body.
- Pregnant women, lactating women, and women during menstruation.
- People suffering from epilepsy, hypertension, severe heart disease, diabetes, tumors or cancer, or thyroid diseases.
- People with open wounds, damage, severe inflammation, eczema, dermatitis on the skin, or telangiectasia (capillary dilation) or rosacea on the face.
- People who have recently (usually within 6 months) received botulinum toxin injections, filler injections, laser treatment, chemical peeling, or plastic surgery.

### Usage Preparation

1. Check the appearance of the slimming blanket to ensure that the power cord is not damaged, the plug is in good contact, and the thermostat displays normally.
2. Clean and dry the skin of the body part to be used. It is recommended to wear light and breathable cotton clothing to avoid direct contact between the skin and the blanket surface (to prevent low-temperature burns).
3. Lay the slimming blanket flat on a flat and hard surface (such as a yoga mat, bed board) and arrange it neatly to avoid wrinkles.

### Operation Steps

- Firmly connect the thermostat plug to the slimming blanket interface, then insert the power plug into a household socket.
- Press the power button on the thermostat, the indicator light will turn on, and it will enter the standby state.
- It is recommended to set the initial temperature to 40-45℃, and gradually adjust it according to your own tolerance. The maximum temperature shall not exceed 55℃.
- After turning on, it is recommended to use it for 20-30 minutes at a time, 2-3 times a week. Do not use it for a long time every day.
- It can be used with body lotion and slimming essential oil to enhance the skin absorption effect.

- After use, first press the power button on the thermostat to turn it off. After the indicator light goes out, unplug the power plug and the blanket interface.
- Wait for the slimming blanket to cool down to room temperature naturally, fold it neatly, put it into the storage bag, and store it in a cool and dry place.

### Additional Precautions

1. Contraindicated Groups: Pregnant women, infants and young children, diabetics, people with sensitive or damaged skin, and people with chronic diseases such as hypertension and heart disease are prohibited from using this product.
2. Safety Warnings: It is strictly prohibited to use it while sleeping. Avoid long-term high-temperature hot compress on the same part to prevent low-temperature burns. If skin redness, stinging or other discomfort occurs during use, turn off the machine and stop using it immediately.
3. Maintenance: The slimming blanket cannot be washed with water or exposed to the sun. Surface stains can be gently wiped with a damp cloth. Avoid bending or pulling the power cord to prevent damage.

### V. Common Troubleshooting

| Fault Phenomenon                   | Possible Cause                                               | Possible Cause                                                    |
|------------------------------------|--------------------------------------------------------------|-------------------------------------------------------------------|
| No response when turning on        | Power plug is not inserted firmly or the socket has no power | Check the plug connection and try replacing with a powered socket |
| Slow temperature rise              | Wrinkles on the blanket surface or loose wrapping            | Arrange the blanket surface to ensure it fits the body            |
| Abnormal display of the thermostat | Poor interface contact                                       | Unplug the thermostat plug and reinsert it firmly                 |