


Instructions Manual

Beauty instrument

F833.N

BODY MASSAGER

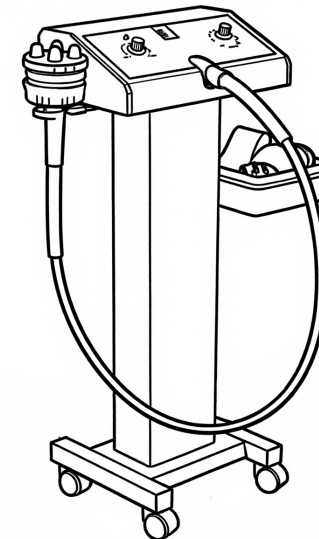

VERSION 1.0

WK
Imported by
WEEKO BARCELONA SL

Pol. Ind. La Vall d'an
Camí de Garreta, 92
08600 Berga
Barcelona - Spain
T. +34 93 824 90 20
E. info@weeko.com

Thank you sincerely for choosing our product. This manual is compiled to help customers use this instrument better. We hope our company's products and services will bring you a distinguished experience.

This user manual contains all general information about the instrument. It is the starting point for new users to begin using the instrument, as it provides a general introduction to the system, explains how it works, and details the instrument's installation and various operation methods.

Usage Precautions

1. Before turning on the machine, be sure to read the user manual carefully. Operators must be familiar with the user manual and operation methods, and the instrument must be managed by a designated person;
2. Do not use a power strip to avoid fire hazards;
3. When replacing the fuse, always turn off the power switch to prevent fire and electric shock;
4. Must use a power supply of 220V±22V, 50HZ±1HZ and an effectively grounded socket;
5. Using a damaged or modified power cord may cause fire and electric shock hazards;
6. If the instrument is loose, falling, or damaged, please contact the company in a timely manner;
7. Do not use the instrument in unstable environments such as tilting, vibrating, or impacting conditions;
8. No water or debris is allowed to enter the interior of the instrument. In case of foreign objects entering, please contact the company promptly;
9. When moving the instrument, the power plug must be pulled out first;
10. When moving the instrument, try to avoid vibration;
11. The instrument should be placed on a stable surface, with at least 5cm of space left around it to ensure air circulation and facilitate heat dissipation of the instrument;
12. Do not cover the instrument with anything to prevent blockage of the air vents;
13. Avoid using the instrument in the following environments: direct sunlight; high humidity environments; near water; vibrating or tilting places; strong magnetic field environments; places where electromagnetic waves or impulse voltages occur; places where chemicals are stored or corrosive gases exist;

Accessory List

Name	Image Quantity
Main Unit	1 pc
Vertical Column	1 pc
Column-type Massage Head	1 pc
Shower-type Massage Head	1 pc
Round Sponge Massage Head	1 pc
Half-ball Type Massage Head	1 pc
U-shaped Sponge	1 pc
Power Cord	1 pc

Product Specifications

- Rated Voltage: AC 220V-240V, 50/60Hz
- Maximum Power Consumption: 75W
- Main Unit Dimensions: 92x39x45CM
- Gross Weight: 16Kg
- Net Weight: 13.5Kg

Notes

1. Be sure to use the instrument according to the correct operation method (as shown in the figure).
2. During operation, the bending angle of the hose shall not be less than 100 degrees. Do not press the hose hard with your hand during operation.
3. Do not excessively pursue high intensity: The appropriate feeling is a comfortable high-frequency vibration on the skin. Higher intensity is not necessarily better.
4. Not a Substitute for Exercise: This is an auxiliary shaping tool. To achieve ideal results, it must be combined with a healthy diet and proper exercise.
5. Pay Attention to Body Signals: If you feel any pain, severe discomfort, or abnormal skin redness during use, stop immediately.
6. Check Instrument Condition: Ensure the instrument is in good condition, especially that the probe or brush head in direct contact with the skin has no cracks.
14. This instrument shall not be placed or used near radios, televisions, copiers, fax machines, etc.;
15. Do not use it together with other instruments (such as microwaves) to avoid mutual interference between instruments and cause misoperation.
16. If an abnormal situation occurs to the instrument or the patient during operation, the operator shall immediately move the operation handle away from the patient, press the stop switch, and stop using the instrument;
17. The instrument shall be correctly used by skilled operators. Personnel other than operators are not allowed to use the instrument.
18. After the operation, place the control handle back on the control handle holder.
19. After using the operation handle, clean the probe for next use.

Precautions for Troubleshooting, Storage, and Inspection

- Do not use the instrument if it malfunctions. In case of smoke, burning smell, or inability to operate, continued use may cause fire or electric shock. Immediately press the stop button, turn off the power, unplug the power plug, and contact the company or distributor immediately;
- Personnel other than professional technicians are strictly prohibited from repairing or modifying the instrument. The interior of the instrument must never be opened;
- The instrument and its components must be inspected regularly;
- Do not use cleaning tools that may damage the instrument's protective coating, such as brushes, brooms, etc.;
- If the instrument is not used for more than 2 weeks, it must be carefully inspected before reuse.

Contraindicated Groups (Do Not Use If You Belong to the Following Groups)

- Pregnant women and lactating women.
- Patients with severe heart disease, hypertension, or epilepsy.
- Individuals with severe skin inflammation, eczema, dermatitis, sunburn, or unhealed wounds.
- Those with implanted cardiac pacemakers or other internal electronic devices.

- Individuals with a history of skin allergies or allergies to metals, silica gel, or cosmetics.

Important Safety and Usage Taboos

- Usage Frequency:** It is recommended to use 3-4 times a week to allow the body time to recover for better results.
- Usage Time:** Avoid using immediately after meals or on an empty stomach. It is recommended to use 1 hour after meals or before bedtime.
- Absolute Contraindicated Groups:** Pregnant women and women during menstruation are prohibited from using; postpartum women (especially those who have undergone cesarean section) are prohibited from using. Individuals with implanted electronic medical devices such as cardiac pacemakers or metal implants are also prohibited from using. Avoid using on areas with excessive acute inflammation, bleeding, skin damage, or malignant lesions. It is not suitable for individuals with thrombophlebitis or lymphatic disorders. Children must use under adult supervision; infants are prohibited from using.
- Do not start the machine when the customer is not ready, otherwise it will increase the wear rate of the vibrating massage head.
- Do not disassemble the centrifugal vibrating head without permission.
- There is no need to add oil to the rotating shaft or the machine. All bearings and rotating shafts are sealed and self-lubricating.
- Clean the exterior of the machine with mild soapy water and lukewarm water. Do not use alcohol or cleaning fluids. Clean accessories in the same way.



Instrument Operation

- Plug in the power supply, and the red indicator light at the top left corner will illuminate.
- Energy Adjustment Knob**  : Adjust the intensity moderately from low to high as needed.
- Timed Start Knob**  : The "0" position indicates the OFF/stop state, and the maximum adjustable timing is 30 minutes.

Operation Precautions for the Vibrating Massage Fat-Removing Machine

Correct Operation Method



Incorrect Operation Method

