
ComfortGlide®

Instructions for Use

- de Gebrauchsanweisung
- fr Mode d'emploi
- nl Gebruiksaanwijzing
- it Istruzioni per l'uso
- es Instrucciones de uso
- pt Instruções para Uso
- sk Návod na použitie
- cs Návod k použití
- sv Instruktioner för användning
- fi Käyttöohjeet
- da Instruktion
- no Bruksanvisning
- el Οδηγίες χρήσης
- pl Instrukcja obsługi
- sl Navodila za uporabo
- ro Instructiuni de folosire
- bg Инструкции за употреба
- et Kasutusjuhised
- hu Használati utasítás
- lt Norādījumi par lietošanu
- hr Upute za uporabu
- tr Kullanim talimatlari

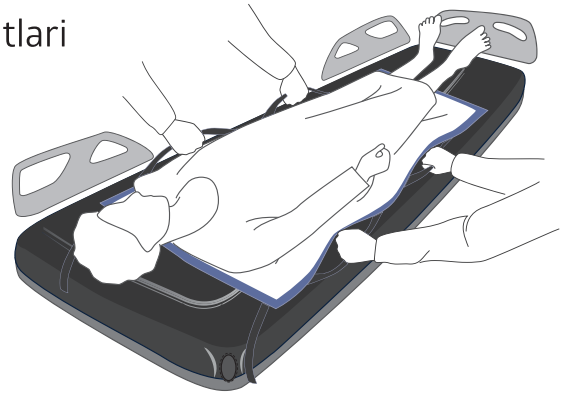REF

CGSHEET
CGSHEETNS
CGSHEETXL
CGSHEETR
CGSHEETRXL
CGWEDGE
CGWEDGEXL
CGDRYPAD
CGDRYPADXL
CGWEDGE-R
CGWEDGE-RXL

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Latex free, Latexfrei, Sans latex, Latexvrij, Senza lattice, Sin látex, Sem látex, Bez latexu, Bez Latexu, Latexfri, Lateksiton, Latexfri, Lateksfri, Χωρίς λάτεξ, Bezlateksowe, Brez lateksa, Fără latex, Несъдържащ латекс, Lateksivaba, Latex mentes, Sudėtyje nėra latekso, Bez lateksa, Lateks içermez

Manufactured in China for: **Medline Industries, LP,** Three Lakes Drive, Northfield, IL 60093 USA

CE EC REP Medline International France SAS,
5 rue Charles Lindbergh, 44110 Châteaubriant, France

CH REP **CH Rep and importer:** Medline International Switzerland SARL, c/o MN & Associés SA,
1 Place de Longemalle, 1204 Genève, Switzerland

UK CA **UK representative and importer in UK:**
Medline Industries Ltd, WA4 6HL, UK

Medline International B.V.,
Nieuwe Stationsstraat 10, Arnhem, 6811 KS, Netherlands

RJ24EMX - 10/2024

en

Instructions for Use:

1. CGSHEET, CGSHEETXL, CGSHEETR, CGSHEETRXL, CGSHEETNS: For Use With: Wedges and Drypads. CGWEDGE, CGWEDGEXL, CGWEDGE-R, CGWEDGE-RXL: For Use With: Sheets and Drypads.
2. Unfold and position the sheet next to patient, align the top of the sheet with the patient's shoulders. The 2 retention straps (if applicable) should be towards the head of the bed.
3. Place a Drypad (if applicable) over the sheet and fold the edges of both long-ways to tuck under the patient.
4. Roll patient onto his/her side following patient handling policies.
5. Tuck the glide sheet under patient's side.
6. On the opposite side of the bed, roll patient away from you and unfold the glide sheet and pad.
7. Make sure to lock the wheels of the bed. Firmly grip the handles and gently slide the patient as needed to center them on the bed.

To Offload the Sacrum:

1. Gently-roll patient onto his/her side following patient handling policies. Place wedges underneath the sheet, with "DOWN" arrows pointing towards the mattress. One wedge should be above the sacrum and the other below.
2. Slowly roll patient back on top of the wedges. Sacrum should be fully off-loaded between both wedges, leaving about a "fist" width of space between them.
3. Check that there are no folds and wrinkles under the patient and smooth out as needed.
4. Adjust the bed to the desired angle and use the dual retention straps (if applicable) to secure the top of the glide sheet to the bed. There are two ways to properly secure the straps.

- Secure the straps to a movable part of the bed frame. Do not attach the straps to any stationary part of the bed. Secure them to an opening in the bed frame, or
- If there is no bed structure within reach, attach the straps to each other around the bottom of the mattress.

To Clean:

Reusable Sheet: 30% nylon, 70% polyester, wash prior to use, wash with like colors, wash warm, rinse well, tumble dry low heat, remove promptly. Do not use bleach or softener. After each washing, check one box of the grid available on the sheet.


Reusable Wedges may be cleaned according to your infection prevention policies. Discard if soiled, damaged, or discolored.


⚠ Caution:

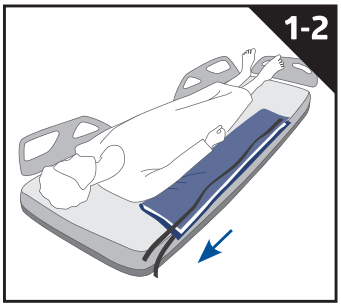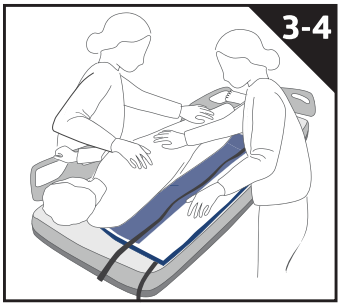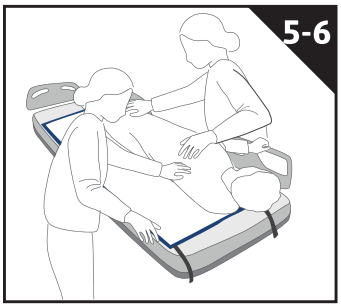
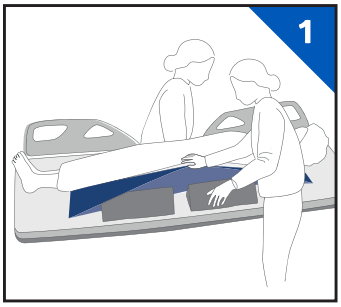
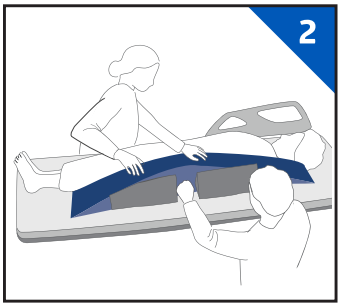
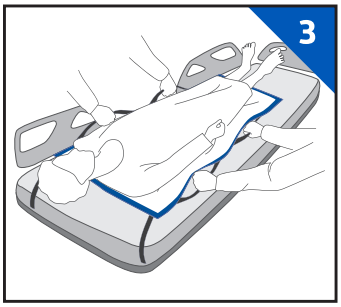
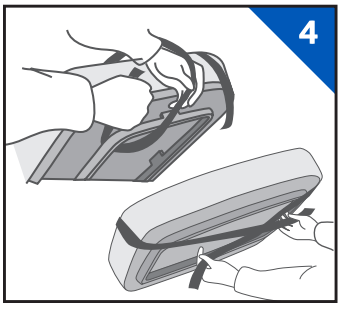
1. Components are not designed for patient lifting.
2. Not to be used in combination with other lifting or positioning devices.
3. If using wedges with Comfort Glide, do not use a sheet to slide the patient onto the wedges. Use log-roll method to place wedges under the patient.
4. Always follow your facility's guidelines for proper patient handling.
5. Two people should be present during use of the repositioning sheet, or bed rails must be in an upward position during reposition.
6. For reusable sheets, examine the sheet prior to use.
7. Reuse of single use devices creates a potential risk for the patient or user. It may cause contamination and/or impairment of function; which may lead to injury, illness and/or death of a patient.
8. If a serious incident has occurred in relation to the device, it should be reported to Medline and the competent authority of the Member State in which the user and/or patient is established.
9. Precautions should be taken when discarding this device and disposal of the device shall be made in accordance with applicable hospital or national regulations for biologically hazardous waste.
10. **CGSHEETR, CGSHEET, CGSHEETNS:** To assist in repositioning patients **up to 261Kg**. **CGSHEETRXL:** To assist in repositioning patients **up to 286 Kg. CGSHEETXL:** To assist in repositioning patients **up to 453Kg**.

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