HUDSON RCI®

# Adaptors and Connectors

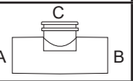
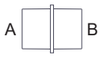
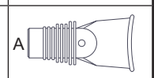
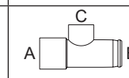
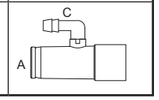
REF **41077, 41409, 41421, 41422, 41565, 41639, 41640, 41642, 41648, 41651, 41659, 41660, 41669**

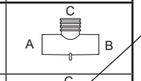
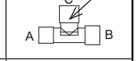
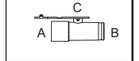
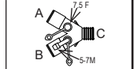
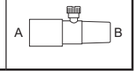
**Rx ONLY**

Teleflex MEDICAL
3015 Carrington Mill Blvd
Morrisville, NC 27560 USA
USA: (866) 246 6990
International: +1 919 544 8000

CE 1639 EC REP Teleflex Medical IDA Business and Technology Park Dublin Rd., Athlone, Co. Westmeath Ireland

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Hudson RCI is a registered trademark of Teleflex Incorporated or its affiliates.

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## Adaptors and Connectors

### Not made with natural rubber latex.

**WARNING:**
- Be sure connections are secure at all times.
- Do not sterilize.
- Where an Adaptor port cap is provided, plug port while not in use to prevent leakage.
- Where an Adaptor port cap is not provided, do not use adaptor without appropriate monitoring equipment. Failure to do so will result in leakage.
- Single patient use: The use of products intended for single patient use on more than one patient creates a potential risk of serious injury and/or infection which may lead to death. Reprocessing of medical devices intended for single patient use only may result in degraded performance or a loss of functionality

**Instructions for Use:**

REF **41659**
**Metered Dose Inhaler Adaptor**
1. Connect the adaptor in-line.
2. Turn the blue arrow in direction medication is to be administered.
3. Insert the canister and dispense medication as prescribed.

REF **41660**
**Universal Y-Piece**
1. Observe the appropriate manufacturer's instructions for circuit installation.
2. If patient pressure sensing line or a temperature probe is not used, place cap on open port.

**Do not use on multiple patients (Single Patient Use Only)**

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© Copyright 2006, 2010, 2011, 2012, 2013, 2016, 2017, 2018 Teleflex L02806 R01