ULTRACOVER™ for TEE
Endocavity Probe Cover
Instructions for Use

en Endocavity Probe Cover

Instructions for Use

Intended Use / Indications for Use
A device intended to be used as a physical barrier for protection against the effects of environmental exposure (e.g., body fluids, gels) and/or to maintain the required hygienic level of various diagnostic or surgical procedures which utilize a transducer (probe).
This is a single use device. [Not to be reprocessed.]
The devices are further indicated to be used in or on:
- Natural body orifices (e.g., vagina, rectum, esophagus)

Intended Patient Population
Probe covers have no limitation for any definition of the patient group and their estimated state/condition. Probe covers can be used at the discretion of the healthcare professional.

Contraindications
There are no known contraindications for these devices.

Warnings
- This product contains Natural Rubber Latex which may cause allergic reactions.
- WARNING: CHOKING HAZARD. Don't turn patient (from stomach to back or reverse) during procedure.
- Do not use if sterile barrier package is damaged. If damage is found, call your Medline-Microtek representative.
- Visually inspect for breaches of packaging integrity (such as tears, cuts, open seals, etc.) prior to use.
- Report serious incidents to legal manufacturer/importer and competent authority.


Precautions
- It is important to spread the gel regularly thoughout the entire length of the cover for easy application and to avoid breaking.
- Secure distal tip to probe with rubber bands or ties to prevent gel from migrating away from tip.
- Use care when handling the Probe Cover in close proximity to sharps, cautery devices, and other medical devices or procedures that may damage the integrity of the Probe Cover.


Product / Device Description
Probe Covers are composed of natural or synthetic materials and come in a variety of sizes and shapes to allow these devices to properly fit equipment with diverse shapes and uses.
Probe Covers are used in a variety of clinical and surgical settings by healthcare professionals. The devices cover a transducer / probe and, in some instances, the cables connected to the transducers / probes. These coverings offer a physical barrier between the transducer and the patient.


How Supplied
- Do not use if package is opened or damaged, store in a cool, dry, dark place.
- This product may be supplied with ultrasound gel, rubber bands and/or ties.

Reuse Warning Statement
For single patient use only. Do not reuse, reprocess or resterilize. Reuse, reprocessing or resterilization may compromise the structural integrity of the device and/or lead to device failure which, in turn, may result in patient injury, illness or death. Reuse, reprocessing or resterilization may also create a risk of contamination of the device and/or cause patient infection or cross-infection, including, but not limited to, the transmission of infectious disease(s) from one patient to another. Contamination of the device may lead to injury, illness or death of the patient.

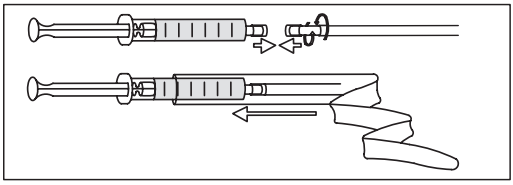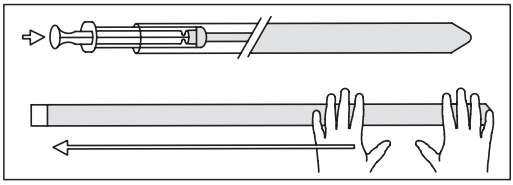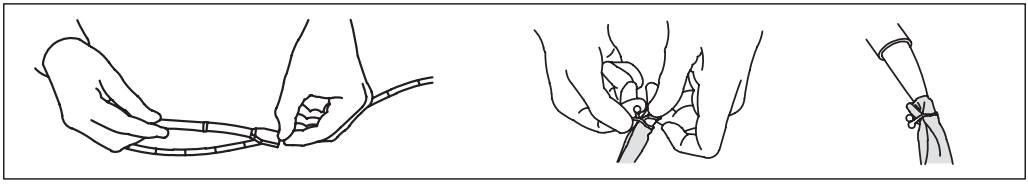Handling and Storage
- Keep away from sunlight
- After use, dispose of product and packaging in accordance with hospital, administrative and/or local government policy.

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©2025 Medline Industries, LP. All rights reserved.
Microtek Medical B.V.