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颜色：BLACK ■ PANTONE 293C

MEDLINE UROSID® Silicone SENSOR 400

Adapter Cable

en	**Adapter Cable** Instruction for Use	cs	**Propojovací kabel** Návod k použití
de	**Adapterkabel** Gebrauchsanweisung	pl	**Kabel adaptera** Instrukcja stosowania
fr	**Câble adaptateur** Notice d'utilisation	sk	**Kábel adaptéra výrobcu** Návod na použitie
nl	**Adapterkabel** Gebruiksaanwijzing	hu	**Adapterkábel** Használati utasítás
it	**Cavo adattatore** Istruzioni per l'uso	sl	**Adapterski kabel** Navodila za uporabo
es	**Cable adaptador** Instrucciones de uso	ro	**Cablu adaptor** Instrucțiuni de utilizare
pt	**Cabo Adaptador** Instruções de utilização	bg	**Адаптерен кабел** Инструкции за употреба
sv	**Adapterkabel** Instruktioner för användning	et	**Adapterikaabel** Kasutusjuhised
fi	**Sovitinkaapeli** Käyttöohje	lt	**Adapterio kabelis** Naudojimo instrukcija
da	**Adapterkabel** Brugsanvisning	hr	**Kabel adaptera** Upute za korištenje
no	**Adapterkabel** Instruksjoner for bruk	tr	**Adaptör Kablosu** Kullanım Talimatları
el	**Καλώδιο προσαρμογέα** Οδηγίες Χρήσης		

REF	Model	Connecting patient monitor brand
646060A	Phone plug	GE, Philips, Dragier, Siemens, Mindry
646070A	HP/ Philips round plug	HP, Philips
646080A	Siemens round plug	Siemens

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Instruction for Use

DESCRIPTION

The Adapter Cable is used to connect the temperature sensor to a patient monitor (YSI 400-series). The main materials of the Adapter Cable are a copper ring, PVC, HDPE.


As shown in figures 1-3 below, there are three different types of Adaptor Cable:

1. Phone plug
2. HP/Phillips round plug
3. Siemens round plug.

For each cable “Plug 1” is connected to patient monitor and “Plug 2” to the temperature sensor.


The product can be used on with patients who are using the All Silicone Foley Catheter with Temperature Sensor.

Figure 1 – Phone plug

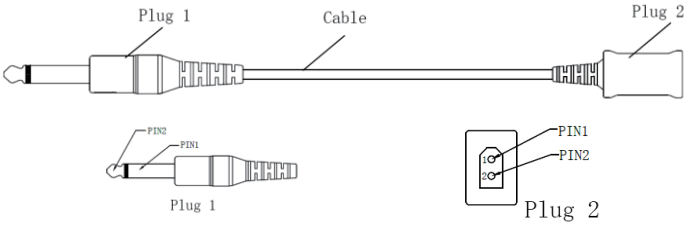
REF646060 Phone plug

Figure 2 - HP/ Philips round plug

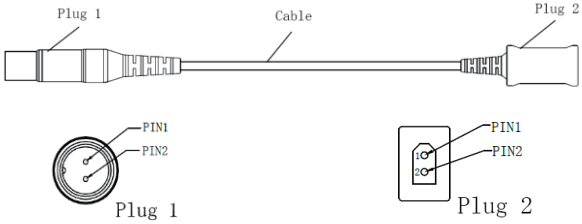REF646070 HP/ Philips round plug

Figure3: Siemens round plug

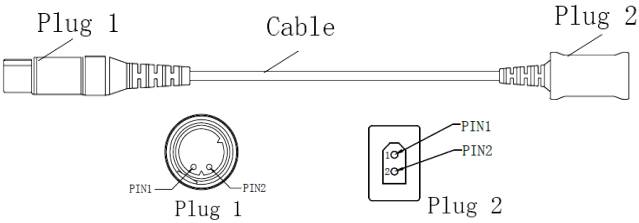REF646080 Siemens round plug

INTENDED USE AND INDICATIONS

The Adapter Cable is used to connect a patient monitor and temperature sensor when using products such as the All Silicone Foley Catheter with Temperature Sensor.

CONTRAINDICATIONS

There are no known contraindications.

POTENTIAL COMPLICATIONS

- Temperature may display incorrectly
- Risk of tissue burns from radio-frequency (RF) current

DIRECTIONS FOR USE

- Select the appropriate Adapter Cable according to TEMP interface type of patient monitor.
- Wipe the cable with an alcohol wipe, however be careful not to wipe the two plugs.

WARNINGS / PRECAUTIONS

- Do not use if the adaptor cable is damaged or irregularly shaped.
- Do not use after the expiry date.
- The product is intended for use by trained healthcare professionals and correct aseptic technique must be practiced at all times.
- Do not expose to temperatures above 49℃.
- Dispose of the product and packaging in accordance with local EBME waste disposal protocol.
- Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the appropriate regulatory authority.

ELECTRICAL CHARACTERISTICS

Input voltage: 5V DC, supplied by patient monitor.

OPERATING CONDITIONS

- Temperature: 25℃ to 45℃
- Humidity: 10% to 95% RH
- Pressure altitude: 70kPa~ 106kPa

STORAGE CONDITIONS

- Store in a dark, cool and dry place.
- Protect product from moisture and excessive heat.
- Avoid prolonged exposure to ultraviolet, sunlight and fluorescent light.
- Avoid storage that may crush the product.
- Stock rotation on first in first out basis.
- Temperature: -15~49℃
- Humidity: ≤95% RH
- Pressure altitude: 70kPa~ 106kPa

MEANING OF SYMBOLS ON PACKAGE

	Separate collection for waste electric and electronic equipment (WEEE) is required.

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Medline International B.V.,
Well Lead Medical Co., Ltd.
EC REP Shanghai International Holding Corp. GmbH(Europe)
CH REP QS Engineering AG