Skin Staple Remover

Instructions for use

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REF C-W

en Skin Staple Remover

Please read the instruction manual carefully before using it!

1. Main Structure

The metal part of the remover is made of 1Cr18Ni9Ti, and the handle of the remover is made of ABS plastic.

2. Working Mechanism

Position the jaw of the skin staple remover on the staple (between the upper and nether cutter). Close the handle of the remover and it will spread the legs of staple. When the legs have completely disengaged from the skin, the staple can be removed.

3. Capability Index

a. The metal parts of the skin staple removers are anti-corrosive.

b. The skin staple remover has been sterilized with ethylene oxide. The residual of the ethylene oxide is less than 10mg/kg.

4. Application area


The skin staple remover is used to remove the staples.

5. How-to-use

a. Position the jaws of the remover on the staple. Put the front part of the nether cutter of the remover under the staple.

b. Close the handle of the remover. The legs of the staple will be disengaged from the skin.

c. Only when the staple has disengaged completely from the skin and the legs are in a vertical position can the staple be removed.

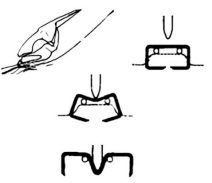6. Important Notice:

a. The remover must only be used by qualified medical person.

b. The remover should not be contaminated after it is removed from the package.

c. The remover should not be used if the packaging has been damaged.

d. Please do not use the skin staple remover after the expiration date.

7. Reserve and Transportation

The relative humidity during storage and transportation should not exceed 80% and ventilation should be good. This product should not be exposed in corrosive gases.

Not to be stored beneath other cargo. Handle with care.

8. Expiration date

5 years following sterilization.

9. After sales service

If you have any complaints about this product, please contact with suppliers or manufactors.

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Ningbo Advan Electrical Co., Ltd. Industrial Development Zone, Fuhai Town, 315332 Cixi City, Ningbo, Zhejiang Province, People's Republic of China

CE 0123

EU REP CMC Medical Devices & Drugs S.L.
C/Horacio Lengo no 18, CP 29006, Málaga, Spain

CH REP CMC Medical Devices GmbH. Rigistrasse 3,
6300 Zug, Switzerland

CH importer: Medline International Switzerland Sàrl, c/o MN & Associés SA, 1 Place de Longemalle, 1204 Genève, Switzerland

UK Representative: SUNGO Certification Company Limited3rd floor, 70 Gracechurch Street, London. EC3V OHR

Importer in UK: Medline Industries Ltd,WA4 6HL, UK

Medline International B.V.,
Catharijnesingel 47, 3511 GC
Utrecht, The Netherlands

Distributed in Australia by:
Medline International, Two Australia Pty. Ltd., Sydney, Australia

RJ25NAE 2025-10