MED-FLEX Suction System - Instructions For Use

Item code	*Item Description*
FSL200S	*MED-FLEX 2L Suction Liner With Filter*
FSL200T	*MED-FLEX 2L Suction Liner Without Filter*

EN **READ INSTRUCTIONS COMPLETELY PRIOR TO USE**


INTENDED USE

The MED-FLEX disposable suction system is used in medical treatments for the hygienic suction of various fluids and/or secretions from patients. The suction system does not come into contact with the patient's body.

Precautions :

-Allow for sufficient length of suction tubing (patient/liner) to prevent canister from tilting during use.
-For proper functioning of suction canisters, use appropriate vacuum tubing between canister and vacuum source that resist from kinking.
-Vacuum should only be introduced after correct device connections are made.
-Ensure vacuum suction remains on until all connecting tubing is disconnected and ports are securely capped.
-Always wear proper personal protective apparel when handling infectious waste.
-Reuse of single use devices creates a potential risk for the patient or user. It may cause contamination and/or impairment of
function; which may lead to injury, illness and/or death of the patient.
-If a serious incident has occurred in relation to the device, it should be reported to Pergo and the competent authority.

SINGLE SET UP

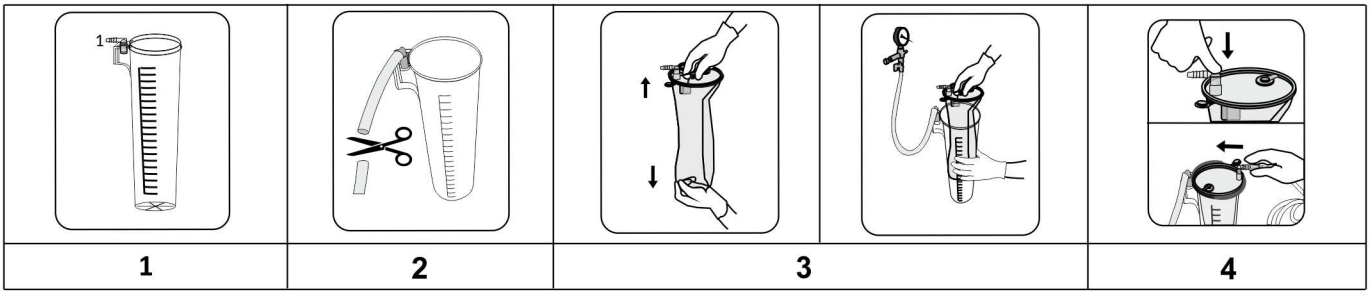Directions for set up and use


1.Make sure the reusable outer canister is in an upright position and well fixed.
2.Connect the vacuum tubing from the vacuum source to the vacuum L-connector on the canister.
3.Install the liner with filter :
 a)Remove elastic, unfold and pull the liner down.
 b)Insert the liner into the appropriately sized reusable outer canister.
4.Turn the vacuum on. Make sure the on/off valve is open (if any).
5.Secure the patient port through the L-shaped connector using your finger, and gently press the lid to seal the suction liner and outer canister assembly.
6.Ensure that the vacuum has been created and the soft flexible liner is fully deployed.
7.Make sure all unused ports are securely capped, prior to initiating suction.
8.Finally, attach the patient tubing to the orange L-Shaped connector of the patient port.
9.The system is ready to be used.

Directions for disposal

1.Once suctioning is finished, leave the vacuum on.
2.Disconnect tregulations he patient tubing and orange L-shaped connector from the patient port.
3.Securely cap the patient port using the green cap on the lid.

4.Make sure that all other ports are securely caped as well.
5.Switch vacuum off and take the liner out of the outer canister.
6.Do not remove the suction liner before the vacuum reaches zero.
7.Dispose of liners, according to your facility's protocols or other applicable for biologically hazardous waste.

TANDEM SET UP

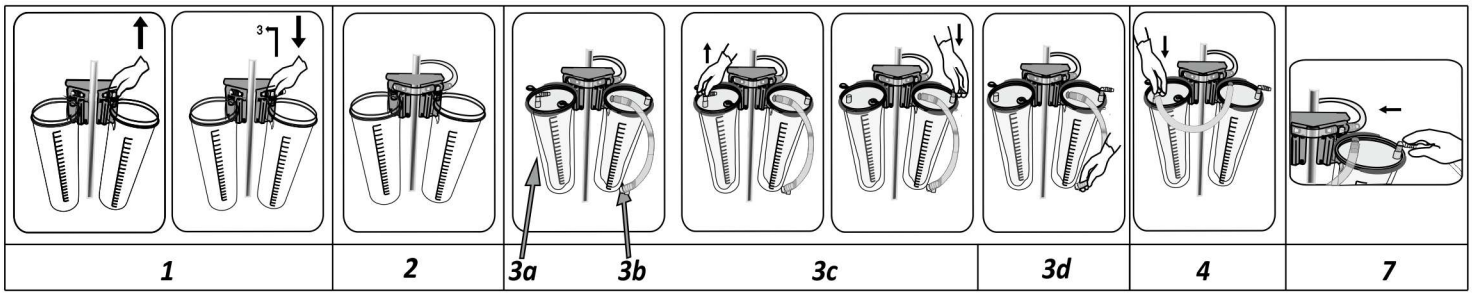
Directions for set up and use

1.Make sure the canisters are in an upright position and well installed on a canister trolley
2.Assemble vacuum circuit :
a)Replace pre-attached vacuum L-connector by T-connector except for one canister
b)Connect the canisters to each other with vacuum tubing.
3.Install liners in series:
a)Liner with filter: remove elastic, unfold and pull the liner down, then place it into the canister with vacuum L-connector.
b)Liners with tandem tubing (without filter): remove elastics, unfold and pull the liners down, then place them into the rest of canisters with vacuum T-connector.
c)Transfer orange patient L-connector from the liner with filter to the patient port of the liner without filter immediately on the right. This becomes the first liner of serial connection to which connect patient tubing.
d)Following anticlockwise pattern, connect the liners to each other with tandem tubing on the liners.
Turn the vacuum on. Make sure the on/off valve is open (if any).
4.Secure the patient port through the L-shaped connector using your finger, and gently press the lids to seal the suction liners and outer canisters assembly.
5.Ensure that the vacuum has been created and the soft flexible liner is fully deployed.
6.Make sure all unused ports are securely capped, prior to initiating suction.
7.Finally, attach the patient tubing to the L-Shaped connector of the patient port.
8.The system is ready to be used.

Directions for disposal
1.Once suctioning is finished, leave the vacuum 'ON'.
2.Disconnect the patient tubing and L-shaped connector
3.Securely cap patient port using the tandem tubing from the same liner
4.Repeat previous instruction for rest of the liners with tandem tubing.
5.Close patient port of the liner with filter using the green cap on the lid.
6.Make sure that all other ports are securely caped as well.
7.Switch vacuum off and take the liners out of the outer canisters.
8.Dispose of the liners, according to your facility's protocols or other applicable regulations for biologically hazardous waste.

For flexible liners in tandem setups, it is recommended that all used and unused disposable soft liners be discarded and subsequently replaced following each procedure as all contaminated by aerosols.

1:3 CONFIGURATION

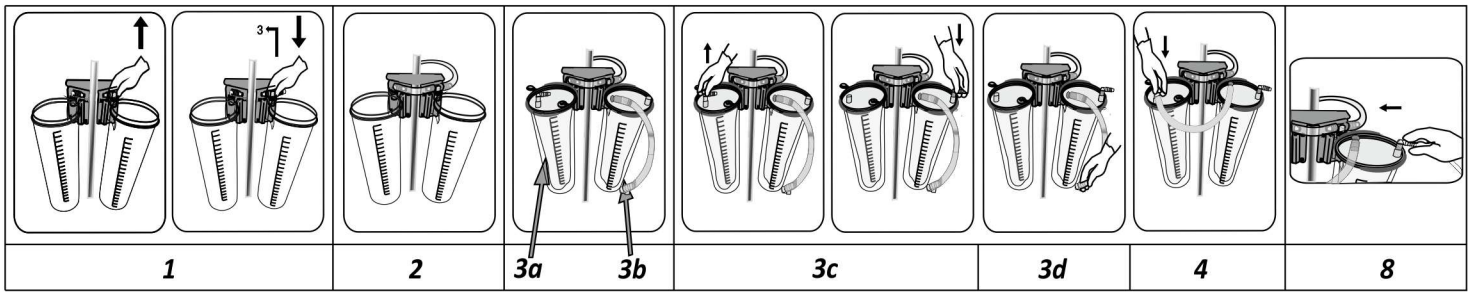1:3 configuration consist of one single and serial suction system installed in parallel. This configuration provides with greater flexibility allowing bettering adjusting the use according to the need and reducing the wastage of liner unused.

Directions for set up and use

1.Make sure the canisters are in an upright position and well installed on a canister trolley
2.Assemble vacuum circuit for 1:3 configuration :
 a)Replace pre-attached vacuum L-connector by T-connector (FSCVTC) in 2 canisters.
 b)Replace the rest of connectors by L and T vacuum connector with on/off valve (FSCVLCOF and FSCVTCOF).
 c)Connect vacuum tubing from the vacuum source to the vacuum T-connector connector with on/off valve
 d)Canisters to each other with vacuum tubing.
3.Install the single set up, select a liner with filter :
 a.Remove elastic, unfold and pull the liner down
 b.Place the liner into the canister with vacuum T-connector on/off valve (FSCVTCOF)
4.Install the tandem set up :
 a.Select a liner with filter, remove elastic, unfold and pull the liner down and place it into the canister with vacuum L-connector on/off valve (FSCVLCOF)
 b.Liners with tandem tubing (without filter): remove elastics, unfold and pull the liners down, then place them into the rest of canisters with vacuum T-connector.
 e)Transfer orange patient L-connector from the liner with filter to the patient port of **the liner without filter** immediately on the right. This becomes the first liner of serial connection to which connect patient tubing.
 f)Following anticlockwise pattern, connect the liners to each other with tandem tubing on the liners. Turn the vacuum on. Make sure the on/off valve is open (if any).
5.Secure the patient port through the L-shaped connector using your finger, and gently press the lids to seal the suction liners and outer canisters assembly.
6.Ensure that the vacuum has been created and the soft flexible liner is fully deployed.
7.Make sure all unused ports are securely capped, prior to initiating suction.
8.Finally, attach the patient tubing to the L-Shaped connector of the patient port.
9.The system is ready to be used.

Directions for disposal

1.Once suctioning is finished, leave the vacuum 'ON'.
2.Disconnect the patient tubing and L-shaped connector
3.Securely cap patient port using the tandem tubing from the same liner
4.Repeat previous instruction for rest of the liners with tandem tubing.
5.Close patient port of the liner with filter using the green cap on the lid.
6.Make sure that all other ports are securely caped as well.
7.Switch vacuum off and take the liners out of the outer canisters.
8.Do not remove the suction liner before the vacuum reaches zero.
9.Dispose of liners, according to your facility's protocols or other applicable for biologically hazardous waste.

For flexible liners in tandem setups, it is recommended that all used and unused disposable soft liners be discarded and subsequently replaced following each procedure as all contaminated by aerosols.

WARNINGS :

W1

-Do not reuse or sterilize. If the same patient undergoes long-term treatment, the flexible liner must be replaced at least every 24 hours.
-Do not subject the canister to active vacuum when not in use.
-Not intended for pleural/thoracic use.
-Inspect the device for the presence of cracks or damaged parts prior to use.
-Monitor vacuum gauge to assess or consistent vacuum pressure prior to and during use.
-The suction liner is not designed to be left under high vacuum for a long time. If using the suction liner, the vacuum source is recommended to be turned off if the procedure is not shortly started. Check that the product is intact before starting suction. Failure to do so can delay starting of the critical operation.
-Do not exceed vacuum level of 736 mmHg (1 bar).
-Do not leave the system under high vacuum for more than 1 hour while the suction tubing is closed using a stop vacuum adapter or a push-activated vacuum stop connector.
-Do not attempt to move the trolley by pulling on the patient tube and do not pull on the patient tube. The patient connector may become detached from the patient tube.
-To ensure proper seating of the lid on the canister, close the patient connector with one finger and gently press the center of the lid with one finger while the vacuum is active. Pressing on the edge of the lid may cause the lid to seat unevenly and prevent full sealing. (W1)

Cleaning of reusable canister
Autoclavable at 121°C for 20 minutes or 134°C for 4 minutes. Remove the vacuum L-connector on the canister before autoclaving.

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UK Importer : Medline Industries Ltd, WA4 6HL, UK.