en - Instructions For Use

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

en - Assembly; es - Ensamblaje; pt - Montagem; de - Versammlung; fr - Assemblée

1.

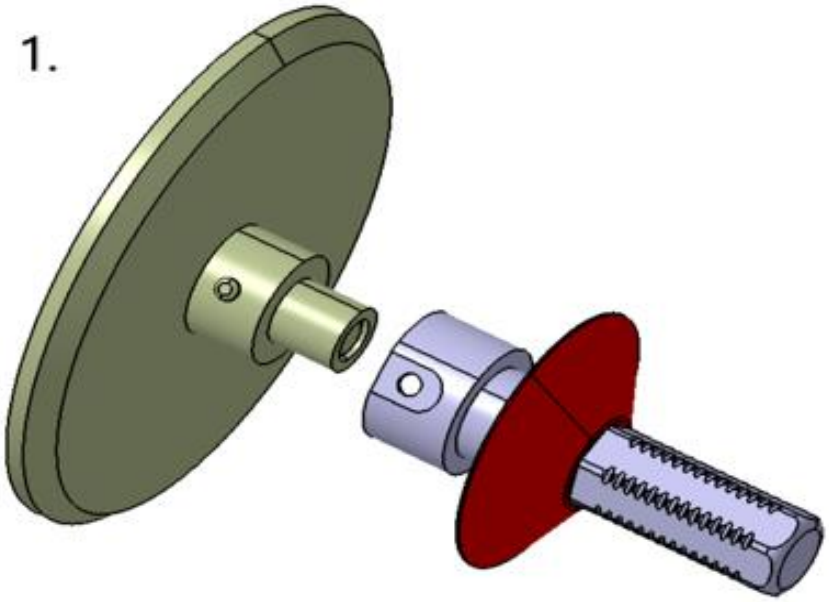2.

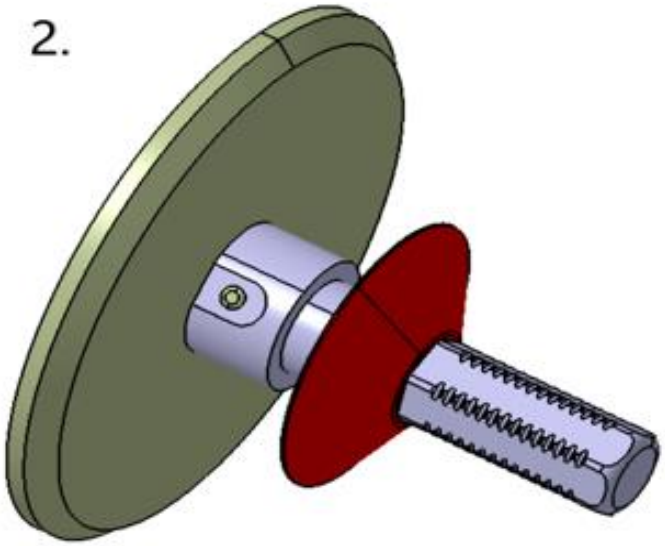en - Instructions For Use

Product Code	Product Description
3600-52S	Universal Light Handle, for use with 'Dr Mach M3, 380/400/500'
3600-54	Universal Light Handle, for use with 'Trumpf Trulight 701'
3600-54D	Universal Light Handle, for use with 'Trulight (DeRoyal)'
3600-55	Universal Light Handle, for use with 'Trumpf Flat Panel Monitor'
3600-56	Universal Light Handle, for use with 'Trumpf ILED 3 & 5'
3600-58	Universal Light Handle, for use with 'Trumpf Laser Plus'
3600-64	Universal Light Handle, for use with 'Amico Pearl & Dr Mach'
3600-74	Universal Light Handle, for use with 'Skytron DoVera'
3600-78	Universal Light Handle, for use with 'Skytron DoVera' (Camera Replacement Handle)
3600-79	Universal Light Handle, for use with 'Skytron TR7/Polaris 50'
3600-80	Universal Light Handle, for use with 'Skytron ANGINIEUX'
3600-104	Universal Light Handle, for use with 'SIEMENS 2004 Heraeus Hanaulux'
3600-107	Universal Light Handle, for use with 'CASTLE Heraeus 2007'
3600-55D	Universal Light Handle, for use with 'Trulight (DeRoyal)'
3600-58D	Universal Light Handle, for use with 'Trumpf ALC Plus Adaptor (DeRoyal)
31146619	Universal Light Handle, for use with 'SIEMENS 2004 Heraeus Hanaulux'
50001203	Universal Light Handle, for use with 'Dr Mach M3, 380/400/500'
3600-52SD	Universal Light Handle, for use with 'Dr Mach, DeRoyal'
3600-56D	Universal Light Handle, for use with 'Trumpf ILED,(DeRoyal)

This document should be thoroughly read before use.
It is designed to provide instructions on how this product is to be fitted and used.

Intended user
The 'Universal Light Handle' accessories are intended for use by clinicians working with Operating Room lighting systems

Intended purpose
The 'Universal Light Handle' accessory is intended to replace the manufacturer fitted 'Light Handle' and allow for the subsequent use of 'Universal Light Handle Cover' barrier devices to be applied to the 'Light Handle' grip.

The use of Sterile 'Light Handle Cover' barrier devices serve to prevent contamination of the sterile field by handling the lighting equipment, thus eliminating the need for detachment of fittings for reprocessing, cleaning, sterilisation, and re-fitting between uses.

Indications for use
Operating Room lighting systems, and associated accessories, are of a general purpose surgical equipment nature and have widespread application that has no condition or treatment based specificity.

How supplied
Supplied Non Sterile. **Warning:** These accessories are not intended to be sterilized prior to, or between, uses.
Supplied individually packaged in a card shelf box of dimensions: 15.24cm x 11.43cm x 11.43cm

Other devices required for use, and interoperability requirements
The 'Universal Light Handle' product codes identified in this Instructions For Use provide fitting compatibility with the lighting system as defined in their product descriptions.

Contraindications
The 'Universal Light Handles' are not intended for use with any lighting systems that are not defined in their product description.

Precautions
The user is to confirm compatibility of the selected on-market barrier cover \ sleeve prior to use, ensuring that the selected cover \ sleeve fits to the handle grip without damage (e.g. splitting) compromising its sterile barrier properties, and is suitably secure so as not to drop free into the sterile field during use.

Disposal
Warning: Due to the clinical environment of use, consideration for decontamination in accordance with relevant national guidelines shall be given. No other disposal precautions apply.

Serious incident reporting
Any serious incident that has occurred in relation to the use of these accessories should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

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