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
en - Assembly; es - Ensamblaje; pt - Montagem; de - Versammlung; fr - Assemblée

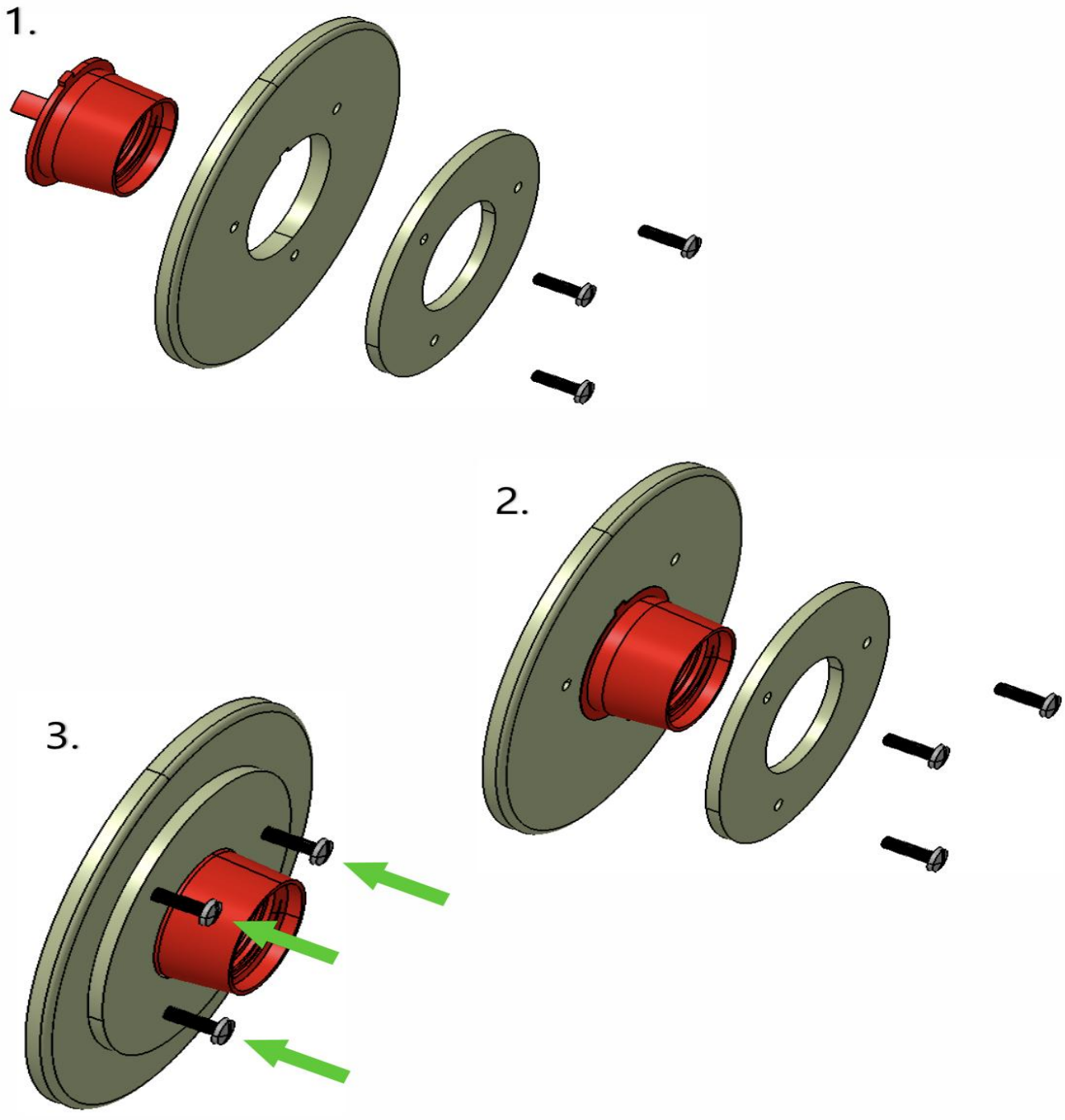en - Instructions For Use

Product Code	Product Description
3660-25	Universal Light Adaptor, for use with 'ALM PRX Adaptor'
3660-25-253	Universal Light Adaptor, for use with 'Maquet (X'Ten)'
3660-25-254	Universal Light Adaptor, for use with 'Maquet Power LED'
3660-25-297	Universal Light Adaptor, for use with 'Maquet PRX 4000/6000'
366—25-298	Universal Light Adaptor, for use with 'Maquet PRX 8000'
3660-25D	Universal Light Adaptor, for use with 'ALM PRX4000 Video'
3660-916	Universal Light Adaptor, for use with 'Volista standop adaptor'
3660-66	Universal Light Adaptor, for use with KLS Martin Disposable Handle for Monitor Adaptor ULA

⚠ **This document should be thoroughly read before use.**
It is designed to provide instructions on how this product is to be fitted and used.

Intended user
The 'Universal Light Adaptor' accessories are intended for use by clinicians working with Operating Room (OR) lighting systems

Intended purpose
The Litex 'Universal Light Adaptor' accessory is intended to permanently replace the fitted Original Equipment Manufacturer 'Light Handle' and allow for the subsequent fitting and use of rigid, disposable, 'light handle grip' devices on the lighting system.

The use of sterile disposable grip devices serve to prevent contamination of the sterile field by handling the lighting equipment, thus eliminating the need for detachment of fittings for reprocessing, cleaning, sterilisation, and re-fitting between uses.

Indications for use
Operating Room lighting systems, and associated accessories, are of a general purpose surgical equipment nature and have widespread application that has no condition or treatment based specificity.

How supplied
Supplied Non Sterile. **Warning:** These accessories are not intended to be sterilized prior to, or between, uses.
Supplied individually packaged in a card shelf box of dimensions: 7.62cm x 7.62cm x 7.62cm.

Other devices required for use, and interoperability requirements
The Litex 'Universal Light Adaptor' product codes identified in this Instructions For Use provide fitting compatibility with the lighting system as defined in their product descriptions.

Contraindications
The Litex 'Universal Light Adaptors' are not intended for use with any lighting systems that are not defined in their product description.

Precautions
The user is to confirm compatibility of the selected on-market barrier screw -on handle prior to use, ensuring that the selected screw-on handle fits securely to the adaptor so as not to drop free into the sterile field during use.

Disposal
Warning: Due to the clinical environment of use, consideration for decontamination in accordance with relevant national guidelines shall be given. No other disposal precautions apply.

Serious incident reporting
Any serious incident that has occurred in relation to the use of these accessories should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

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Litex Manufacturing Inc.
Medical Device Management Ltd