
CENTURION®
Medical Products

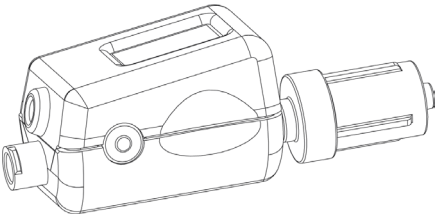Compass® Universal Hg

English....2
German....10
French....18
Dutch....26
Italian....34
Spanish....42
Portuguese....50
Swedish....58
Czech....66
Slovak....74
Slovenian....82
Norwegian....90
Polish....98
Finnish....106
Danish....114
Greek....122
Romanian....130
Bulgarian....138
Estonian....146
Hungarian....154
Latvian....162
Lithuanian....170
Croatian....178
Russian....186
Turkish....194
Korean....202
Arabic....210
Hebrew....218
Malay....226
Serbian....234

English

CENTURION®
Medical Products

Compass® Universal Hg

CAUTION: Federal (USA) law restricts this device to sale by or on the order of a physician.

Description
The Compass® Universal Hg is a self-calibrating disposable pressure transducer with integrated digital display.

Intended Use
The Compass® Universal Hg disposable pressure transducer with integrated digital display is intended for direct measurement and monitoring of physiological pressure, including during the infusion of fluids and therapeutic and diagnostic agents. Compass devices are intended for use by a licensed healthcare professional within the general population including pediatrics.

Indications for Use
The Compass® Universal Hg is indicated for use in the measurement and monitoring of compartment pressure, abdominal compartment pressure, and central venous pressure.

How Supplied
Intended for one-time use. Sterilized using ethylene oxide. Sterile if package is unopened or undamaged. Do not use the product if there is doubt as to whether the product is sterile. Avoid extended exposure to light. Upon removal from package, inspect the product to ensure no damage has occurred. If a serious incident has occurred in relation to the device, it should be reported to Medline and competent authority of the Member State in which the user and/or patient is established.

Suggested Procedure - CVC Insertion


1. Remove Compass® Universal Hg from its sterile package.
2. Activate and calibrate Compass® Universal Hg by depressing power button for approximately one second until display reads “00 mm” (Figure 1).


Warning: Turn on Compass® Universal Hg prior to attaching accessories to avoid incorrect baseline pressure.


3. Remove the proximal luer cap (Figure 2).
4. Attach introducer needle and syringe to Compass® Universal Hg (Figure 3).
5. Insert introducer needle into central vein or other vessel.
6. Verify vessel access using pressure reading. Pressure is displayed on Compass® Universal Hg LCD (Figure 4).
7. Insert guidewire (0.025˝, 0.035˝, or 0.038˝ diameter) through Compass® Universal Hg guidewire port and into vessel (Figure 5). Use of a guidewire insertion tool or “cheater” may be necessary for smaller devices.


Alternate Technique: Remove syringe and insert guidewire through syringe port on Compass® Universal Hg.


8. Hold guidewire in place and remove introducer needle and Compass® Universal Hg assembly.
9. Perform catheter insertion per catheter manufacturer’s instructions.

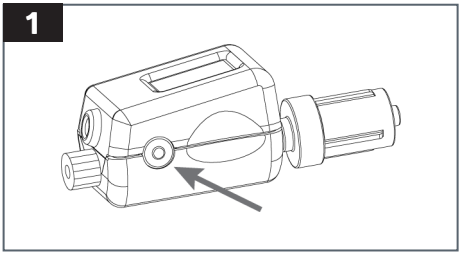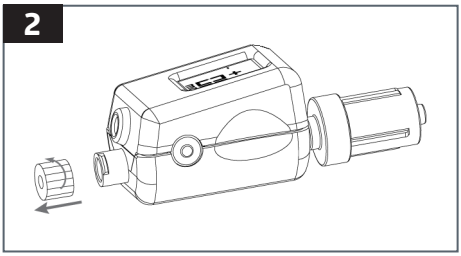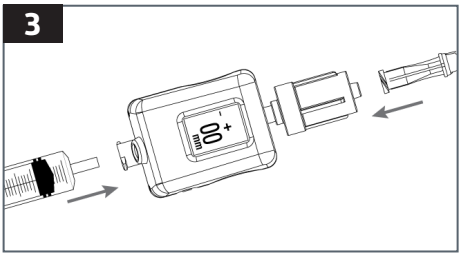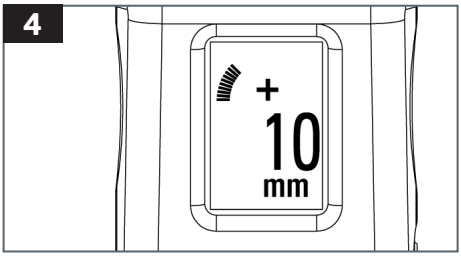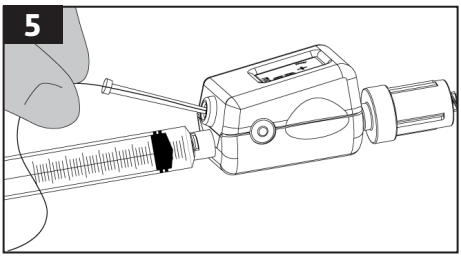General Use

1. Remove Compass® Universal Hg from its sterile package. If accessories (e.g., stopcock or syringe) are to be attached to the proximal end of the Compass® Universal Hg, remove the proximal luer cap.

2. Activate and calibrate Compass® Universal Hg by depressing power button for approximately one second until display reads “00 mm” (Figure 1).


Warning: Turn on Compass® Universal Hg prior to attaching accessories to avoid incorrect baseline pressure.


3. Connect Compass® Universal Hg to accessory monitoring equipment (e.g. needle, catheter, stopcocks, flush devices, extension sets, etc.) per hospital protocol. (Figure 2).


Note: Attach distal (male luer) side of Compass® Universal Hg to needle or catheter in contact with patient.


4. Pressure is displayed on Compass® Universal Hg LCD (Figure 3.)

If desired, insert device (maximum 0.038” diameter) through Compass® Universal Hg guidewire port (Figure 4). Use of guidewire insertion tool or “cheater” may be necessary for smaller devices.

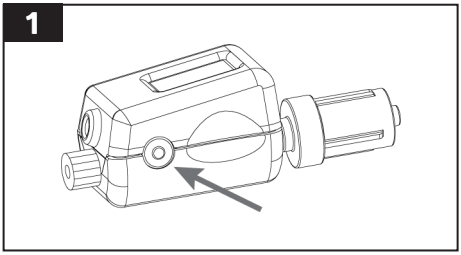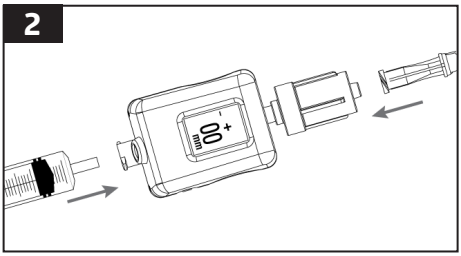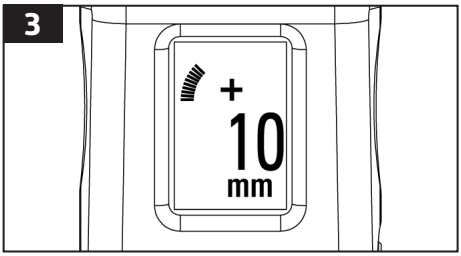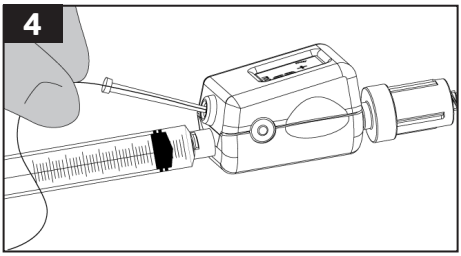Suggested Procedure - CVP Monitoring

1. Remove Compass® Universal Hg from its sterile package.
2. Activate and calibrate Compass® Universal Hg by depressing power button for approximately one second until display reads “00 mm” (Figure 1).


Warning: Turn on Compass® Universal Hg prior to attaching accessories to avoid incorrect baseline pressure.


3. Remove the proximal luer cap, then prime Compass® Universal Hg with saline and attach Compass to distal hub of central venous catheter (Figure 2).


Note: Extension tubing (sold separately) may be needed to position Compass® Universal Hg at phlebostatic axis.


4. Attach infusion tubing (sold separately) to rear of Compass® Universal Hg (Figure 3), or close rear of Compass® Universal Hg with needleless connector.
5. Position Compass® Universal Hg at phlebostatic axis, then measure central venous pressure. Pressure is displayed on Compass® Universal Hg LCD (Figure 4).
 - If Compass® Universal Hg is connected to infusion fluid, close clamp (Figure 5).
 - If Compass® Universal Hg is capped, flush to ensure patency before CVP measurement.


Note: When clamp is open during fluid infusion, the Compass® Universal Hg will report the pressure of the infusion fluid, NOT the central venous pressure (Figure 6).

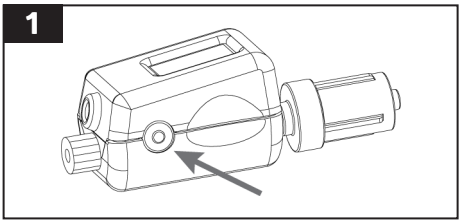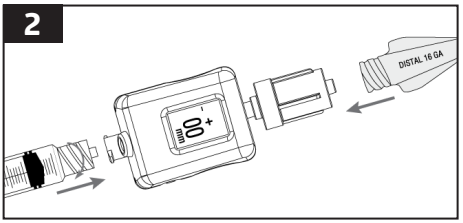
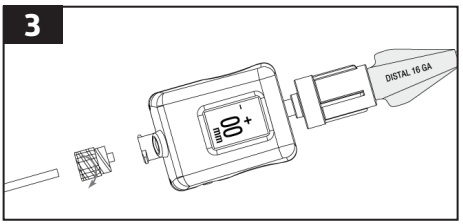
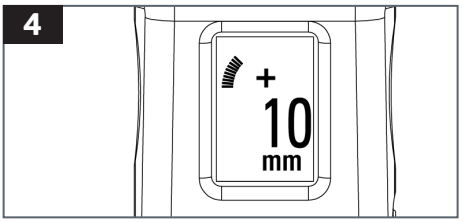
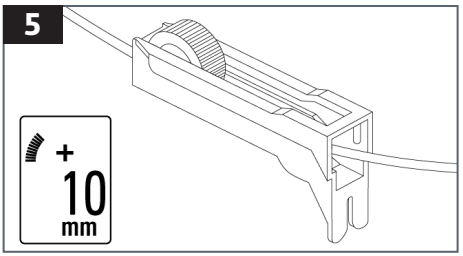
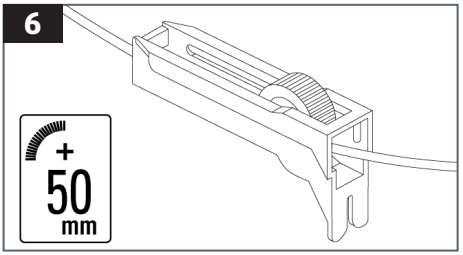
Suggested Procedure - Compartment Pressure


1. Remove Compass® Universal Hg from its sterile package.
2. Remove the proximal luer cap and connect Compass® Universal Hg to18G needle and saline filled syringe (Figure 1), then prime Compass® Universal Hg and needle with saline.
3. With the Compass® Universal Hg and needle positioned at the insertion angle (Figure 2), activate and calibrate Compass® Universal Hg by depressing power button for approximately one second until display reads “00 mm” (Figure 3).
4. Insert needle tip per hospital protocol (Figure 4) and inject 0.3 cc of saline.
5. Pressure is displayed on Compass® Universal Hg LCD (Figure 5).

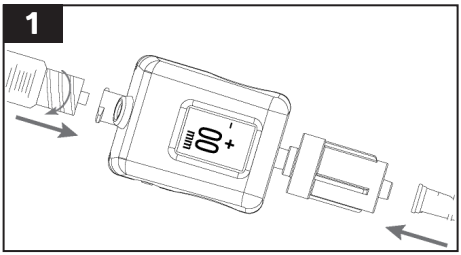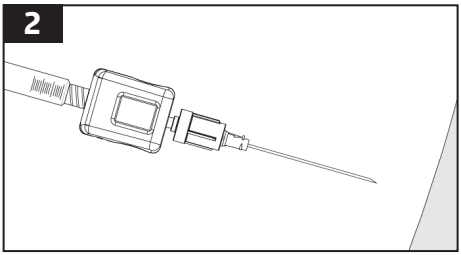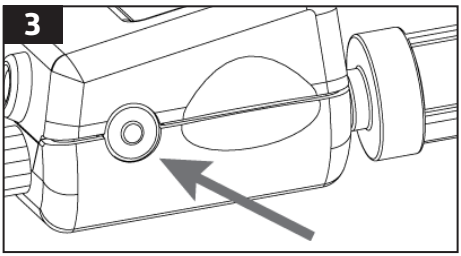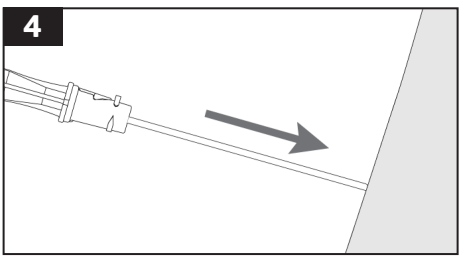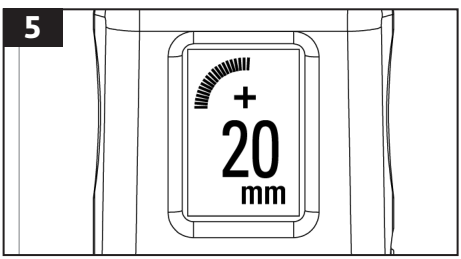Suggested Procedure - Intra-abdominal Pressure

1. Remove Compass® Universal Hg from its sterile package.
2. Activate and calibrate Compass® Universal Hg by depressing power button for approximately one second until display reads “00 mm” (Figure 1).

Warning: Turn on Compass® Universal Hg prior to attaching accessories to avoid incorrect baseline pressure.


3. Remove the proximal luer cap and prime Compass® Universal Hg and accessories prior to attaching to sampling port of Foley catheter.
4. Attach needleless valve connector to the female luer lock connection of the Compass® Universal Hg (Figure 2).
5. Connect Compass® Universal Hg to the female Luer lock connection of the Foley catheter (Figure 3A or Figure 3B).
 5a. CLAMP Foley drain bag PRIOR to instilling saline into the bladder.


NOTE: Extension tubing (sold separately) may be needed to position the Compass® Universal Hg at the mid-axillary line (level of the urinary bladder). For patients in an inclined position, place Compass® Universal Hg at the iliac crest at the mid-axillary line.


6. Instill saline per hospital protocol.


NOTE: Scrub the needleless connector with a sterile alcohol prep pad between attachment of each sterile 10 mL saline syringe.


7. Position Compass® Universal Hg at mid-axillary line, then measure intra-abdominal pressure. Pressure is displayed on Compass® Universal Hg LCD (Figure 4).

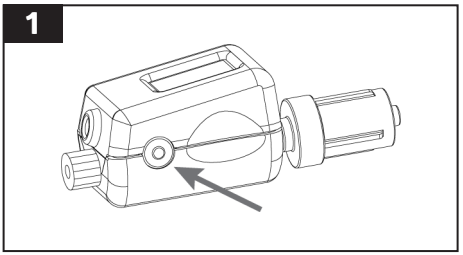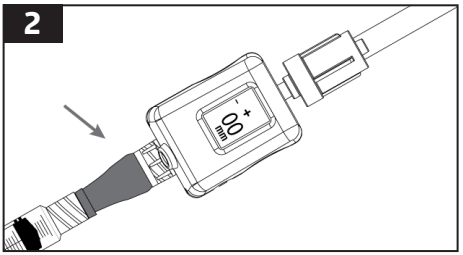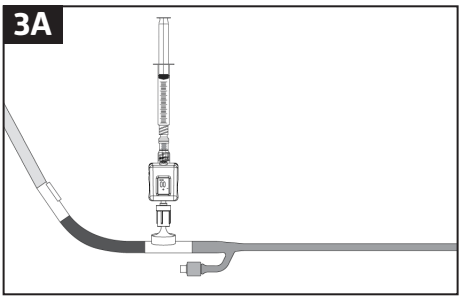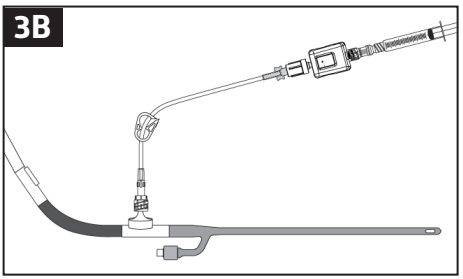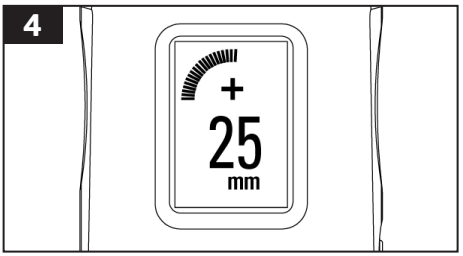Warnings

- Do not reinfuse blood or other bodily fluid into the patient after any such fluid is drawn into the Compass® Universal Hg.
- Do not reuse Compass® Universal Hg for more than one procedure in the same patient. Use a new Compass® Universal Hg device for each procedure. Reuse may lead to patient infection or cross contamination due to lack of sterility.
- To avoid a non-physiological pressure reading, do not move syringe plunger while measuring pressure.
- Pressure reading is not physiological if either front or rear of the device is open to atmospheric pressure.
- Not for use in oxygen-rich environments.


WARNING: No modification of this equipment is allowed.

The Compass® device contains a Lithium Manganese button cell battery that cannot be removed from the device housing. The whole Compass device should be disposed of as medical waste as facilitated by the healthcare facility's protocol and local member state and/or government requirements.

Device Power Off and Recalibration

Press button three times in a "PUSH - PUSH - HOLD" sequence (i.e. quickly push button twice and then hold third button push until device turns off). To recalibrate and zero the device, open the device to atmosphere and depress power button until display reads "00 mm".

Display

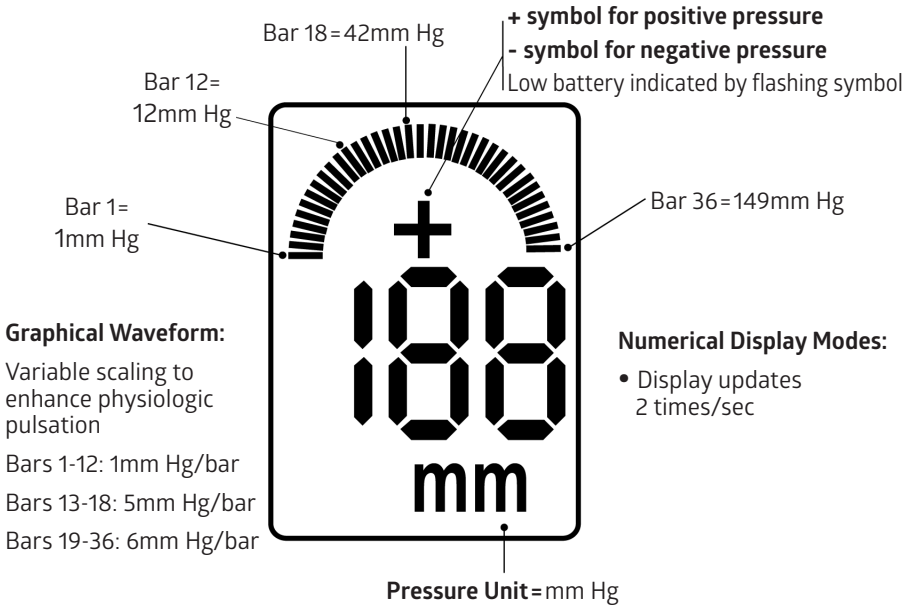Compass® Universal Hg Specifications

Characteristics	Value
Pressure Range	-199 mmHg to +999 mmHg
Accuracy	±2 mm Hg below 50 mm Hg ±3% of reading above 50 mm Hg
Zero Drift	±1 mm Hg per 4 hours
Operating Life	8 hours
Operating Temperature	15 to 40° C (59 to 104° F)
Storage Temperature	-20 to 49° C (-4 to 120° F)
Light Sensitivity	< 1 mm Hg
Lithium Battery Voltage	3 Volts

Software Error Codes

Code	Meaning
SF	Software Failure
HF	Hardware Failure
PF	Pressure Sensor Failure
bF	Battery Failure
EE	Environmental Error-Device outside operating temperature range
CE	Calibration Error-Pressure values outside operating range

©2024 Medline Centurion
Centurion Medical Products, LP
100 Centurion Way
Williamston, MI 48895 USA
www.medline.com
1-800-MEDLINE

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EC REP Medline International France SAS
Medline International B.V.
BSI Assurance UK Ltd,
Medline Industries Ltd,