尺寸：210x145mm
颜色：BLACK

MED-LAP

Asid-Bonz Extraction Bags (Operation Use)

EN	Extraction Bag（Operation Use） Instruction for Use	CS	Extrakční sáček（Použití v provozu） Návod k použití
DE	Bergebeutel（für den OP） Gebrauchsanweisung	PL	Worek do ekstrakcji（działanie i zastosowanie） Instrukcja stosowania
FR	Sac d'extraction (usage opératoire） Mode d'emploi	SK	Extrakčný Vak（Použitie pri zákrokoch） Návod na použitie
NL	Extractiezak (Operatief gebruik) Gebruiksaanwijzing	HU	Elszívó zsák（Működtetés Használat） Használati utasítás
IT	Sacca di estrazione（Uso operativo） Istruzioni per l'uso	SL	Ekstrakcijska vrečka（Uporaba pri operacijah） Navodilo za uporabo
ES	Bolsa de extracción（Uso en operaciones） Instrucciones de uso	RO	Pungă pentru extracție（Operație Utilizare） Instrucțiuni de utilizare
PT	Saco de Extração（Uso Operativo） Instruções de utilização	BG	Сак за екстракция（Употреба при операции） Инструкции за употреба
SV	Extraktionspåse（Operationsanvändning） Användningsinstruktioner	ET	Ekstraktsiooni kott (Operatsiooni-aegne kasutus) Kasutusjuhend
FI	Poistopussi（Toimenpide Käyttö） Käyttöohjeet	LT	Surinkimo maišas（naudojamas operacijos metu） Naudojimo instrukcija
DA	Udtagningspose（Operationsbrug） Brugsanvisning	HR	Vrećica za ekstrakciju（Korištenje pri operacijama） Upute za korištenje
NO	Ekstraksjonspose（Bruk under kirurgi） Instruksjoner for bruk	TR	Ekstraksiyon Torbası（Ameliyat Kullanımı） Kullanım Talimatı
EL	Τσάντα εξαγωγής (Χρήση λειτουργίας) Οδηγίες χρήσης		

REF	Product	Volume	Trocar
100200	Med-Lap Extraction Bag, Type A	200ml	10mm/12mm
100400	Med-Lap Extraction Bag, Type A	400ml	10mm/12mm
100800	Med-Lap Extraction Bag, Type A	800ml	10mm/12mm
0110310	Med-Lap Extraction Bag with Grip, Type C	300ml	10mm/12mm
0110320	Med-Lap Extraction Bag with Grip, Type C	550ml	10mm/12mm
0110380	Med-Lap Extraction Bag with Grip, Type C	800ml	10mm/12mm

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EN Extraction Bag (Operation Use)
Instruction for Use


INTENDED USE / INDICATIONS

The Extraction Bag is intended to use for removal of resected tissues of abdominal cavity during laparoscopy.

CONTRAINDICATIONS

Do NOT use if the endoscopic techniques are contraindicated.


Type A-Self-Opening Extraction Bag, for minimally invasive surgery
DESCRIPTION

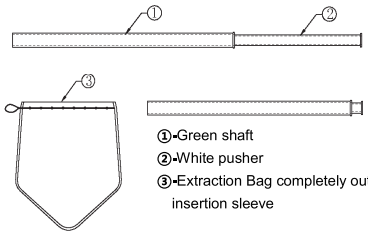DIRECTIONS FOR USE

1. Inspect packaging to ensure it is not broken or damaged.
2. Open packaging using aseptic technique.
3. Insert the green shaft ① of the system into an appropriately sized trocar and into the abdominal cavity.
4. Press the white pusher ② until the Extraction Bag is completely out of the insertion sleeve. As soon as the Extraction Bag is completely out of the insertion sleeve, the bag opens automatically ③.
5. Remove the insertion sleeve ① and the pusher ② through the trocar.
6. Use an appropriate laparoscopic instrument to insert the specimen into the Extraction Bag. (The Extraction Bag must be protected against contact with sharp objects.)
7. Grip the filled Extraction Bag at the closing wire with forceps and pull it through the trocar. This will close the opening of the bag.
8. Remove the Extraction Bag from the abdominal cavity with the trocar and forceps. If the contents of the bag are too large, enlarge the trocar incision to make it easier to extract the bag.

Type C-Self-Opening Extraction bag, detachable
DESCRIPTION

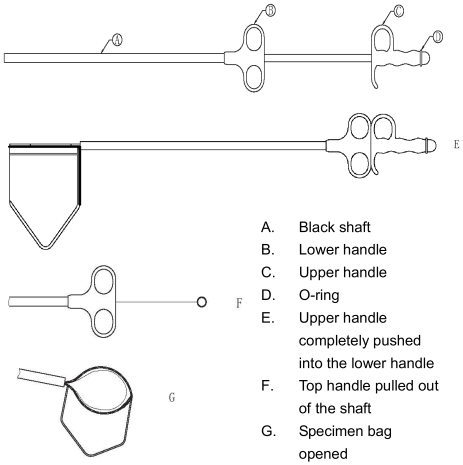DIRECTIONS FOR USE

1. Inspect packaging to ensure it is not broken or damaged.
2. Open packaging using aseptic technique.
3. Line up the arrow on the device upper handle (C) on the lower handle (B) (allowing the Extraction Bag facing upward after deployment).
4. Insert the black shaft (A) of the system (10mm) into an appropriately sized trocar and into the abdominal cavity.
5. Deploy the specimen bag by pushing the upper handle (C) with arrow towards the lower handle (B) until the flexible metal ring of the specimen bag is completely open (E) (G).
6. Use an appropriate laparoscopic instrument to insert the specimen into the extraction bag. (The Extraction Bag must be protected against contact with sharp objects).
7. Release the small red O-Ring (D) on the upper handle (C) and put it aside.
8. Pull back the upper handle (C) while holding the lower handle (B), to close and finally push out the bag.
9. Pull out the upper handle (C) from the shaft (A). (F).
10. Remove the shaft (A) and make sure that the O-ring (D) is outside of the trocar.
11. Remove the Extraction bag from the abdominal cavity once securely locked. If the contents of the bag are too large, enlarge the trocar incision to make it easier to extract the bag.

IMPORTANT

- Single Use Only
- Contents are sterile unless packaging was opened or is damaged.
- Do not use if it's damaged or irregularly shaped.
- Any serious incident that has occurred in relation to device should be reported to Medline and competent authority of the Member State in which the user and/or patient is established.
- Do NOT reuse, reprocess and/or re-sterilize. Reuse, reprocess and/or re-sterilize of single-use devices creates a potential risk for the patient or user. It may cause contamination and/or impairment of function, which may lead to injury, illness and/or death of the patient.

PRECAUTIONS

- The product mush be used exclusively by qualified medical staff.
- Thorough understanding of principles, instrumentation and techniques involving minimal invasive procedures is required to avoid possible injuries to both patient and medical staff.

DISPOSAL INSTRUCTIONS

Discard after use. Precautions should be taken when discarding this device. Disposal of the device should be done in accordance with applicable hospital or national regulations for biologically hazardous waste.

STORAGE CONDITIONS

- Store under light-proof and clean conditions. Avoid exposure to sunlight and relative humidity levels over 80%.
- Store in manner preventing crushing.

MEANING OF SYMBOLS ON PACKAGE

	Date of manufacture		**Use-by-date**	LOT	**Batch code**
	Do not re-use		**Consult instructions for use**		**Do not use if package is damaged**
STERILE EO	**Sterilized using ethylene oxide**	2 STERILIZE	**Do not resterilize**		**Single sterile barrier system with protective packaging inside**
MD	**Medical Device**		**Manufacturer**	EC REP	**Authorized representative in the European Community**
CE 0123	**CE Marked Product**		**Keep dry**		**Keep away from sunlight**
	Importer		**This way up**		**Fragile, handle with care**
REF	**Catalogue number**	49℃ -15℃	**Temperature limits: -15℃~49℃**	LATEX	**Latex Free**

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ASID BONZ GmbH, Hertzstraße 3
Well Lead Medical Co., Ltd.
EC REP **Shanghai International Holding Corp. GmbH(Europe)**
CH REP **QS Engineering AG**