Item code	Item Description	Note
FSL100E	MED-FLEX Essential 1L Suction Liner	Single Use
FSL200E	MED-FLEX Essential 2L Suction Liner	
FSL300E	MED-FLEX Essential 3L Suction Liner	
FSL100EPG	MED-FLEX Essential 1L Suction Liner (Pre-Gelled)	
FSL200EPG	MED-FLEX Essential 2L Suction Liner (Pre-Gelled)	
FSL300EPG	MED-FLEX Essential 3L Suction Liner (Pre-Gelled)	
FSL100EV	MED-FLEX Essential 1L Suction liner with valve	
FSL200EV	MED-FLEX Essential 2L Suction liner with valve	
FSL300EV	MED-FLEX Essential 3L Suction liner with valve	
FSL100EK	MED-FLEX Essential 1L Liner with 6mm x 1.8m tubing	
FSL200EK	MED-FLEX Essential 2L Liner with 6mm x 1.8m tubing	
FSL300EK	MED-FLEX Essential 3L Liner with 6mm x 1.8m tubing	
FSLTT	MED-FLEX Essential Tandem Tubing	
FSC1000	MED-FLEX 1L Reusable Outer Canister	Reusable
FSC2000	MED-FLEX 2L Reusable Outer Canister	
FSC3000	MED-FLEX 3L Reusable Outer Canister	


EN **READ INSTRUCTIONS COMPLETELY PRIOR TO USE**

Precautions :

-Allow for sufficient length of suction tubing (patient/liner) to prevent canister from tilting during use.
-For proper functioning of suction canisters, use appropriate vacuum tubing between canister and vacuum source that resist from kinking.
-Vacuum should only be introduced after correct device connections are made.
-Ensure vacuum suction remains on until all connecting tubing is disconnected and ports are securely capped.
-Always wear proper personal protective apparel when handling infectious waste.
-Reuse of single use devices creates a potential risk for the patient or user. It may cause contamination and/or impairment of function; which may lead to injury, illness and/or death of the patient.
-If a serious incident has occurred in relation to the device, it should be reported to Pergo and the competent authority.
-Using pre-gelled liners:

- As solidifier expands when absorbing fluids, measure accuracy is impacted,
- Pre-gelled liners should not be used in tandem configuration.

SINGLE SET UP

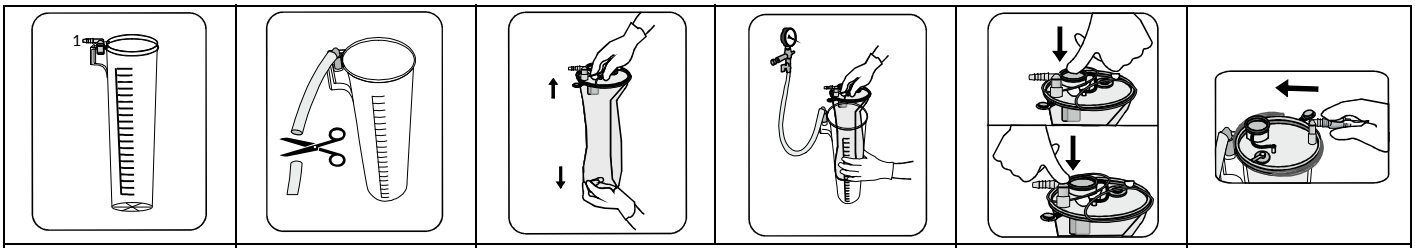1	2	3	4	5

Directions for set up and use


1.Make sure the reusable outer canister is in an upright position and well fixed.
2.Connect the vacuum tubing from the vacuum source to the vacuum L-connector on the canister.
3.Install the liners:
a)Remove elastic, unfold and pull the liner down.
b)Insert the liner into the appropriately sized reusable outer canister.
4.Make sure all unused ports are securely capped, prior to initiating suction.
5.Turn the vacuum on. Make sure the on/off valve is open (if any).
6.Secure the patient port through the L-shaped connector using your finger, and gently press the lid to seal the suction liner and outer canister assembly.
7.Ensure that the vacuum has been created and the soft flexible liner is fully deployed.

8.Attach the patient tubing to the orange L-Shaped connector of the patient port.
9.The system is ready to be used.

Directions for disposal

1.Once suctioning is finished, leave the vacuum 'ON'.
2.Securely cap the patient port using the green cap on the lid.
3.Make sure that all other ports are securely caped as well.
4.Switch vacuum off and take the liner out of the outer canister.
5.Do not remove the suction liner before the vacuum reaches zero.
6.Dispose of liners, according to your facility's protocols or other applicable for biologically hazardous waste.

TANDEM SET UP

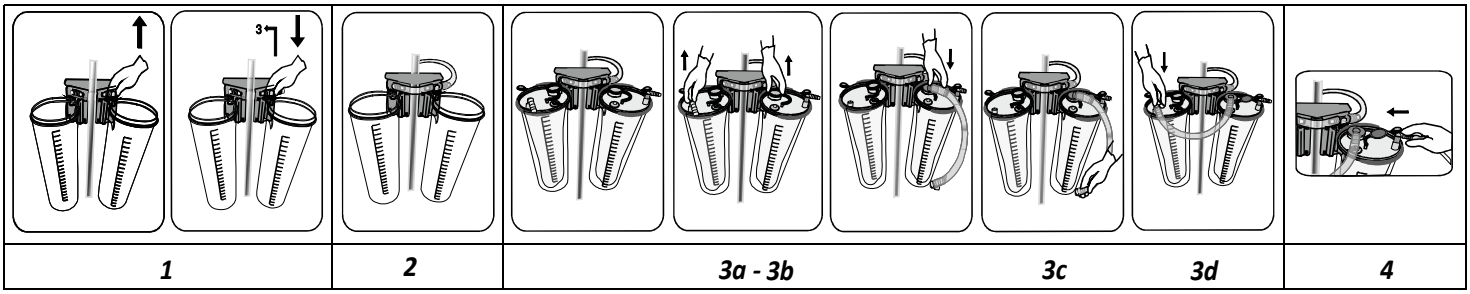
Directions for set up and use

1.Make sure the canisters are in an upright position and well installed on a canister trolley.
2.Assemble vacuum circuit :
a) Replace pre-attached vacuum L-connector by T-connector except for one canister.
b) Connect the canisters to each other with vacuum tubing.
3.Use tandem tubing (FSLTT) to link all the liners as follows:
a) Starting from the first liner, insert the orange connector of the tandem tubing (FSLTT) into the accessory port and press it firmly.
b) Then, following a clockwise pattern, remove the patient L-shaped connector from the patient port of the next liner and connect the white connector the of tandem tubing to the patient port of the second liner.
c) Repeat the tandem tube instructions (a and b) above, following a clockwise pattern until all the liners are connected.
d) No tandem tubing should be inserted to the accessory port of the last liner.
4.Make sure all unused ports are securely capped.
5.Turn the vacuum on. Make sure the on/off valve is open (if any).
6.Secure the patient port through the L-shaped connector using your finger, and gently press the lids to seal the suction liners and outer canisters assembly.
7.Ensure that the vacuum has been created and the soft flexible liners are fully deployed.
8.Attach the patient tubing to the orange L-Shaped connector of the patient port.
9.The system is ready to be used.

Directions for disposal
1.Once suctioning is finished, leave the vacuum 'ON'.
2.Disconnect the patient tubing and L-shaped connector.
3.Disconnect the tandem tubing for all liners.
4.Close patient ports of the liners using the green cap on the lid. Close all accessory ports with orange cap.
5.Make sure that all other ports are securely caped as well.
6.Switch vacuum 'OFF' and take the liners out of the canisters.
7.Do not remove the suction liner before the vacuum reaches zero.
8.Dispose of liners, according to your facility's protocols or other applicable for biologically hazardous waste.

For flexible liners in tandem setups, it is recommended that all used and unused disposable soft liners be discarded and subsequently replaced following each procedure as all contaminated by aerosols.

WARNINGS :

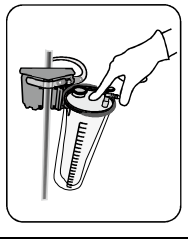W1

-Do not reuse or sterilize. If the same patient undergoes long-term treatment, the flexible liner must be replaced at least every 24 hours.
-Do not subject the canister to active vacuum when not in use.
-Not intended for pleural/thoracic use.
-Inspect the device for the presence of cracks or damaged parts prior to use.
-Monitor vacuum gauge to assess or consistent vacuum pressure prior to and during use.
-The suction liner is not designed to be left under high vacuum for a long time. If using the suction liner, the vacuum source is recommended to be turned off if the procedure is not shortly started. Check that the product is intact before starting suction. Failure to do so can delay starting of the critical operation.
-Recommended Dosage for Pre-Gelled Liners: Dosage is dependent on the composition of the fluids. A high percentage of electrolytes in the solution typically requires a higher dosage. Gallbladder fluids can partially reduce the effectiveness (fat).
-Do not exceed vacuum level of 736 mmHg (1 bar).
-Do not leave the system under high vacuum for more than 1 hour while the suction tubing is closed using a stop vacuum adapter or a push-activated vacuum stop connector.
-Do not attempt to move the trolley by pulling on the patient tube and do not pull on the patient tube. The patient connector may become detached from the patient tube.
-To ensure proper seating of the lid on the canister, close the patient connector with one finger and gently press the center of the lid with one finger while the vacuum is active. Pressing on the edge of the lid may cause the lid to seat unevenly and prevent full sealing. (W1)

Cleaning of reusable canister
Autoclavable at 121°C for 20 minutes or 134°C for 4 minutes. Remove the vacuum connectors on the canister before autoclaving.

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UK Importer : Medline Industries Ltd, WA4 6HL, UK.