
Invisishield™ Wound Protector

de **Wundrandschutzfolie**
fr **Protecteur de plaie**
nl **Wondbeschermer**
it **Telo ad anello in PE con foro centrale**
es **Protector Incisión Quirúrgica**
pt **Protetor Incisão Cirúrgica**
sk **Protektor rany**
cs **Ochrana rány**
sv **Sårskydd**
fi **Haavanreunansuojus**
da **Sårafdækning**
no **Sårbeskyttelse**
el **Προστατευτικό τραύματος**
pl **Ochraniacz na rany**
sl **Zaščita za rane**
ro **Protector plaga**
bg **Предпазител за рани**
et **Haavakaitse**
hu **Sebszél védő gyűrűs fólia**
lv **Brūču aizsargs**
lt **Žaizdos apsaugojimo priemonė**
hr **Zaštita za ranu**
tr **Yara koruyucu örtü**
ar **واقي الجروح**


REF	Ø
IDE3707	7cm
IDE3712	12cm
IDE3718	18cm
IDE3723	23cm
IDE3728	28cm

1

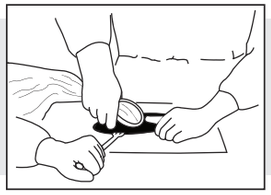en Squeeze the malleable ring together and insert into the wound
de Drücken Sie den verformbaren Ring zusammen und führen Sie ihn in die Wunde ein.
fr Presser l'anneau malléable et l'insérer dans l'incision.
nl Knijp de kneedbare ring samen en plaats deze in de wonde.
it Stringere l'anello malleabile e inserirlo nella ferita.
es Apriete el anillo maleable e introdúzcalo en la herida.
pt Apriete el anillo maleable e insértelo en la herida
sk Stlačte k sebe tvárny krúžok a vložte ho do rany.
cs Stiskněte k sobě tvárný kroužek a vložte do rány.
sv Kläm samman den formbara ringen och för in den i såret.
fi Purista muotoutuva rengas kasaan ja aseta haavaan.
da Klem den formbare ring sammen og indsæt i såret.
no Press den myke ringen sammen og før den inn i såret.
el Πιέστε τον εύπλαστο δακτύλιο και εισάγετε στην πληγή
pl Ścisnąć plastyczny pierścień i wprowadzić go do rany.
sl Upogibljiv obroč stisnite skupaj in ga vstavite v rano.
ro Strângeți complet inelul maleabil și introduceți-l în plagă.
bg Стиснете гъвкавия пръстен и го поставете в раната
et Suruge deformeeritav rõngas kokku ja sisestagesee haavasse
hu Nyomja össze a formálható gyűrűt, majd helyezze a sebbe.
lv Saspiediet kopā elastīgo gredzenu un ievietojiet to brūcē
lt Suspauskite lankstų žiedą ir įkiškite jį į žaizdą
hr Stisnite kovani prsten i umetnite ga u ranu
tr Genleşebilir bileziği sıkarak küçültün ve yaranın içine yerleştirin
ar اضغط على الحلقة المرنة وأدخلها في الجرح.

2

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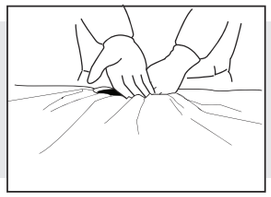en
Intended Use and Clinical Benefits:
Wound protector drapes are specifically designed for abdominal and thoracic surgical procedures, giving health care professionals the benefit of a larger surgical view by retracting the abdominal wall. They keep the wound edges moist, reducing trauma to the incision site.
The wound protector has also the advantage of delivering maximum exposure of the abdominal cavity with minimum incision size and protect against wound contamination during laparoscopic and open surgeries
Directions for Use:
These directions are general guidelines intended for use by qualified medical personnel.
Disposal Instructions:
Precautions should be taken when discarding this device and disposal of the device shall be made in accordance with applicable hospital or national regulations for biologically hazardous waste.
Contraindications and Warnings:
1. Do not use if package is opened or damaged
2. Do not expose to a flame or high intensity heat sources.
3. Reuse of single use devices creates a potential risk for the patient or user. It may cause contamination and/or impairment of function, which may lead to injury, illness and/or death of a patient.
4. If a serious incident has occurred in relation to the device, it should be reported to Medline and the competent authority of the Member State in which the user and/or patient is established.

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Medline International France SAS,
Medline International, Two Australia Pty. Ltd., Sydney, Australia
Medline Industries Ltd, WA4 6HL, UK