MEDLINE

REF S3500P REF S3513P
REF S3510P REF S3581P
REF S3540P

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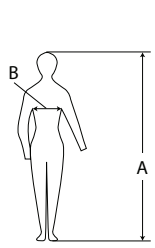
EU declaration of conformity available on: www.medline.eu/document-library

CE 2797 Medical Device class I sterile under Regulation 2017/745,
PPE category III under Regulation (EU) 2016/425,
certified by BSI Group The Netherlands B.V. 1066 EP, Amsterdam, Netherlands.

Medline International France SAS, 5 rue Charles Lindbergh, 44110 Châteaubriant, France
Distributed in Australia by: Medline International, Two Australia Pty. Ltd., Sydney, Australia

UK CA 0086 **UK approved body:** BSI Assurance UK Ltd, MK5 8PP, UK
UK representative and importer in UK: Medline Industries Ltd, WA4 6HL, UK

CH REP **CH Rep and importer:** Medline International
Switzerland SARL, c/o MN & Associés SA,
1 Place de Longemalle, 1204 Genève, Switzerland

LAB190419 - 2024.02

en Personal Protective Equipment of Complex Design Category III, in compliance with Regulation (EU) 2016/425 – Type PB [4]-B and PPE Regulation (EU) 2016/425 as amended to apply in GB. Only the gown's reinforced areas protect against chemicals and infective agents; Class I Medical Device in compliance Regulation 2017/745 and UK Medical device regulation.
Information on performance level is available in the corresponding Technical Data Sheet. Compliance with EN 13795-1 High Performance in reinforced area only.
The shelf life is given on the product label. The product is not known to cause skin reactions.
These gowns are made from materials with good breathability properties to reduce heat stresses and designed for comfortability.
Storage conditions: 5 – 30°C.
Precautions and warnings: For single use only, Do not reuse. Flammable material. Keep away from fire.
Reuse of single use devices creates a potential risk for the patient or user. It may cause contamination and/or impairment of function; which may lead to injury, illness and/or death of a patient.
If a serious incident has occurred in relation to the device, it should be reported to Medline and the competent authority of the Member State in which the user and/or patient is established.
Disposal:
The used protective gowns have to be disposed as hazardous waste. Special marked containers have to be used. Local regulations must be considered.
Gowns must remain in their primary packaging until required.
Declaration of conformity: www.medline.eu/document-library

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