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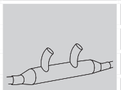
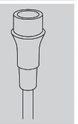
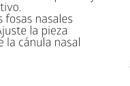
**Distributed in Australia by:** Medline International
Two Australia Pty. Ltd., Sydney, Australia

RC26SWY 2026-03

| Softech Plus | Prong[8] | | | Oxygen tube[11] | | | |
|---|---|---|---|---|---|---|---|
| Reference[1] | Maximum flow rate[2] | Size[3] | Type[9] | Length[10] | Type[9] | Colour[12] | End connector[15] |
| 1870CE | 6 LPM | Adult[4] |  | 2.1 m | Star Lumen | Green[13] |  |
| 1871CE | 4 LPM | Pediatric[5] | | 2.1 m | | | |
| 1872CE | 2 LPM | Infant[6] | | 2.1 m | | | |
| 1873CE | 2 LPM | Neonatal[7] | | 2.1 m | | | |
| 1877CE | 6 LPM | Adult[4] | | 4.2 m | | | |

| Softech | Prong[8] | | | Oxygen tube[11] | | | |
|---|---|---|---|---|---|---|---|
| Reference[1] | Maximum flow rate[2] | Size[3] | Type[9] | Length[10] | Type[9] | Colour[12] | End connector[15] |
| 1820CE | 6 LPM | Adult[4] |  | 2.1 m | Star Lumen | Transparent[14] |  |
| 1826CE | 4 LPM | Pediatric[5] | | 2.1 m | | |  |
| 1828CE | 2 LPM | Infant[6] | | 2.1 m | | |  |

1 Reference; Referenz; Référence; Artikelnummer; Articolo; Referencia; Referência; Referens; Viite; Reference; Referanse; Αναφορά; Reference; Informacje referencyjne; Referenca; Referință; Справка; Referencia; Referencia; Ürün Kodu

2 Maximum flow rate; Maximale Flussrate; Débit maximum; Maximale doorstroomsnelheid; Portata del flusso massimo; Tasa de flujo máximo; Taxa de fluxo máxima; Maximalt flöde; Suurin virtausnopeus; Maksimal strømningshastighed; Maksimal strømningshastighet; Μέγιστος ρυθμός ροής; Maximální průtok; Maksymalne natężenie przepływu; Največja stopnja pretoka; Debit maxim; Максимален дебит; Maximális áramlási sebesség; Maximálny prietok; Maksimum akış hızı

3 Size; Größe; Taille; Afmeting; Taglia; Tamaño; Tamanho; Storlek; Koko; Størrelse; Μέγεθος; Velikost; Wymiary; Velikost; Dimensiune; Размер; Méret; Velkosť; Boyut

4 Adult; Erwachsene; Adulte; Volwassen; Adulto; Adulto; Adulto; Vuxen; Aikuinen; Voksen; Voksen; Ενηλίκων; Dospělí; Dla dorosłych; Odrasli; Adult; Възрастен; Felnőtt; Dospelý; Yetişkin

5 Pediatric; Pädiatrisch; Pediatrique; Pediatrisch; Pediatrico; Pediatrica; Pediátrico; Pediatrisk; Lapsi; Pædiatrisk; Barn; Παιδιατρική; Pediatrický; Dla dzieci; Otroci; Pediatric; Педиатричен; Gyermek; Pediatrický; Pediatrik

6 Infant; Säugling; Nourrisson; Zuigeling; Lattante; Infantil; Infantil; Spædbarn; Pikkulapsi; Spædbarn; Spedbarn; Βρεφική; Kojenec; Dla niemowląt; Dojenčki; Bebélus; Дете; Csecsemő; Kojenec; Bebek

7 Neonatal; Neugeborene; Nouveau-né; Neonataal; Neonatale; Neonatal; Neonatal; Nyfödd; Vastasyntynyt; Nyfødt; Neonatal; Νεογνική; Novorozenec; Dla noworodków; Novorojenčki; Neonatal; Новородено; Újszülött; Novorodenec; Yenidoğan

8 Prong; Nasenstück; Embout nasal; Pinnen; Rebbi; Punta; Ponta; Nässtycke; Piikki; Næseklemme; Tapp; Ρινική Διάταξη; Hrot; Wąsy tlenowe; Držalo; Mărțară nazală; Зъбец; Orrszonda; Hrot; Burun aparatı

9 Type; Typ; Type; Type; Tipo; Tipo; Tipo; Typ; Τyppi; Type; Type; Τύπος; Typ; Typ; Vrsta; Tip; Тип; Tipus; Typ; Tip

10 Length; Länge; Longueur; Lengte; Lunghezza; Longitud; Comprimento; Längd; Pituus; Længde; Lengde; Μήκος; Délka; Długość; Dolžina; Lungime; Дължина; Hossz; Dĺžka; Uzunluk

11 Oxygen tube; Sauerstoffschlauch; Tubulure à oxygène; Zuurstofslang; Tubo d'ossigeno; Tubo de oxigeno; Tubo de oxigénio; Syrgasslang; Happiletku; Iltslange; Oksygenslange; Σωλήνας οξυγόνου; Kyslíková trubice; Rurka tlenowa; Cev za kisik; Tub de oxigen; Кислородна тръба; Oxigén cső; Kyslíková trubica; Oksijen tüpü

12 Colour; Farbe; Couleur; Kleur; Colore; Color; Cor; Färg; Väri; Farve; Farge; Χρώμα; Barva; Kolor; Barva; Culoare; Цвят; Szín; Farba; Renk

13 Green; Grün; Vert; Groen; Verde; Verde; Verde; Grön; Vihreä; Grøn; Grønn; Πράσινο; Zelená; Zielony; Zelena; Verde; Зелен; Zöld; Zelená; Yeşil

14 Transparent; Transparent; Transparent; Transparant; Trasparente; Transparente; Transparente; Tranşparent; Läpinäkyvä; Transparent; Gjennomsiktig; Διάφανο; Průhledný; Przezroczysty; Prozoren; Transparent; Прозрачен; Átlátszó; Priehľadný; Şeffaf

15 End connector; Endanschluss, Connexion; Aansluiting; Connettore; Conector final; Conexão; Anděkontakt; Päätyliitin; Endetilslutning; Endekobling; Χώνεψους; Koncový konektor; Złącze końcowe; Končni priključek; Conector de capăt; Краен конектор; Végcsatlakozó; Koncový konektor; Uç konnektör

## en

### Intended purpose:

Nasal Cannulas are used to deliver a low flow of oxygen to the patient through the nostrils. These devices are intended for use by medical personnel, lay persons, or patient that have been prescribed oxygen. Intended to be used for a maximum of 30 days.

### Directions for Use:

Prior to use, ensure compatibility with gas outlet. Standard connectors are compatible with 5-7mm male fittings and universal connectors are compatible with male DISS 1240 fittings or 5-7mm male fittings.
Attach nasal cannula by pushing end connector directly onto compatible gas outlet.
Set oxygen to prescribed flow rate, as directed by clinician.
Check for secure connections by confirming flow of gas through device.
Insert nasal prongs into nostrils and position lariat tubing over the ears. Adjust bolo slide under chin until nasal cannula is secure.

### Warnings:

- Non-sterile. Do not sterilize prior to use.
- Do not use in the presence of an ignition source.
- Do not use if device is soiled or damaged.
- Prior to use, ensure compatibility with gas outlet.
- Prior to use, ensure all connections are secure. Do not over tighten or force connections. Do not alter device. Improper connections may result in disconnection or leakage during use, which may result in performance degradation and increase the risk of respiratory distress.
- Single Use Only: Do not reuse, reprocess, or re-sterilize. Reuse of device creates a potential risk of serious injury and/or infection which may lead to death. Reprocessing of medical devices intended for single use only may result in degraded performance or a loss of functionality.
- If a serious incident has occurred in relation to the device, it should be reported to Medline and the competent authority of the Member State in which the user and/or patient is established.

### Disposal instructions:

Precautions should be taken when discarding this device and disposal of the device shall be made in accordance with applicable hospital or national regulations for biologically hazardous waste.

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