HUDSON RCI®

REF 1852, 1853, 1870, 1871, 1872, 1873, 1874, 1876, 1877, 1878


NON STERILE · Rx Only · single use only · caution

Not Made with DEHP

Not made with natural rubber latex.

SOFTECH® Plus

Nasal Cannula w/ Green Star Lumen Tubing

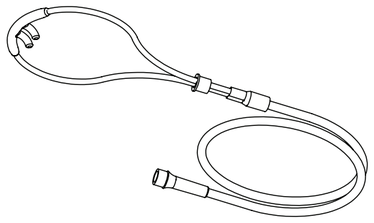Manufactured for:

Teleflex MEDICAL

3015 Carrington Mill Blvd
Morrisville, NC 27560 USA
USA: (866) 246 6990
International: +1 919 544 8000

© 2011, 2012, 2016, 2017 Teleflex Incorporated
L07099R01

CE 0120 EC REP Teleflex Medical, IDA Business and Technology Park, Dublin Rd., Athlone, Co. Westmeath, Ireland

SOFTECH, STAR LUMEN and HUDSON RCI are registered trademarks of Teleflex Incorporated in the United States of America and in selected other countries

English (EN)

SOFTECH® Plus

Nasal Cannula w/ Green Star Lumen Tubing

Directions for Use:
1. Attach the oxygen supply tubing to the gas source and set oxygen to prescribed flow.
2. Check for gas flow throughout the device.
3. Insert the nasal tips into the nostrils with the two plastic tubes over the ears and under the chin.
4. Gently adjust plastic slide until the cannula is secure.
5. If necessary, secure the cannula to the patient's face with tape.

Cautions:
- Do not use if package is open or damaged.
- Do not store at extreme temperature and humidity. Store in a cool and dry place.
- Ensure that all connections are secure and air circulates freely through the tubing.

Warnings:
- Re-processing of products intended for single use only may result in degraded performance or a loss of functionality. Re-use of single use only products may result in exposure to viral, bacterial, fungal, or prionic pathogens. Validated cleaning and sterilization methods and instructions for re-processing to original specifications are not available for these products. This product is not designed to be cleaned, disinfected, or re-sterilized.
- This nasal cannula is intended to supply supplemental oxygen only. Take care to avoid continuous positive airway pressure (CPAP) while using this device.

SIZE	REF	SIZE	REF
Adult	1852, 1870, 1874, 1877	Pediatric	1853, 1871, 1876, 1878
Neonate	1873	Infant	1872

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