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


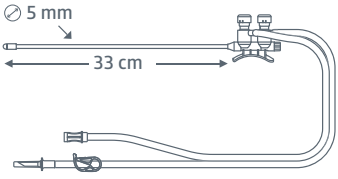
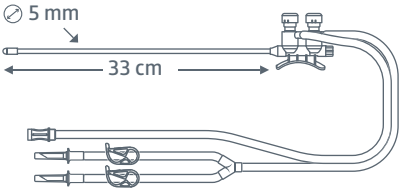
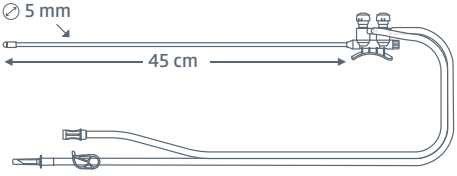
REF ORC45


en 1- Trumpet Valve

Indication for use:
The trumpet valve for suction irrigation (single spike or double spike) is used to deliver irrigation fluids to surgical sites during laparoscopic procedures and to remove fluid waste, tissue debris and smoke from the surgical site.

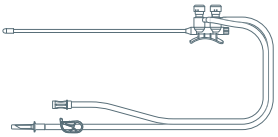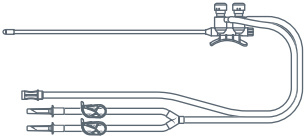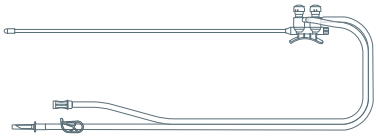OR1200S : Trumpet Valve with 33 cm cannula for suction and irrigation - Single Spike	OR1201D : Trumpet Valve with 33 cm cannula for suction and irrigation – Double Spikes	OR1245S : Trumpet valve with 45 cm cannula for suction and irrigation
		

Description:
A – The single use trumpet valve for suction & irrigation, comprises of a double piston hand piece, a cannula and tubings :

1. The cannula is screwed on the front of the double piston hand piece. Cannula distal tip is protected with a cap.
2. A cap is fitted (screwed) at the rear of the double piston hand piece
3. The blue valve (back of the hand piece) is for irrigation. Connected next to the blue piston, the irrigation line comes with a clamp and a spike to connect to a saline bag. For dual spikes trumpet valves, the irrigation line comes with a "Y" connector with 2 spikes to connect to two saline bags.
4. The red valve (front of the hand piece) is for suction
5. The extremity of the tube comes with a female connector to connect to a suction liner/canister

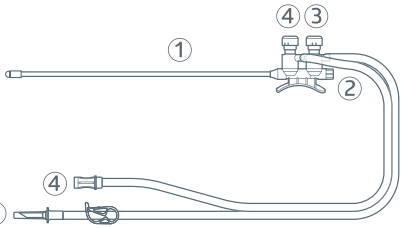B – Using the fingers :

- By pressing the blue valve (3),
- it will control the flow of irrigation fluid from the saline bag through the irrigation tube, the hand piece, the cannula and be delivered on to the surgical site.

- The red valve (4) is for suction. The extremity of the tubes comes with a female connector to connect to a suction liner/canister.

Note: If it is necessary to separate the two tubes, simply pull in opposite directions.

Instruction for use:

1. Check that the package is not compromised or the content is not expired. Do not use the device if the sterility has been compromised.
2. Open the pouch and transfer onto sterile field
3. Inspect the device for overall condition and physical integrity
4. Remove the spike(s) cap(s), connect the spike shape connector(s) to the irrigation bag(s)
5. Connect the suction line to the suction liner/canister and adjust vacuum pressure
6. Before beginning the procedure, prime lines, test suction irrigator, do not use if defects are noted.
7. Press blue valve for irrigation, press red valve for suction.

Note: When the suction valve is released and suction does not discontinue, depress the irrigation valve in order to stop suction. After applying the alternate valve, rapidly press and release the suction valve. Repeat process as many times as necessary in order to fix suction valve. If these efforts fail, remove the cannula from the surgical site, press and hold the suction valve, and wait ten seconds. If problem continues, stop the use of this irrigator and replace with a new device.

Precaution:
When not using suction, ensure that the trumpet does not remind depressed so as to minimize loss of pneumoperitoneum.

2- Suction Irrigation Cannula

ORC45
Trumpet valve with 45 cm cannula for suction and irrigation

Indication for use:
Cannula compatible with Medline trumpet valves : OR1200S, OR1245S and OR1201D.

Description:
Single-use and sterile suction-irrigation cannula to be connected to Medline suction-irrigation for enlarging cannula length on OR1200S/OR1201D or replacing cannula on OR1245S.

Instructions for use :
In order enlarge or replace existing trumpet valve cannula, remove old cannula and screw the new one:

Important:
This device was designed tested and manufactured for single use only. Any attempt to reuse or reprocess may lead to failure and potential patient injury. Reprocessing and/or re-sterilization of this device may give rise to patient infection through contamination. Do not reuse, reprocess or re-sterilize this device.

1. Refer to your suction/irrigation systems operator manual lor proper setup.
2. Assure that all suction/irrigation manufacturer's safety precautions are followed.
3. This device is intended to be used by healthcare professionals. Thorough understanding and training of the operating principles and the hazards involved in utilizing an laparoscopic approach is necessary to avoid possible injury to the user and or patient.
4. Device is provided sterile. Do not use the device if the sterility has been compromised or the packaging is damaged. Do not reuse and or re-sterilize. Discard after use.
5. Carefully inspect instrument for overall integrity before use.
6. Reuse of single use devices creates a potential risk for the patient or user. It may cause contamination and/or impairment of function; which may lead to injury, illness and/or death of a patient.
7. If a serious incident has occurred in relation to the device, it should be reported to Medline and the competent authority of the Member State in which the user and/or patient is established.

Storage & handling:

1. Device must be stored in a clean, cool and dry place.
2. It must be handled with care to avoid damage to the packaging and its content during transportation and storage.
3. Device must not be in contact with acids, organic solvents or other corrosive liquids.
4. Avoid exposure of device to sunlight.

Disposal Instructions:
Precautions should be taken when discarding this device and disposal of the device shall be made in accordance with applicable hospital or national regulations for biologically hazardous waste.

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GCMEDICA ENTERPRISE LTD.(WUXI),
EU REP Kingsmead Service B.V. Zonnehof 36, 2632 BE,
CH REP CMC Medical Devices GmbH
Medline International B.V.,