0.25 in

0.25 in

MEDLINE

HUDSON RCI®

Tracheostomy Mask

Instructions For Use

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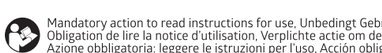
Manufactured in Mexico by:
Medline Industries, LP, Three Lakes Drive, Northfield, IL 60093 USA.

Distributed in Australia by: Medline International Two Australia Pty. Ltd., Sydney, Australia.

CE 2797 EC REP Medline International France SAS, 5 rue Charles Lindbergh, 44110 Châteaubriant, France

CH REP **CH Rep and importer:** Medline International Switzerland SARL, c/o MN & Associés SA, 1 Place de Longemalle, 1204 Genève, Switzerland

UK CA 0086 **UK representative and importer in UK:** Medline Industries Ltd, WA4 6HL,UK

Medline International B.V., Nieuwe Stationsstraat 10, Arnhem, 6811 KS, Netherlands

REF 1075CE REF 1076CE

RB25145 2025-02
HNL02175

en Tracheostomy Mask

Intended use:
Aerosol Masks are intended to deliver aerosol therapy to the patient for inhalation of nebulized medications or for humidification of the respiratory tract.
These devices are intended for use by medical personnel, lay persons, or patients that have been prescribed aerosol therapy.
1075CE: Intended for use on adults only.
1076CE: Intended for use on pediatrics only.
Intended to be used for a maximum of 30 days.

Instructions for use:
1. Connect aerosol tubing or a nebulizer to the mask inlet.
2. Check for secure connections by confirming flow of aerosol through device.
3. Place the elastic strap behind the neck.
4. Gently pull the ends of the strap until the mask is secure.
5. When using suction, loosen the mask and drop the mask out of the work area.

Note: The mask inlet swivels 360° to position the tubing for supine or upright patients.

Disposal Instructions:
Precautions should be taken when discarding this device and disposal of the device shall be made in accordance with applicable hospital or national regulations for biologically hazardous waste.

Warnings:
Single use: Do not reuse, reprocess or re-sterilize. Reuse of device creates a potential risk of serious injury and/or infection which may lead to death. Reprocessing of medical devices intended for single use only may result in degraded performance or a loss of functionality.
Non-Sterile. Do not sterilize prior to use.
Do not use if device is soiled or damaged.
If a serious incident has occurred in relation to the device, it should be reported to Medline and the competent authority of the Member State in which the user and/or patient is established.

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