TRIDENT
Sales Corporation

Your Source for Quality Medical Disposable Products

en — Non-Conductive, Single Use Coaxial Rebreathing Circuit


While this product has been carefully designed and manufactured, it is the user's responsibility to determine the suitability of this product in any intended clinical application.

CAUTION: While this product has been A.Q.L. tested by the manufacturer, qualified personnel must functionally check the product prior to use. Patient should be monitored constantly while product is in use.

WARNING: Do not soak, rinse, gas or steam sterilize these products. Attempts to clean or sterilize these products may result in harmful residual gases or dangerously increase breathing resistance.

Instructions for Use

1. Install breathing bag on anaesthesia machine connection.
2. Install circuit, noting the inspiratory and expiratory connections.
3. Prior to use, examine circuit for damage that may have occurred during shipment. Pressure test circuits according to hospital or machinery pre-check procedures.

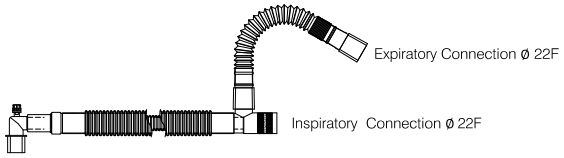Additional Cautionary Notes pertaining to Phthalates:
Products marked with PHT/DEHP symbol contains traces of phthalates. Phthalates are widely used as plasticizers in medical products and have undergone extensive testing. However, when using the device, it should be taken into consideration that possible health effects when used for long durations in pregnant women, nursing mothers and children cannot be entirely discounted.

SYMBOL LEGEND (Note: not all symbols may be included in the labeling of this product)

Manufacturer	Catalogue number	Batch (lot) number	Use by	Do not reuse	Caution				
Consult instruction for use	Contains phtalates	Contains DEHP	Phtalates Free	DEHP Free	Latex free				

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