PHT

LATEX

DEHP

en This product contains the phthalate DEHP (applicable for tubing only). Phthalates are widely used as plasticizers in medical products and have undergone extensive testing. However, when using the device, it should be taken into consideration that possible health effects when used for long durations in pregnant women, nursing mothers and children cannot be entirely discounted. Not Made With Natural Rubber Latex.

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

Medline is a registered trademark of Medline Industries, Inc.
Manufactured for | fabriqué pour | fabricado para:
Medline Industries, Inc., Three Lakes Drive, Northfield, IL 60093 USA.
Distributed in Australia by: Medline International Two Australia Pty. Ltd., Sydney, Australia
Made in China - **www.medline.com** RL20SJP 2020-12


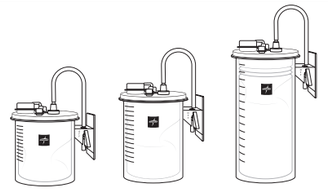
TANDEM
REF DYNDTND

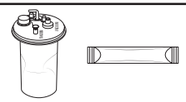		REF
	1000 cc/mL 1500 cc/mL 3000 cc/mL	DYNDSCL1000 DYNDSCL1500 DYNDSCL3000 OR1930
	1000 cc/mL 1500 cc/mL 3000 cc/mL	OR1910PG OR1920PG OR1930PG
 PHT DEHP **6 mm (1/4 in) x180 cm (6 ft)**	1000 cc/mL 1500 cc/mL 3000 cc/mL	OR916K OR926K OR936K
 PHT DEHP **6 mm (1/4 in) x300 cm (10 ft)**	1000 cc/mL 1500 cc/mL 3000 cc/mL	- OR929K OR939K
 PHT DEHP **6 mm (1/4 in) x180 cm (6 ft)**	1000 cc/mL 1500 cc/mL 3000 cc/mL	OR53916 OR53926 -
 PHT DEHP **6 mm (1/4 in) x300 cm (10 ft)**	1000 cc/mL 1500 cc/mL 3000 cc/mL	- OR53929 -
 PHT DEHP **6 mm (1/4 in) x180 cm (6 ft)**	1000 cc/mL 1500 cc/mL 3000 cc/mL	OR54916 OR54926 -
	1000 cc/mL 1500 cc/mL 3000 cc/mL	DYNDCLO1000 DYNDCLO1500 DYNDCLO3000
	1000 cc/mL 1500 cc/mL 3000 cc/mL	DYNDCLO10T DYNDCLO15T DYNDCLO30T


SETUP AND USE I SINGLE

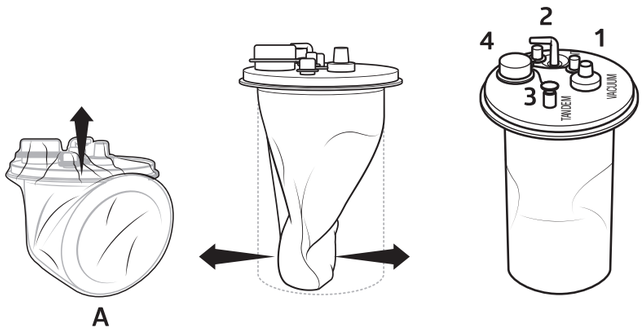A


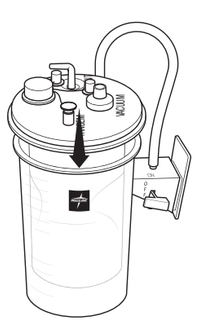
B, C

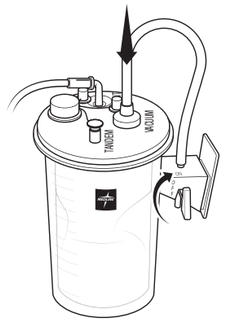D, E, F, G, H, I

SETUP AND USE I TANDEM

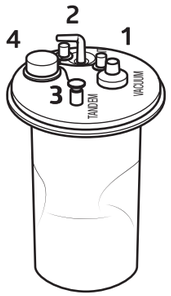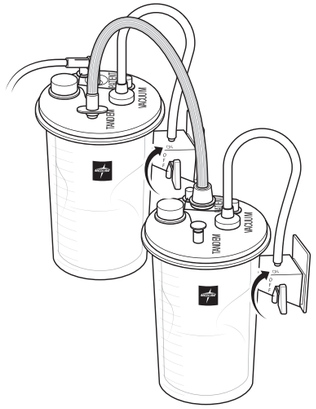en

SOFT LINERS

READ INSTRUCTIONS COMPLETELY PRIOR TO USE.

Precautions:

- Avoid tilting, jarring, or sudden movements of the canister which may activate the valve and can cause premature shut-off.
- Allow for sufficient length of suction tubing (patient/liner) to prevent canister from tilting during use.
- For proper functioning of suction canisters, use appropriate vacuum tubing between canister and vacuum source that resist from kinking.
- Vacuum should only be introduced after correct device connections are made.
- Ensure vacuum suction remains on until all connecting tubing is disconnected and ports are securely capped.
- Always wear proper personal protective apparel when handling infectious waste.
- Reuse of single use devices creates a potential risk for the patient or user. It may cause contamination and/or impairment of function; which may lead to injury, illness and/or death of a patient.
- If a serious incident has occurred in relation to the device, it should be reported to Medline and the competent authority of the Member State in which the user and/or patient is established.

SETUP AND USE for Soft Liner : SINGLE

Port identification key

(1) Vacuum Port (2) Patient Port (3) Tandem Port (4) Accessory Port

Directions for set up and use

A. Remove liner from its packaging and then pull the liner down in order to extend and straighten out.
B. Insert proper sized soft liner into the coinciding sized hard outer reusable canister. Connect tubing from the suction regulator (or other vacuum source) to the port underneath the canister bracket.
C. Press down on the liner lid in order to make sure that you have securely sealed the liner and outer canister.
D. Attach red tubing from vacuum source to the Vacuum Port (1) on the soft liner lid.
E. Securely cap off Tandem Port (3) and Accessory Port (4).
F. Attach patient tubing to the Patient Port (2) on the soft liner lid.
G. Confirm all unused ports are securely capped off prior to initiating suction.
H. Ensure that the vacuum has been created and the soft liner is fully inflated.
I. The system is ready for use

Using the soft liners

- When the suction is on, the liquid flows inside the soft liner. The soft liner can be used throughout the procedure or until it is full. When the soft liner is full, the shut-off valve cuts the suction and prevents the flow of liquid goes to the suction source.
- Shut-off valve is angle sensitive and if canister is bumped or moved suddenly away from a vertical position there is potential for premature shut-off. If suction flow is stopped prematurely, disconnect the vacuum tube from Vacuum Port (1) and allow the shut-off valve to move back into the open position. Then reconnect the vacuum tube to Vacuum Port (1).

Using the soft liners with solidifier

If adding solidifier, follow manufacturer's instructions. Solidifiers should be added before shutting off the vacuum source to prevent risk of fluid leaking from the accessory port when the solidifier is added.
As solidifier expands when absorbing fluids, measure accuracy is impacted.

Directions for disposal

Note: Vacuum should remain on until steps J to L are completed

J. Disconnect the connective tubing from Patient Port (2).
K. Ensure that all ports are capped off and sealed other than Vacuum Port (1).
L. Disconnect red vacuum tube from Vacuum Port (1) and cap that port securely.
M. Turn vacuum off.
N. Remove suction soft liner from outer canister.
O. Comply with facility's protocols and local, state and federal requirements when disposing of soft liners and contents.

Precautions should be taken when discarding this device and disposal of the device shall be made in accordance with applicable hospital or national regulations for biologically hazardous waste.

SETUP AND USE for Soft Liners : TANDEM

Port identification

(1) Vacuum Port (2) Patient Port (3) Tandem Port (4) Accessory Port

Directions for set up and use

In a tandem application canister orientation should follow a counter clockwise pattern.

A. Place soft liners in outer reusable canister on the canister carrier.
B. Ensure that the red manifold tubing is attached to each of the reusable outer canisters.
C. Connect tubing from the suction regulator (or other vacuum source) to the red manifold tubing underneath the holder.
D. Insert the male tandem tube connector into Tandem Port (3). Be certain connector is inserted fully into Tandem Port and secured in place.
E. Position the female tandem tube valve over the Patient Port (2) of the next canister (removing elbow connector) that the overflow will enter. In order to ensure a secure seal, press down firmly.
F. Repeat the tandem tube instruction above (Step E) until all the liners are connected.
G. Attach connecting patient suction tubing to the Patient Port (2) of first canister and attach adapted suction device (eg. Yankauer) to other end of the patient suction tubing.
H. Confirm all unused ports are securely capped off prior to initiating suction.
I. Ensure that the on/off valves for all canisters desired for use are in the "ON" position prior to initiating suction.
J. Once setup as instructed, ensure that the vacuum has been created and the soft liners are fully inflated.
K. The system is ready for use.

Caution: When using soft liners in tandem set-up, in order to avoid potential leakage during disconnect do not completely fill the last liner.

Using the soft liners with solidifier

Follow manufacturer's instructions.
It is not recommended to use solidifier prior to suction during tandem configuration, as fluid will solidify, volume expands and suction is altered in the tandem tube during procedures.
For best results, solidifiers should be added after the procedure and before shutting off the vacuum source to prevent risk of fluid leaking from the accessory port when the solidifier is added.

Directions for disposal

Note: Vacuum should remain on until steps L to Q are completed. Securely cap all ports prior to turning off suction.

L. Disconnect the patient tubing from the Patient Port (2) on the first soft liner (removing elbow connector) and cap the port.
M. Carefully remove tandem tube valve from the Patient Port (2) on the tandem port of the liner.
N. Secure the cap over the end of the tandem tube valve and then cap that Patient Port (2).
O. Repeat steps M and N for all remaining liners.
P. Disconnect the red vacuum tube coming from the outer canister from the Vacuum Port (1) and cap the port. Repeat this step until all red vacuum tubes are disconnected.
Q. Confirm all unused ports are securely capped off prior to turning off suction.
R. Turn off the vacuum.
S. Remove the liner from the outer canister carefully. **Note:** The liner/tandem tube assembly remain together.
T. Dispose according to hospital protocol. **Note:** The liner and tandem tube assembly are single use only.

Note: Liner contents are considered potentially hazardous and proper caution should be taken when handling.
For soft liners tandem setups, it is recommended that all used and unused disposable soft liners be discarded and subsequently replaced following each procedure as all contaminated by aerosols.

⚠ WARNINGS

- Do not reuse or sterilize. If the same patient undergoes long-term treatment the soft suction liner must be replaced at least every 24 hours.
- Do not subject the canister to active vacuum when not in use.
- Not intended for pleural/thoracic use.
- Inspect the device for the presence of cracks or damaged parts prior to use.
- Monitor vacuum gauge to assess for consistent vacuum pressure prior to and during use.
- Causes of premature shut-off include:
 - ◊ Bumping the canister or tilting at an angle can activate the shut-off mechanism inside the lid.
 - ◊ Suctioning agents that release gas, including hydrogen peroxide or povidone iodine may lead to the release of foam that can activate the shut off mechanism or even eject fluid from the canister.

Cleaning of reusable canister

Reusable outer canister can be wiped with disinfectant.

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