
REF ORVN12
REF ORVN15

en INSTRUCTIONS FOR USE
de GEBRAUCHSANWEISUNG
fr NOTICE D'UTILISATION
nl GEBRUIKSAANWIJZING
it ISTRUZIONI PER L'USO
es INSTRUCCIONES DE USO
pt INSTRUÇÕES PARA USO

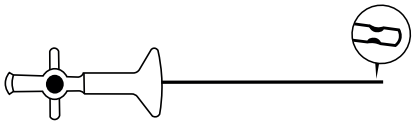STERILE EO

1639

Latex Free - Latexfrei - Sans latex -
Latexvrij - Senza lattice -
Sin Látex - Sem látex

Medical Device - Medizinprodukt - Dispositif médical -
Medisch hulpmiddel - Dispositivo Medico -
Producto Sanitario - Dispositivo médico

Wickimed (Huizhou) Medical Equipment Manufacturing Co.,
Ltd Tang Jiao Xing Wang Street, Lilin Town,
Zhongkai Hi-Tech Zone, Huizhou, Guangdong,
516035, P.R CHINA.

EC REP Shanghai International Holding Corp.
GmbH (Europe) Eiffestrasse 80,
20537 Hambourg, Germany

Distributed in Europe by: Medline International France SAS,
5 rue Charles Lindbergh, 44110 Châteaubriant, France

RG20WIC 2020-07

en VERESS NEEDLE
DISPOSABLE INSUFFLATION NEEDLE

Indication for use:
Veress Needles are intended to establish pneumoperitoneum prior to laparoscopy procedures.

Description:
Medline Veress Needles is a stainless steel needle with 14G outside diameter.
Veress Needles have a spring-loaded blunt stylet mechanism and comprise an ergonomic handle with a stopcock and a luer lock connector for insufflation tubing.

The inner stylet cannula retracts as the needle is pushed through the peritoneum and extends once the abdominal cavity has been reached.

Sterilised by ethylene oxide (EO).

Medline Veress Needle is available in two different lengths: 120 mm (ORVN12) and 150 mm (ORVN15).

Instructions for use:
1. Inspect packaging to ensure it is not broken or damaged.
2. Open packaging using aseptic technique, remove needle protection, and inspect overall condition and physical integrity.
3. Close and open the stopcock to ensure it works properly. Ensure it remains closed for insertion.
4. Make a small incision to insert the needle and inspect the handle to ensure that the safety red indicator slides back towards the stopcock. This indicates the inner stylet is retracting.
5. Once inserted, verify that the safety red indicator is back to its initial position (fully away from the stopcock), meaning the inner stylet is in full distal position; ready for insufflation.
6. Connect Veress Needle to the insufflation tubing for establishing pneumoperitoneum.

Contraindications:
Do **NOT** use if the endoscopic techniques are contraindicated.

Important:
- Single Use Only.
- Contents are sterile unless packaging was opened or is damaged.
- Any serious incident that has occurred in relation to the device should be reported to Medline and the competent authority of the Member State in which the user and/or patient is established.
- **Do NOT reuse and/or re-sterilise. Reuse of single-use devices creates a potential risk for the patient or user. It may cause contamination and/or impairment of function, which may lead to injury, illness and/or death of the patient.**

Precautions:
Do **NOT** attempt to insert the Veress Needle if the safety red indicator does not slide back.
Do **NOT** use the device if sterile conditions have been compromised or the packaging is damaged.
The product must be used exclusively by qualified medical staff. Thorough understanding of principles, instrumentation and techniques involving minimal invasive procedures is required to avoid possible injuries to both patient and medical staff.

Disposal instructions:
Discard after use. Precautions should be taken when discarding this device. Disposal of the device should be done in accordance with applicable hospital or national regulations for biologically hazardous waste.

Storage:
Store under light-proof and clean conditions. Avoid exposure to sunlight and relative humidity levels over 80%.

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