REF OR84

en **SPECIMEN CUP**

de **PROBENGEFÄSS**

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nl **SPECIMENCUP**

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pt **COPO PARA AMOSTRAS**

sv **PROVKOPP**

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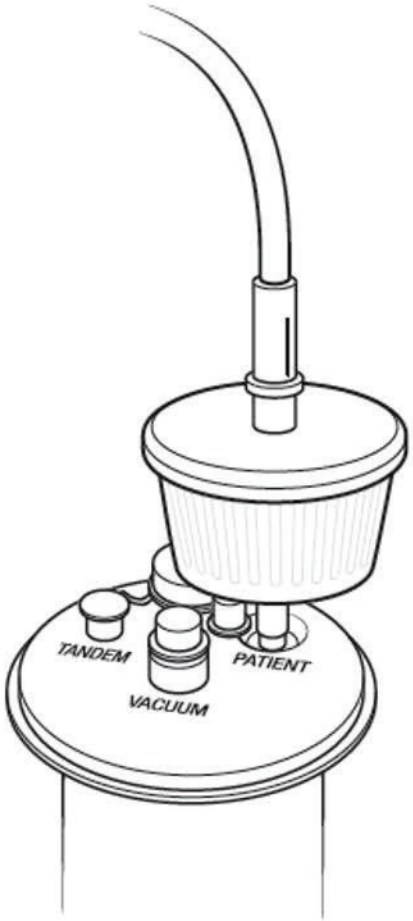
no **PRØVEKOPP**

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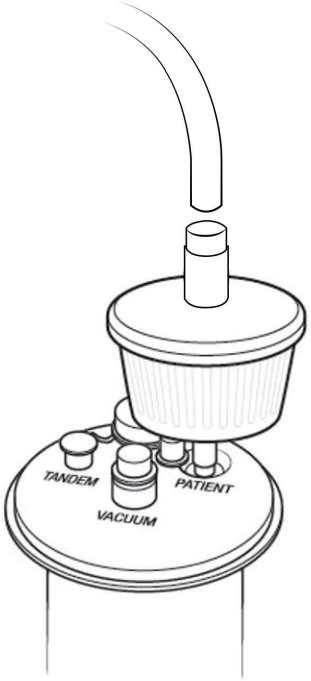
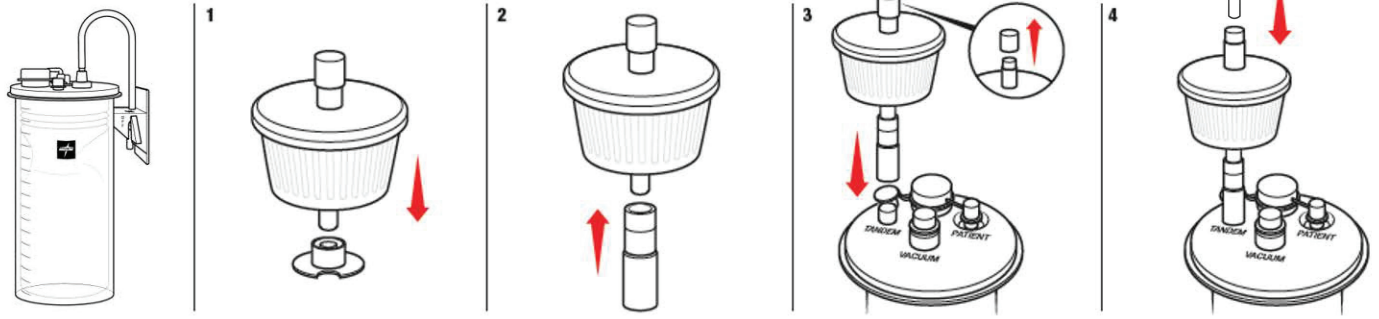
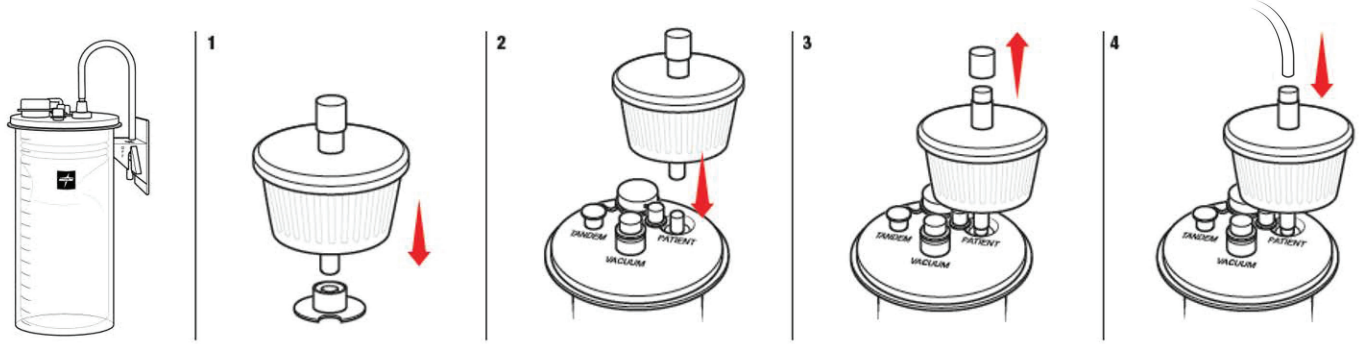
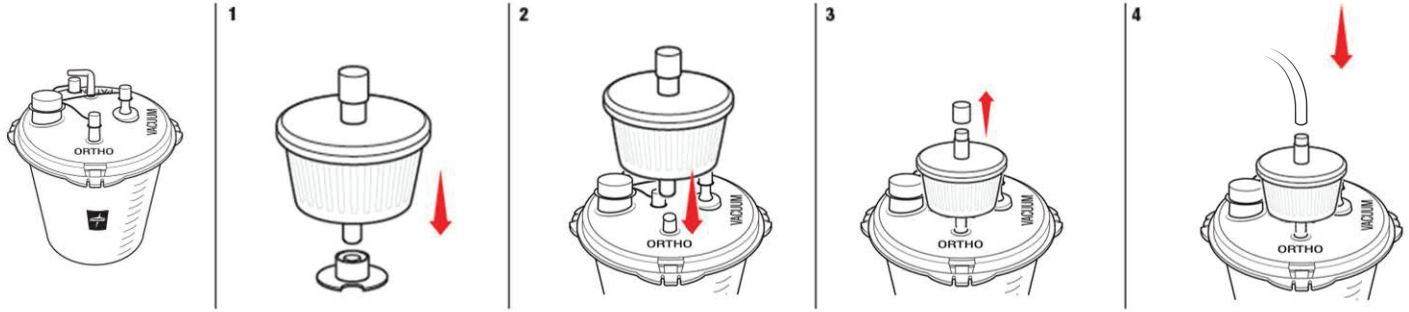
1
2
3
4
VACUUM
ORTHO

en INSTRUCTIONS FOR USE

There are two different adaptors for this product

- **For Medline canisters:** use the blue adaptor to connect the specimen cup.
- **For non-Medline canisters:** use the transparent adaptor to connect the specimen cup.
 NOTE: the transparent adaptor might not be compatible with all non-Medline canisters.

SET UP

1. Remove the two blue caps from the specimen cup. Hold onto the caps until the end of the procedure because they will be needed for the specimen cup removal process.
2. Connect the port on the bottom of the specimen cup directly to the port of the suction canister:
 - **SET 1: for US Soft Liners/EU Med-Soft Liners/AU Soft liners:** connect one end of the adaptor to the specimen cup and connect the other end of the adaptor to the port on the suction canister labeled “TANDEM” by pushing it into place.
 - **SET 2: for EU Med-Soft Liners/EU Flex One Liners/AU Soft Liners/ AU Semi-Rigid Liners:** connect the specimen cup directly to the port on the suction canister labeled “PATIENT” by pushing it into place.
 - **SET 3: for US Rigid Canisters/US Semi-Rigid Canisters/EU Med-Rigid Canisters/EU Med-Semi Rigid Canisters/AU Rigid Canisters:** connect the specimen cup directly to the port on the suction canister labeled “ORTHO” by pushing it into place.
3. Connect the patient tubing (10mm diameter) to the shorter port on the top of the specimen cup.
4. Check that all ports of the suction canister are closed.
5. Attach the vacuum tubing to the port of the suction canister labeled “VACUUM”.
6. Turn on the vacuum.

REMOVAL

1. Keep the vacuum on.
2. Disconnect the patient tubing from the shorter port on the top of the specimen cup.
3. Close the port on the top of the specimen cup with the blue cap.
4. Remove the specimen cup from the port of the suction canister.
5. Close the longer port on the bottom of the specimen cup with the blue base cap.
6. Close all of the ports on the suction canister except the one labeled ‘VACUUM’.
7. Disconnect the vacuum tubing from the port of the suction canister labeled “VACUUM” and close the port.
8. Turn off the vacuum.
9. Ensure both ports of the specimen cup are properly closed before sending to the laboratory.

Reuse of single use devices creates a potential risk for the patient or user. It may cause contamination and/or impairment of function; which may lead to injury, illness and/or death of a patient.

If a serious incident has occurred in relation to the device, it should be reported to Medline and the competent authority of the Member State in which the user and/or patient is established.

Directions for disposal: Precautions should be taken when discarding this device and disposal of the device shall be made in accordance with applicable hospital or national regulations for biologically hazardous waste.

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Medline Industries, Inc., Three Lakes Drive, Northfield, IL 60093 USA.
Medline International France SAS,