1. Page size : 8.5x11 inches 1. Margins: 0.25 in

HUDSON RCI®

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Specifications
Driving Gas: Compressed Air
Driving Gas Pressure: 345 kPa / 50 PSI
Driving Gas Flow: 6-8 LPM
Angles of Operation: 0-45° from vertical
Maximum Fill Volume: 6mL
Empty Mass of Nebulizer: 0.093kg

Specifications; Spezifikationen; Caractéristiques; Specificaties; Specifiche; Especificaciones; Especificações; Špecifikácie; Specifikace; Tekniska specifikationer; Tekniset tiedot; Specifikationer; Spesifikasjoner; Προδιαγραφές; Dane techniczne; Specifikacije; Specificații; Спецификации; Specifikációk; Teknik Özellikler; المواصفات

Driving Gas: Compressed Air; Treibgas: Druckluft; Entraînement du gaz : Air comprimé; Gebruikt gas: Gecomprimeerd gas; Guida a gas: Aria compressa; Gas impulsor: aire comprimido; Gás acionador: Ar comprimido; Hnací plyn: Stlačený vzduch; Hnací plyn: Stlačený vzduch; Drivgas: Tryckluft; Ajokaasu: Paineilma; Drivluft: Trykluft; Kjøregass: Trykkluft; Οδήγηση αερίου: Πεπιεσμένος αέρας; Gaz napędowy: Sprężone powietrze; Pogonski plin: Stisnjen zrak; Gaz de acționare: Aer Comprimat; Задвижващ газ: Въздух под налягане; Vezető gáz: Sűrített levegő; Sürücü Gaz: Basınçlı Hava; غاز التشغيل: هواء مضغوط


Driving Gas Pressure: 345 kPa / 50 PSI; Treibgasdruck: 345 kPa / 50 PSI; Pression d'entraînement du gaz : 345 kPa / 50 PSI; Druk gebruikt gas: 345 kPa / 50 PSI; Pressione guida a gas: 345 kPa / 50 PSI; Presión del gas impulsor: 345 kPa / 50 PSI; Pressão do gás acionador: 345 kPa / 50 PSI; Tlak hnacieho plynu: 345 kPa/50 PSI; Tlak hnacího plynu: 345 kPa / 50 PSI; Gastryck vid drift: 345 kPa / 50 PSI; Ajokaasupaine: 345 kPa / 50 PSI; Drivlufttryk: 345 kPa/50 PSI; Drivgasstrykk: 345 kPa / 50 PSI; Πίεση οδήγησης αερίου: 345 kPa / 50 PSI; Ciśnienie gazu w układzie napędowym: 345 kPa / 50 PSI; Tlak pogonskega plina: 345 kPa/50 PSI; Presiune de acționare a gazului: 345 kPa / 50 PSI; Налягане на задвижващия газ: 345 kPa / 50 PSI; Vezető gáznyomás: 345 kPa / 50 PSI; Sürücü Gaz Basıncı: 345 kPa / 50 PSI; ضغط غاز التشغيل: 345 كيلو باسكال/ رطل لكل بوصة مربعة

Driving Gas Flow: 6-8 LPM; Treibgasfluss: 6 - 8 LPM; Débit d'entraînement du gaz : 6-8 LPM; Stroom gebruikt gas: 6-8 LPM; Flusso guida a gas: 6-8 LPM; Caudal de gas impulsor: 6-8 LPM; Fluxo do gás acionador: 6-8 LPM; Prietok hnacieho plynu: 6 – 8 l/min.; Průtok hnacího plynu: 6–8 l/min; Drivgasflöde: 6-8 LPM; Ajokaasuvirtaus: 6- 8 LPM; Drivluftstrømning: 6-8 l/m; Strømmen av drivgass: 6-8 LPM; Ροή οδήγησης αερίου: 6-8 LPM; Przepływ gazu napędowego: 6-8 LPM; Pretok pogonskega plina: 6–8 LPM; Debit de acționare a gazului: 6-8 LPM; Поток на задвижващ газ: 6-8 LPM; Vezető gázáramlás: 6-8 LPM; Sürücü Gaz Akışı: 6-8 LPM; تدفق غاز التشغيل: 6-8 لتر/دقيقة

Angles of Operation: 0-45° from vertical; Betriebswinkel: 0-45° von der Senkrechten; Angles d'utilisation : 0-45° Par rapport à a verticale; Werkingshoeken: 0-45° verticaal; Angoli di funzionamento: 0-45° da posizione verticale; Ángulos de operación: 0-45° de la vertical; Ângulos de Operação: 0-45° em relação a vertical; Uhly prevádzky: 0-45° z vertikálnej polohy; Provozní úhel: 0-45° od svislé polohy; Driftvinklar: 0-45° från vertikal; Käyttökulmat: 0-45° pystysuorasta; Betjeningsvinkler: 0-45° fra lodret; Driftsvinkler: 0-45° fra vertikal; Γωνίες λειτουργίας: 0-45° από την κάθετη; Kąty działania: 0-45° w pionie; Koti delovanja: 0-45° navpično; Unghiuri de Operare: 0-45° din poziție verticală; Ъгли на работа: 0-45° от вертикала; Működési szögek: 0-45° a függőlegeshez képest; Çalışma Açıları: 0-45° dikeye göre; زوايا التشغيل: 0-45 درجة من الوضع الرأسي

Maximum Fill Volume: 6mL; Maximales Füllvolumen: 6mL; Volume maximum de remplissage 6 mL; Maximaal vulvolume: 6mL; Volume massimo di riempimento: 6mL; Volumen de llenado máximo: 6 mL; Volume de enchimento máximo: 6 mL; Maximálny objem náplne: 6 ml; Maximální objem naplnění: 6 ml; Maximal påfyllningsvolym: 6mL; Suurin täyttötilavuus: 6mL; Maksimal fyldningsvolumen: 6 ml; Maksimalt fyllvolum: 6mL; Μέγιστος όγκος πλήρωσης: 6mL; Maksymalna objętość napełnienia: 6 ml; Največji volumen polnjenja: 6mL; Volum Maxim de Umplere: 6mL; Максимален обем за пълнене: 6mL; Maximális töltési térfogat: 6mL; Maksimum Dolum Hacmi: 6 mL; أقصى حجم تعبئة: 6 مل

Empty Mass of Nebulizer: 0.093kg; Leere Masse des Verneblers: 0,093 kg; Masse à vide du nébuliseur : 0,093 kg; Lege massa van vernevelaar: 0,093 kg; Massa vuota del nebulizzatore: 0,093 kg; Masa del nebulizador vacío: 0,093 kg; Massa do nebulizador em vazio: 0,093 kg; Hmotnosť prázdneho nebulizéra: 0,093 kg; Hmotnost prázdného nebulizátoru: 0,093 kg; Nebulisatorns tomvikt: 0,093kg; Sumuttimen tyhjä massa: 0,093kg; Forstøverens tomme masse: 0,093kg; Tom masse av forstøveren: 0,093kg; Κενή μάζα νεφελοποιητή: 0,093kg; Pusta masa nebulizatora: 0,093 kg; Masa praznega nebulizatorja: 0,093 kg; Greutatea netă a Nebulizatorului: 0,093 kg; Празна маса на небулизатора: 0,093kg; A porlasztó üres tömege: 0,093kg; Nebulizatörün Boş Ağırlığı: 0,093 kg; الكتلة الفارغة لجهاز الرذاذ: 0.093 كجم

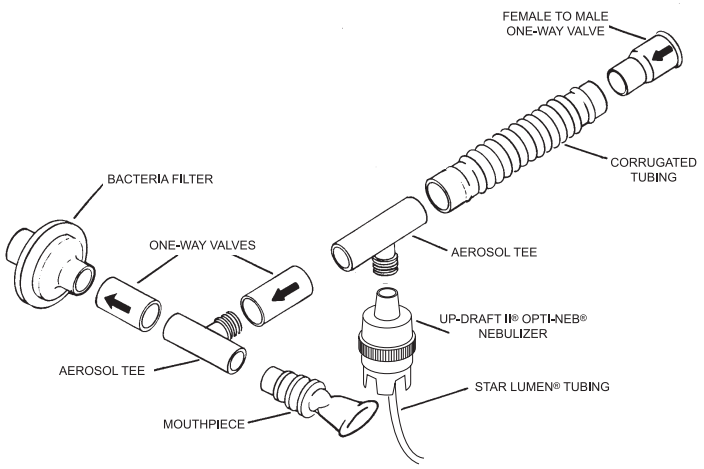BACTERIA FILTER; BAKTERIENFILTER; FILTRE À BACTÉRIES; BACTERIEFILTER; FILTRO PER BATTERI; FILTRO BACTERIANO; FILTRO DE BACTÉRIAS; FILTER BAKTERIÍ; BAKTERIÁLNÍ FILTR; BAKTERIEFILTER; BAKTEERISUODATIN; BAKTERIEFILTER; BAKTERIEFILTER; ΦΙΛΤΡΟ ΒΑΚΤΗΡΙΩΝ; FILTR BAKTERYJNY; BAKTERIJSKI FILTER; FILTRU DE BACTERII; АНТИБАКТЕРИАЛЕН ФИЛТЪР; BAKTÉRIUMSZŰRŐ; BAKTERİ FİLTRESİ; مرشح البكتيريا

ONE-WAY VALVES; EINWEGVENTILE; VALVES ANTI-RETOUR; EENRICHTINGSKRAAN; VALVOLE UNIDIREZIONALI; VÁLVULAS UNIDIRECCIONALES; VÁLVULAS UNIDIRECCIONAIS; JEDNOSMERNÉ VENTILY; JEDNOCESTNÉ VENTILY; ENVÄGSVENTILER; YKSITIEVENTTIILIT; ENVEJSVENTILER; ENVEISVENTILER; ΜΟΝΟΔΡΟΜΕΣ ΒΑΛΒΙΔΕΣ; ZAWORY JEDNOKIERUNKOWE; ENOSMERNI VENTILI; SUPAPE UNIDIRECȚIONALE; ЕДНОПОСОЧНИ КЛАПАНИ; EGYUTAS SZELEPEK; TEK YÖNLÜ VANALAR; صمامات أحادية الاتجاه

AEROSOL TEE; T-AEROSOLVERNEBLER; TEE AÉROSOL; AEROSOL STUK; MASCHERA PER AEROSOL; MASCARILLA PARA AEROSOL TEE; AEROSSOL T; AEROSÓLOVÝ T; AEROSOLOVY TEE ADAPTER; AEROSOL TEE; AEROSOLIN T-LIITIN; AEROSOL-T-STYKKE; AEROSOL T-SKJORTE; ΑΕΡΟΖΟΛ ΤΕΕ; KOSZULKA AEROZOLOWA; NASTAVEK ZA AEROSOL; CONECTOR T AEROSOLI; Т-ОБРАЗЕН АДАПТОР ЗА АЕРОЗОЛ; AEROSZOLOS T-ADAPTER; AEROSOL TE BORU; وصلة ثلاثية الشكل مزودة برذاذ

MOUTHPIECE; MUNDSTÜCK; EMBOUT BUCCAL; MONDSTUK; BOCCAGLIO; PIEZA BUCAL; PIEZA BUCAL; NÁUSTOK; NÁUSTEK; MUNSTYCKE; SUUKAPPALE; MUNDSTYKKE; MUNNSTYKKE; ΕΠΙΣΤΟΜΙΟ; USTNIK; USTNIK; MUȘTIUC; ПРИСТАВКА ЗА УСТА; CSUTORA; AĞIZLIK; قطعة فموية

FEMALE TO MALE ONE-WAY VALVE; WEIBLICH ZU MÄNNLICH EINWEGVENTIL; VALVE ANTI-RETOUR FEMELLE À MÂLE; VROUWELIJKE NAAR MANNELIJKE EENRICHTINGSKRAAN; VALVOLA UNIDIREZIONALE FEMMINA-MASCHIO; VÁLVULA UNIDIRECCIONAL HEMBRA-MACHO; VÁLVULA UNIDIRECCIONAL FÊMEA - MACHO; JEDNOSMERNÝ VENTIL ŽENSKÝ K MUŽSKÉMU; JEDNOCESTNÝ VENTIL ŽENSKÝ - MUŽSKÝ; ENVÄGSVENTIL HANE TILL HONA; NAARAS-UROS YKSITIEVENTTIILI; HUN TIL HAN ENVEJSVENTIL; HUNN TIL HANN ENVEISVENTIL; ΘΗΛΥΚΗ ΣΕ ΑΡΣΕΝΙΚΗ ΜΟΝΟΔΡΟΜΗ ΒΑΛΒΙΔΑ; ZAWÓR JEDNOKIERUNKOWY ŻEŃSKO-MĘSKI; ŽENSKI DO MOŠKI ENOSMERNI VENTIL; SUPAPĂ MAMĂ-TATĂ UNIDIRECȚIONALĂ; ЕДНОПОСОЧЕН КЛАПАН, ЖЕНСКИ КЪМ МЪЖКИ ПРЕХОД; BELSŐ-KÜLSŐ MENETES EGYUTAS SZELEP; DİŞİ ERKEK TEK YÖNLÜ VALF; صمام أحادي الاتجاه من أنثى إلى ذكر

CORRUGATED TUBING; WELLENSCHLAUCH; TUBULURE ANNELÉE; GEGOLFDE SLANG; TUBO CORRUGATO; TUBO CORRUGADO; TUBOS CORRUGADOS; VLNITÉ HADIČKY; VLNITÁ HADIČKA; KORRUGERAD RÖRLEDNING; AALLOTETTU PUTKI; KORRUGEREDE SLANGER; KORRUGERTE SLANGER; ΚΥΜΑΤΟΕΙΔΗΣ ΣΩΛΗΝΩΣΗ; RURKA KARBOWANA; VALOVITA GIBKA CEVKA; TUBULATURA CORUGATĂ; ГОФРИРАНИ ТРЪБИ; HULLÁMOSÍTOTT CSŐ; OLUKLU BORU; أنبوب ممو ج

AEROSOL TEE; T-AEROSOLVERNEBLER; TEE AÉROSOL; AEROSOL STUK; MASCHERA PER AEROSOL; MASCARILLA PARA AEROSOL TEE; AEROSSOL T; AEROSÓLOVÝ T; AEROSOLOVY TEE ADAPTER; AEROSOL TEE; AEROSOLIN T-LIITIN; AEROSOL-T-STYKKE; AEROSOL T-SKJORTE; ΑΕΡΟΖΟΛ ΤΕΕ; KOSZULKA AEROZOLOWA; NASTAVEK ZA AEROSOL; CONECTOR T AEROSOLI; Т-ОБРАЗЕН АДАПТОР ЗА АЕРОЗОЛ; AEROSZOLOS T-ADAPTER; AEROSOL TE BORU; وصلة ثلاثية الشكل مزودة برذاذ

UP-DRAFT II® OPTINEB® NEBULIZER; UP-DRAFT II OPTINEB VERNEBLER; NÉBULISEUR UP-DRAFT II® OPTINEB®; UP-DRAFT II® OPTINEB® VERNEVELAAR; NEBULIZZATORE UP-DRAFT II® OPTINEB®; NEBULIZADOR UP-DRAFT II® OPTINEB; NEBULIZADOR UP-DRAFT II® OPTINEB; UP-DRAFT II® OPTINEB® NEBULIZÁTOR; ROZPRAŠOVAČ UP-DRAFT II® OPTINEB®; UP-DRAFT II® OPTINEB® NEBULISATOR; UP-DRAFT II® OPTINEB® -SUMUTIN; UP-DRAFT II® OPTINEB®-FORSTØVER; UP-DRAFT II® OPTINEB®-FORSTØVER; ΝΕΦΕΛΟΠΟΙΗΤΗΣ OPTINEB® UO-DRAFT II®; NEBULIZATOR UP-DRAFT II® OPTINEB®; NEBULIZATOR UP-DRAFT II® OPTINEB®; NEBULIZOR UP-DRAFT II® OPTINEB®; ИНХАЛАТОР UP-DRAFT II® OPTINEB®; UP-DRAFT II® OPTINEB® PORLASZTÓ; UP-DRAFT II® OPTINEB® NEBULİZATÖR; ®UP-DRAFT II® OPTINEB جهاز استنشاق/إصدار الرذاذ

STAR LUMEN® TUBING; SCHLAUCH MIT STERNLUMEN; TUBULURE STAR LUMEN; STAR LUMEN® SLANG; TUBO CON LUME A STELLA; TUBO STAR LUMEN; TUBAGEM STAR LUMEN; HADIČKY STAR LUMEN®; HADIČKA STAR LUMEN; STAR LUMEN® RÖRLEDNING; STAR LUMEN® -LETKUT; STAR LUMEN®-SLANGE; STAR LUMEN®-TUBE; ΣΩΛΗΝΩΣΗ STAR LUMEN®; RURKA STAR LUMEN®; CEVI STAR LUMEN®; TUBULATURĂ STAR LUMEN®; ТРЪБА STAR LUMEN®; STAR LUMEN® CSŐ; STAR LUMEN® HORTUM; ®STAR LUMEN أنبوب

en Intended purpose

Medline Hudson RCI Small Volume Nebulizers are intended to be used for intermittent aerosol therapy to deliver prescribed solutions for inhalation. Intended for use on pediatric or adult patients. Intended for use by medical personnel, lay persons, or patients that have been prescribed aerosol therapy. For use in acute and non-acute care settings. Intended for single patient, single use.

Directions for Use

1. Remove the ISO-NEB from its package and be sure all components are intact.
2. Rotate the nebulizer into upright position for right- or left-handed use.
3. Verify that the arrows on one-way valves are oriented correctly: two on the inspiratory limb toward the mouthpiece and one on the expiratory limb toward the particle filter.
4. Remove the nebulizer from the tee.
5. Add prescribed medication to the nebulizer jar through the opening at the top of the cap (6 mL maximum).
6. Reattach the nebulizer to the tee.
7. Attach tubing by pushing end connector directly onto compatible gas outlet. Push the other end onto the fitting under the nebulizer.
8. Adjust the gas source to prescribed flow rate (6 to 8 LPM).
9. Check for secure connections by confirming flow of aerosol through device.
10. Proceed with treatment as prescribed. Exhale through the mouthpiece.
11. Cover mouthpiece to divert flow toward filter when taking breaks during treatment.

Disposal instructions

Discard after single application. Precautions should be taken when discarding this device and disposal of the device shall be made in accordance with applicable hospital or national regulations for biologically hazardous waste.

Warnings

- Do not use with Anaesthetic or Ventilator breathing systems
- Using a solution other than 0.083% Albuterol Sulfate, in particular a suspension and/or high-viscosity solution, can alter the particle size distribution curve, the mass median aerodynamic diameter (MMAD), aerosol output, and/or aerosol output rate, as compared to values disclosed above.
- Materials used in the nebulizer may not be compatible with solutions, suspensions, or emulsions that have not been evaluated.
- Non-sterile. Do not sterilize prior to use.
- Do not use if device is soiled or damaged.
- Prior to use, ensure compatibility with gas outlet. Standard connectors are compatible with 5-7mm male fittings and universal connectors are compatible with male DISS 1240 fittings or 5-7mm male fittings.
- Prior to use, ensure all connections are secure. Do not over tighten or force connections. Do not alter device. Improper connections may result in disconnection or leakage during use, which may result in performance degradation and increase the risk of respiratory distress.
- Do not attach the outlet of the nebulizer directly to an endotracheal tube, tracheostomy tube or other device that does not offer an alternate pathway for exhalation.
- Single Use Only: Do not reuse, reprocess, or sterilize. Reuse of device creates a potential risk of serious injury and/or infection which may lead to death. Residuals from medication or cleaning agents retained within the device may result in contamination, dangerously increased breathing resistance and impaired valve function.
- If a serious incident has occurred in relation to the device, it should be reported to Medline and the competent authority of the Member State in which the user and/or patient is established.

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Manufactured in Mexico by:
Medline Industries, LP, Three Lakes Drive, Northfield, IL 60093 USA.

2797 EC REP Medline International France SAS,
5 rue Charles Lindbergh,
44110 Châteaubriant, France

CH REP **CH Rep and importer:** Medline International Switzerland SARL, c/o MN & Associés SA, 1 Place de Longemalle, 1204 Genève, Switzerland

UK CA 0086 **UK representative and importer in UK:** Medline Industries Ltd, WA4 6HL, UK

Medline International B.V.,
Nieuwe Stationsstraat 10,
Arnhem, 6811 KS, Netherlands

Distributed in Australia by:
Medline International Two Australia Pty. Ltd., Sydney, Australia

RA25145 2025-1
HNL02196