en

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en HEELMEDIX® heel protector is single patient use only. Product is used on patients while in bed to help prevent pressure ulcers.
Warning: Patient must not walk or stand while wearing the HEELMEDIX® heel protector. Reuse on multiple patients may lead to infection.
For quick cleaning: Use a damp cloth to wipe down the inside of the HEELMEDIX® boot, and air dry. **For a thorough cleaning:** Secure the HEELMEDIX® boot straps and put it in its mesh laundry bag. Washable in machine (38°C/100°F maximum). Tumble dry on low/delicate cycle (38 °C/100 °F maximum).
If a serious incident has occurred in relation to the device, it should be reported to Medline and the competent authority of the Member State in which the user and/or patient is established.
Disposal instructions: Precautions should be taken when discarding this device and disposal of the device shall be made in accordance with applicable hospital or national regulations for biologically hazardous waste.

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Medline International France SAS, 5 rue Charles Lindbergh, 44110 Châteaubriant, France
Distributed in Australia by: Medline International Two Australia Pty. Ltd.,Sydney, Australia

CH REP **CH Rep and importer:** Medline International Switzerland SARL, c/o MN & Associés SA, 1 Place de Longemalle, 1204 Genève, Switzerland

UK CA **UK representative and importer in UK:** Medline Industries Ltd, WA4 6HL, UK.


RA25JRH 2025-01

HeelMedix®

Heel Protector, Instructions for use

de **Fersenschutz,** Gebrauchsanweisung
fr **Protecteur de talon,** Mode d'emploi
nl **Hielbeschermer,** Gebruiksaanwijzing
it **Stivaletto,** Istruzioni per l'uso
es **Protector de talón,** Instrucciones de uso
pt **Protetor de calcanhar,** Instruções para uso
sk **Chránič päty,** Návod na použitie
cs **Chránic paty,** Návod k použití
sv **Hälskydd,** Instruktioner för användning
fi **Kantapääsuojus,** Käyttöohjeet
da **Hælbeskytter,** Instruktion
no **Beskytter, Bruksanvisning,** Bruksanvisning
el **Προστατευτικό φτέρνας,** Οδηγίες χρήσης
pl **Ochraniacz na pięty,** Instrukcja obsługi
sl **Ščitnik pete,** Navodila za uporabo
ro **Protecție pentru călcâi,** Instructiuni de folosire
bg **Протектор за пети,** Инструкции за употреба
et **Kannakaitse,** Kasutusjuhised
hu **Sarokvédő,**

1

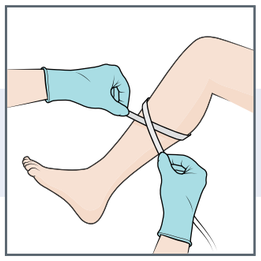en
Measure calf circumference, using measuring tape and sizing chart provided.

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