Probe Covers

Instructions for Use

Intended Use / Indications for Use
A device intended to be used as a physical barrier for protection against the effects of environmental exposure (e.g., body fluids, gels) and/or tomaintain the required hygienic level of various diagnostic or surgical procedures which utilize a transducer (probe).
This is a single use device. [Not to be reprocessed.]
The devices are further indicated to be used in or on:

- Intact or compromised skin
- Natural body orifices (e.g., vagina, rectum, esophagus)
- Intra-operatively (e.g., on internal organs)

Intended Patient Population
Probe Covers have no limitation for any definition of the patient group and their estimated state/condition. Probe Covers can be used at the discretion of the healthcare professional.

Contraindications

- Intraoperative Probe Covers are not specifically intended for use in direct contact with the heart or central circulatory system or the central nervous system.

Warnings

- Visually inspect for breaches of packaging integrity (such as tears, cuts, open seals, etc.) prior to use.
- Do not use if package is opened or damaged. If damage is found, call your Medline-Microtek representative.
- Report serious incidents to legal manufacturer/importer and competent authority.

Precautions

- Secure distal tip to probe with rubber bands or ties to prevent gel from migrating away from tip.
- Use care when handling the Probe Cover in close proximity to sharps, cautery devices, and other medical devices or procedures that may damage the integrity of the Probe Cover.


Product / Device Description
Probe Covers are composed of natural or synthetic materials and come in a variety of sizes and shapes to allow these devices to properly fit equipment with diverse shapes and uses. Probe Covers are used in a variety of clinical and surgical settings by healthcare professionals. The devices cover a transducer / probe and, in some instances, the cables connected to the transducers / probes. These coverings offer a physical barrier between the transducer and the patient.

This product may be supplied with ultrasound gel, rubber bands and/or ties.

Reuse Warning Statement
For single patient use only. Do not reuse, reprocess or resterilize. Reuse, reprocessing or resterilization may compromise the structural integrity of the device and/or lead to device failure which, in turn, may result in patient injury, illness or death. Reuse, reprocessing or resterilization may also create a risk of contamination of the device and/or cause patient infection or cross-infection, including, but not limited to, the transmission of infectious disease(s) from one patient to another. Contamination of the device may lead to injury, illness or death of the patient.

Handling and Storage

- Keep away from sunlight
- After use, dispose of product and packaging in accordance with hospital, administrative and/or local government policy.

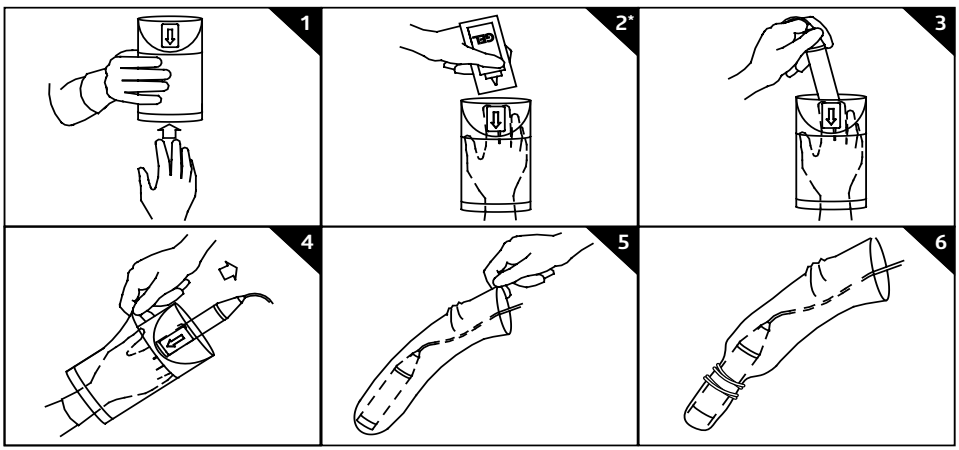**Step #2 is optional when gel is not included. / Стъпка #2 е опционална, когато гелът не е включен. / Krok č. 2 je volitelný, pokud není součástí balení gel. / Trin nr. 2 er valgfrit, hvis der ikke medfølger gel. / Schritt 2 ist optional, wenn kein Gel enthalten ist. / Το βήμα αρ. 2 είναι προαιρετικό όταν δεν χρησιμοποιείται γέλη. / El paso 2 es opcional si no se incluye gel. / 2. sämm on valikuline, kui geel ei ole tootega kaasas. / Vaihe 2 on valinnainen, kun geeliä ei ole toimitettu tuotteen osana. / L'étape 2 est facultative si le gel n'est pas inclus. / Korak br. 2 nije obavezan kada gel nije priložen. / A 2. lépés opcionális, amikor nincs mellékelve gél. / Skref 2 er valkvæt þegar hlaup fylgir ekki. / Il passaggio n. 2 è facoltativo quando il gel non è in dotazione. / Jei gelis nepridedamas, 2 žingsnis yra pasirinktinis. / Ja gels nav iekļauts, 2. darbība nav obligāta. / Stap 2 is optioneel wanneer er geen gel meegeleverd is. / Trinn 2 er valgfritt når gel ikke følger med. / Krok 2 jest opcjonalny, gdy żel nie należy do zestawu. / O passo n.º 2 é opcional quando o gel não está incluído. / Pasul nr. 2 este opțional atunci când gelul nu este inclus. / Если гель не входит в комплект, шаг № 2 необязателен. / KrokKrok č. 2 je voliteľný, keď gél nie je dodaný. / Če gel ni priložen, 2. korak ni obvezen. / Korak br. 2 je opcioni kada se gel ne koristi. / Steg nr 2 är valfritt när gel inte medföljer. / Jélin ürünle birlikte verilmediği durumlarda 2. adım isteğe bağlıdır.*

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©2025 Medline Industries LP. All rights reserved.
Microtek Medical LLC
Northfield, IL 60093 U.S.A.
Microtek Medical B.V.
Medline Industries Ltd