Endocavity Probe Cover

de **Ultrschall- sondenbezug endokavitär**
fr **Housse de sonde endocavitaire**
nl **Endocavity sonde cover**
it **Coprisonda endocavitaria**
es **Funda de ecógrafia endocavitaria**
pt **Cobertura da sonda de endocavidade**
sk **Návlek na endokavitálnu sondu**
cs **Krytí endokavitární sondy**
sv **Endoskopiskt probeskydd**
fi **Endoskooppisen anturin suoja**
da **Endoskopisk probeafdækning**
no **Endoskopisk probebeskyttelse**
el **Κάλυμμα ενδοσκοπικού καθετήρα**
pl **Osłona do sondy**
sl **Pokrov za endokavitetno sondo**
ro **Huse pentru sonde endocavitate**
bg **Предпазно покритие за ендокавитарна сонда**
et **Endoskoopilise sondi kate**
hu **Endokavitális szonda védőburok**
lv **Endokavitālās zondes pārvalks**
lt **Endokavitalinio zondo apdangalas**
hr **Navlaka za endokavitalnu sondu**
tr **Endokavite prob kapağı**
ar **غطاء مسبار الفحص التصويري داخل تجويف الجسم**

Karex Industries Sdn. Bhd.
PTD 7906 & 7907, Taman Pontian Jaya,
Batu 34, Jalan Johor, 82000 Pontian, Johor, Malaysia

EU REP Advena Limited Tower Business Centre 2nd Flr.,
Tower Street, Swatar, BKR 4013 Malta

CH REP **CMC Medical Devices GmbH**
Rigistrasse 3, 6300 Zug, Switzerland


CH Importer : Medline International Switzerland SARL,
c/o MN & Associés SA, 1 Place de Longemalle,
1204 Genève, Switzerland

UK importer: Medline Industries Ltd WA4 6HL UK

Medline International B.V.,
Catharijnesingel 47, 3511 GC Utrecht, The Netherlands

REF ICE28295L
REF ICE33230L

K-IFU-03/2 (v2)
LAB252629 2025-06

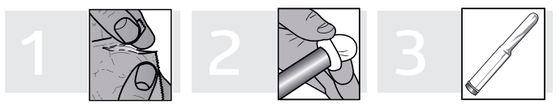en

It is intended to be used by healthcare provider trained in the use of ultrasound.
There are no unique skills or user abilities required for protective transducer cover use.
User groups may include, but not limited to:
Anesthesiologists, Department Nursing staff, Emergency Medicine physicians, Endocrinologists, Interventional Radiologists & Interventional Radiology Technologists, Radiologists, Sonographers, Surgeons, Surgical Technicians, and other clinicians.
Prevent cross-contamination of probes and consequently reduce the risk of cross of transducer probe during endocavity. (i.e transvaginal and transrectal) ultrasound procedure.
Transducer cover is a medical device used by healthcare provider. It is intended for barrier protection during invasive medical examination and diagnosis using transducer probe. Transducer cover use to prevent cross-contamination of Transducer probe during endocavity (i.e. transvaginal or transrectal) ultrasound procedure.
Before use, you should be trained in ultrasonography.
DISPOSE of single use components as infectious waste.
DO NOT REUSE transducer covers. Reuse or reprocessing may create a risk of contamination of the device and cause patient infection or cross-infection.
USE ONLY WATER BASED LUBRICANT, avoid petroleum or mineral oil based lubricant that may harm the covers.
Disinfect the transducer probe after each examination could reduce the probability of cross-contamination.
Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority.
This Product Contains Natural Rubber Latex Which May Cause Allergic Reactions.
The device was sterilized using irradiation and it is packed in a single sterile barrier system.
DO NOT USE the transducer covers if the individual package is obviously damaged or unintentionally open before use.
Discard that covers and use a new one from an undamaged package.
1. Carefully tear off the foil pack at the location where the package can be opened. Make sure not to damage the cover.
2. Remove the cover using proper sterile technique.
3. During removal from foil, avoid damage to the cover by sharp edges, fingernails and etc.
4. Place the cover over the transducer. Pull cover tightly over the transducer to remove wrinkles and air bubbles, taking care to avoid puncturing the cover.
5. Inspect the cover to ensure no holes or tears.

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