MICROTEK

Video Camera Cover

Instruction for use

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
Medline International France SAS,
5 rue Charles Lindbergh, 44110 Châteaubriant, France
Distributed in Australia by:
Medline International, Two Australia Pty. Ltd., Sydney, Australia

CE 2797 UKCA 0086 **UK Representative and importer in UK:**
Medline Industries Ltd, WA4 6HL, UK

CH REP **CH Rep and importer:** Medline International Switzerland Sàrl, c/o MN & Associés SA, 1 Place de Longemalle, 1204 Genève, Switzerland

LAB260128-2026-01

REF ICE4240

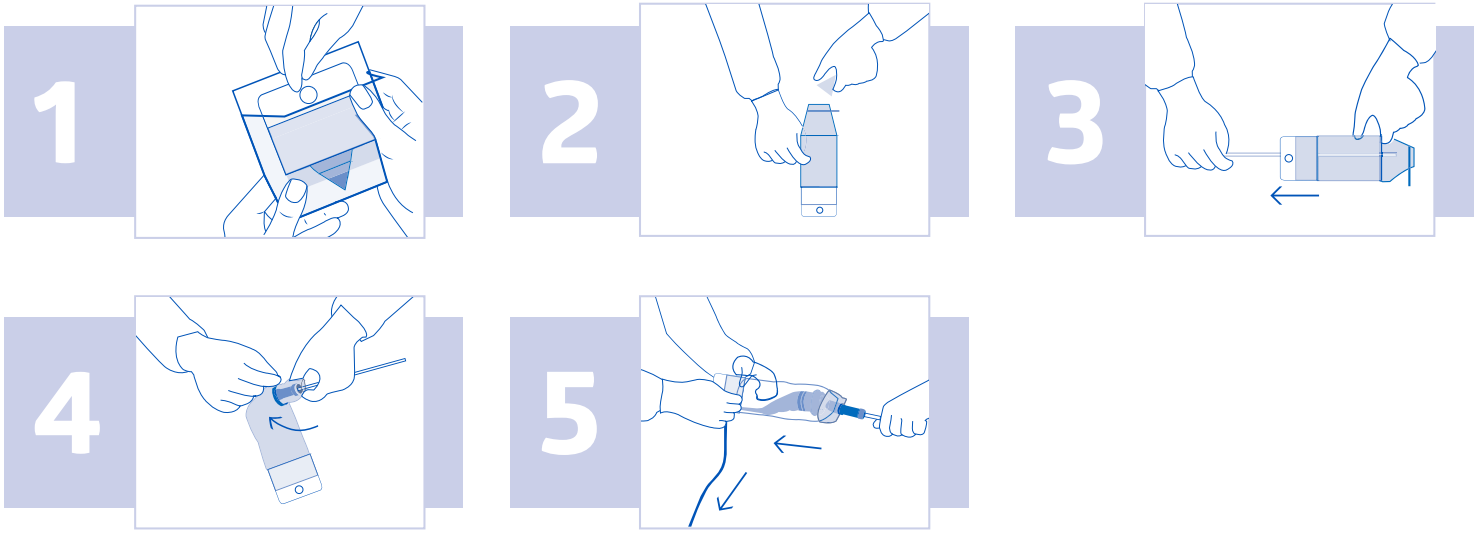en

Intended Use and Clinical Benefits:

Video camera covers are designed to cover cameras that are present in the operating room. These covers that are sterilized prior use, help to maintain the sterile field and prevent potential contamination for the patient undergoing a surgery.

Directions for Use:

Figure 1: Peel up individual pouch and remove the cover.
Figure 2 and 3: Insert endoscope through rigid applicator of folded drape and triangle cut tip up to its camera connection side.
Figure 4: Wrap the portion of drape that covers camera connection side of endoscope with adhesive tape pre-attached.
Figure 5: After connecting the camera to the endoscope, pull and apply the protective cover over the remaining camera cabling. Do not hesitate to ask for help.

Disposal Instructions:

Precautions should be taken when discarding this device and disposal of the device shall be made in accordance with applicable hospital or national regulations for biologically hazardous waste.

Contraindications and Warning:

- Do not use if package is opened or damage.
- Do not expose to a flame or high intensity heat sources.
- Reuse of single use devices creates a potential risk for the patient or user. It may cause contamination and/or impairment of function, which may lead to injury, illness and/or death of a patient.
- If a serious incident has occurred in relation to the device, it should be reported to Medline and the competent authority of the Member State in which the user and/or patient is established.

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