SureSite® SELECT Dressing


REF **VAD1685-INT**

(EN) SureSite® SELECT Dressing
(BG) Обяснение на символите на опаковката
(CS) Legenda značek na obalu
(DA) SureSite® SELECT-forbinding
(DE) SureSite® SELECT-Verband
(EL) Επίδεσμος SureSite® SELECT
(ES) Apósito SureSite® SELECT
(ET) Lehel olevate sümbolite selgitus
(FI) SureSite® SELECT -sidos
(FR) Pansement SureSite® SELECT
(GA) Míniú na Siombailí ar an bPacáiste
(HR) Objašnjenje simbola na pakiranju
(HU) A csomagoláson található szimbólumok magyarázata
(IT) Cerotto per medicazione SureSite® SELECT
(LT) Ant pakuotės esančių simbolių paaiškinimas
(LV) Uz iepakojuma esošo simbolu skaidrojums
(MT) SureSite® SELECT Dressing
(NL) SureSite® SELECT-verband
(PL) Opatrunek SureSite® SELECT
(PT) Penso SureSite® SELECT
(RO) Pansament SureSite® SELECT
(SK) Obväz SureSite® SELECT
(SL) Obveza SureSite® SELECT
(SV) SureSite® SELECT förband
(AR) SureSite® SELECT ضمادة
(JA) SureSite® SELECT ドレッシング材
(NB) SureSite® SELECT kompress
(TR) SureSite® SELECT Pansuman Örtüsü

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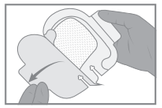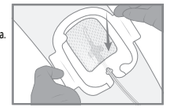
(EN) **Indications:** As a primary dressing, to secure vascular access line sites. Can be left in place, if clean and dry, for up to 7 days. Keep dry during use. Remove/change at first sign of soiling or moisture.
If a serious incident has occurred in relation to the device, it should be reported to Medline and the competent authority of the Member State in which the user and/or patient is established.
General Precautions & Observations: Can be used as a primary or secondary dressing. Should any signs of redness or irritation appear, discontinue use and consult a healthcare professional.
Used devices should be disposed of in proper waste containers per hospital protocol.
Caution: Federal (USA) law restricts this device to sale by or on the order of a physician.
Distributed and imported by

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MD

UK CA 0086 **UK approved body:** BSI Assurance UK Ltd, Kitemark Court Davy Avenue, Knowlhill, Milton Keynes MK5 8PP, UK (Approved body number 0086)
UK representative and importer in UK: Medline Industries Ltd, WA4 6HL, UK

Australian Sponsor:
Emergo Australia Level 20 Tower II,
Darling Park 201 Sussex Street
Sydney, NSW 2000 Australia

Manufactured by;
Centurion Medical Products, LP
100 Centurion Way
Williamston, MI 48895 USA

Manufactured for Medline Industries, LP
Medline and SureSite are registered trademarks of Medline Industries, LP
Three Lakes Dr. Northfield, IL 60093 USA

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EC REP Medline International France SAS