Empty Humidifier 27.6 and 41.4kPa (4 and 6PSI) Pressure Relief Valve

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
REF 3260CE (2-15LPM)

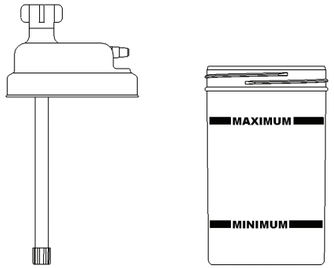en

Intended use:
Adds humidity in water vapor form to respiratory gases delivered to the patient. Medline Hudson RCI Empty Humidifiers are intended for use on general populations, which may include pediatric and adult. To be used with oxygen as their main gas source. It is a single patient-multiple use device used by qualified medical personnel, home healthcare providers, or patients under the general supervision of a qualified medical personnel in both acute and non-acute care setting such as hospitals or homecare environments. These products are intended to be used for a maximum of 30 days.

Instructions for Use:
1. Fill the humidifier jar to the maximum fill line with only sterile water for inhalation, distilled, or purified drinking water. Tightly secure the lid to the jar.
2. Connect the humidifier lid wing nut to the oxygen flowmeter outlet.
3. Attach oxygen tubing (not supplied) to the humidifier lid outlet stem. Ensure all connections are secure.
4. Turn the oxygen source to the prescribed flow.
5. Test the pressure relief valve by occluding oxygen tubing at 3LPM. Listen for an alarm condition. If alarm does not trigger, discard the device.

Cleaning Instructions:
NOTE: Clean when emptied after 3 days of use.
1. Unscrew the lid from the jar.
2. Dilute a ratio of 1mL of dishwashing detergent into 1 L of warm tap water in a clean tub or bowl.
3. Immerse the lid and jar in the soapy solution and while soaking gently wipe the surface of the components with a non-linting wipe without disassembling the components.
4. Rinse the lid and jar under warm water.
5. Dry the lid and jar with a non-linting wipe.
6. Visually inspect the lid and jar to ensure it is clean and that the components are properly attached. Repeat cleaning steps if necessary.
7. Screw the lid into the jar.

NOTE: if scaling (mineral deposit) occurs, dispose of humidifier.
CAUTION: An alarm whistle from the pressure relief valve will activate if the jar accumulates too much pressure. Inspect the system for excessive blockage or a kinked tube. To test the pressure relief valve, occlude oxygen tubing at 3 LPM. Listen for an alarm condition.

Warnings:
· Do not use in the presence of an ignition source.
· If a serious incident has occurred in relation to the device, it should be reported to Medline and the competent authority of the Member State in which the user and/or patient is established.
· Oxygen supplies of 68.9kPa (10PSI) or less using low flow (less than 2LPM) may affect the audible alarm.
· Do not use if device is soiled or damaged.

Disposal instructions:
Precautions should be taken when discarding this device and disposal of the device shall be made in accordance with applicable hospital or national regulations for biologically hazardous waste.

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