
HEEL RAISER PRO®
Heel Protector

Instructions for use

- de **Fersenschutz,** Gebrauchsanweisung
- fr **Protecteur de talon,** Mode d'emploi
- nl **Hielbeschermer,** Gebruiksaanwijzing
- it **Proteggi-tallone,** Istruzioni per l'uso
- es **Protector de talón,** Instrucciones de uso
- pt **Protetor de calcanhar,** nstruções para Uso
- sk **Chránič päty,** Návod na použitie
- cs **Chrániče na paty,** Návod k použití
- sv **Hälskydd,** Instruktioner för användning
- fi **Kantapäänsuoja,** Käyttöohjeet
- da **Hælbeskytter,** Instruktion
- no **Hælbeskytter,** Bruksanvisning
- el **Προστατευτικό φτέρνας,** Οδηγίες χρήσης
- pl **Ochraniacz na pięty,** Instrukcja obsługi
- sl **Ščitnik za peto,** Navodila za uporabo
- ro **Protecţie pentru,** Instructiuni de folosire
- bg **Протектор за пети,** Инструкции за употреба
- et **Kannakaitse,** Kasutusjuhised
- hu **Sarokvédő,** Használati utasítás
- lt **Apsauginis kulno įtaisas,** Naudojimosi instrukcijos
- hr **Štitnik za petu,** Upute za uporabu
- tr **Topuk Koruyu,** Kullanim talimatlari

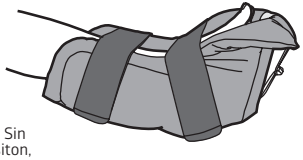MD

Latex Free, Latexfrei, Sans latex, Latexvrij, Senza lattice, Sin látex, Sem látex, Bez latexu, Bez Latexu, Latexfri, Lateksiton, Latexfri, Lateksfri, Χωρίς λάτεξ, Bezlateksowe, Brez lateksa, Fără latex, Несъдържащ латекс, Lateksivaba, Latex mentes, Sudėtyje nėra latekso, Bez lateksa, Lateks içermez

REF MDTE823200CS

en

HEEL RAISER PRO® is designed to off-load the heel for protection and treatment of heel pressure ulcers.

Benefits:
Helps control foot drop
Offers pressure relief by offloading the heel
Open design allows for easy observation and helps promote air-flow
Top folds down for more air-flow to toes or smaller feet

Directions for use:
1. Place lower leg on the long side of the boot with heel above the opening.
2. Secure straps around the boot. Do not wrap so tightly that it restricts circulation. (make sure heel is elevated).
3. Optional: Top can be folded down (see illustration) for more airflow or smaller feet.

DVT Tubing should be run through top two straps, not through the heel opening.
For quick cleaning: Use a damp cloth to wipe down and allow to air dry. For a thorough cleaning: Secure straps and put into mesh laundry bag (included). Washable in machine (38°C/100°F maximum). Tumble dry on low/delicate cycle (38°C/100°F maximum).
Warning: Patient must not walk or stand while wearing the HEEL RAISER PRO®.

For single patient use only. Reuse on multiple patients may lead to infection. If a serious incident has occurred in relation to the device, it should be reported to Medline and the competent authority of the Member State in which the user and/or patient is established.

Disposal instructions: Precautions should be taken when discarding this device and disposal of the device shall be made in accordance with applicable hospital or national regulations for biologically hazardous waste.

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Medline International France SAS, 5 rue Charles Lindbergh, 44110 Châteaubriant, France
Medline Industries Ltd, WA4 6HL, UK.