

HUDSON RCI®

REF 41650

Water Trap

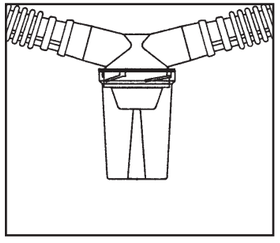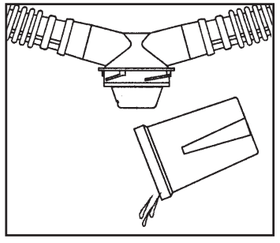Rx ONLY

Teleflex MEDICAL

3015 Carrington Mill Blvd
Morrisville, NC 27560 USA
USA: (866) 246 6990
International: +1 919 544 8000

CE 1639

EC REP
Teleflex Medical
IDA Business and Technology Park
Dublin Rd.,
Athlone, Co. Westmeath
Ireland

© Copyright 2006, 2010, 2011, 2012, 2017, 2018 Teleflex

L02804 R02 (2021-04)

Hudson RCI is a registered trademark of Teleflex Incorporated or its affiliates.

en

Water Trap

Not made with natural rubber latex.

Do not use on multiple patients
(Single Patient Use Only)

Reference	Lid Connection	Jar Size
41650	22 mm M	70 cc

Instructions for Use:

1. Thread the jar onto the lid assembly with a quarter-turn clockwise.
2. Connect the water trap in-line with the appropriate size tubing. The water trap functions best when placed at the lowest point in the breathing circuit to assist drainage.
3. If water is not draining, tilt and/or tap the jar to break surface tension.
4. When the jar is full, remove with a quarter turn anticlockwise and empty the accumulated moisture. The valve in the lid assembly will seal the breathing circuit while the jar is emptied.

Warning:

Single patient use: The use of products intended for single patient use on more than one patient creates a potential risk of serious injury and/or infection which may lead to death. Reprocessing of medical devices intended for single patient use only may result in degraded performance or a loss of functionality.

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