MedEnvision
EFFICIENCY THROUGH STANDARDIZATION

EsySuit Supine

English – Instruction for use

1. Place the approximation field under the patient's legs and keep legs up. Open the EsySuit and put both hands in the apertures (label 1). Slide both hands over the labels and insert as far as possible.
2. Slide with the palm of your hands over the toes and the instep of the feet. Turn your hands 180° and grab the ankles.
 ! Don't grab the toes!
3. Open the side panels and hold with both hands (label 2). Pull towards the patient's head unto the level of the pelvis (label 3).
4. Open the EsySuit (label 4) and unfold it under the patient's legs (label 5). Bring the legs down.
5. Unfold the EsySuit upwards (label 6) until the incision window is entirely exposed.
6. Keep the EsySuit tight and remove the protective backing of the incision window (label 7).
7. Position the adhesive incision foil on the surgical field. Unfold the EsySuit further to the patient's head and both sides of the table.
 Tip: when marking the anatomical landmarks before draping, it is easier to apply the incision foil in the right place.
8. Put sterile jerseys or bandages on both legs (label 8). To prevent perforation, stick the reinforcement to create a safe zone.

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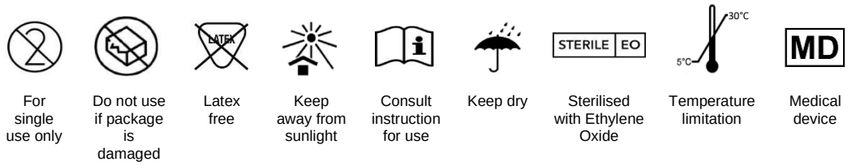
For single use only | Do not use if package is damaged | Latex free | Keep away from sunlight | Consult instruction for use | Keep dry | Sterilised with Ethylene Oxide | Temperature limitation | Medical device

CE 1014

Manufactured by:
DINA-HITEX, spol. s r.o.
Ždánská 987
685 01 Bučovice,
Czech Republic

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