Suresite®
Window Dressing

Instructions for use

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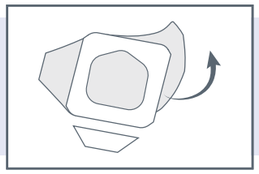
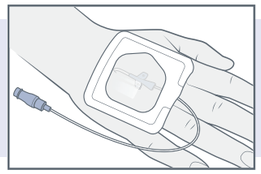
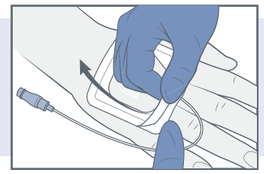
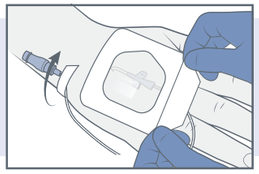
en

Indications:
Intended to secure and maintain catheters. Can be left in place for up to 7 days.

General Precautions & Observations:
Can be used as a primary (stand alone) or secondary dressing (with one fixation method) .
Do not stretch the film during application.
Should any signs of redness or irritation appear, discontinue use and consult a healthcare professional.

Contraindications:
· Do not use on damaged skins or wounds
· Those devices are not wound dressing
· Must not be used on neonates or burn patients

Warning:
Using Suresite Transparent Film Dressing to cover a drug delivery patch could lead to faster absorption, increasing the risk of overdose. Reuse of single use devices creates a potential risk for the patient or user. It may cause contamination and/or impairment of function; which may lead to injury, illness and/or death of a patient.
If a serious incident has occurred in relation to the device, it should be reported to Medline and the competent authority of the Member State in which the user and/or patient is established.

Disposal instructions:
Precautions should be taken when discarding this device and disposal of the device shall be made in accordance with applicable hospital or national regulations for biologically hazardous waste.

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Medline International France SAS,
UKCA 0086 **UK representative and importer in UK:** Medline Industries Ltd, WA4 6HL, UK