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
400mm

180mm

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FFP2 NR
FILTERING HALF MASK

- en Instructions for Use
- fr Mode d'emploi
- nl Gebruiksaanwijzingen
- de Gebrauchsanweisung
- it Istruzioni per l'uso
- es Instrucciones de uso
- pt Instruções de utilização

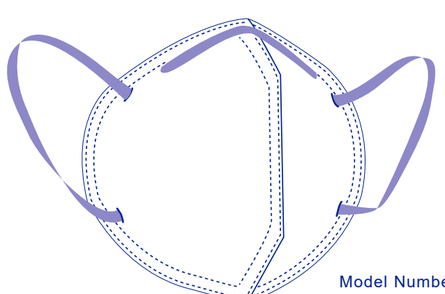Model Number:
ZD-001

Notified Body: UNIVERSAL CERTIFICATION and SURVEILLANCE SERVICES Trade Co. NB:2163
Address: Necip Fazıl Bulvarı Keyap Sitesi E2 Blok No:44/84 Yukarı Dudullu Ümraniye/İSTANBUL – TURKEY

Manufacturer:
Zhende Medical Co., Ltd.
Tel: +86-0575-88771820
Gaobu Town, 312035, Shaoxing City, Zhejiang Province, China
EU Declaration of Conformity Internet download address: www.zhende.com

The Manufactured date: Month/ Year
The Expiry time: 3 years
Version: 02
Rev. Date: 2020.08.11

Instructions for Use

Product Name: FILTERING HALF MASK
Model: ZD-001
The filtering half mask is mainly composed of spunbonded nonwoven, melt-blown nonwoven, hot-air nonwoven, nose clamp, ear loop, and adjustable buckle. When inhaling, the low air pressure inside the mask causes airflow to enter the mask. The air pressure in the mask is higher than the ambient air pressure when exhaling.
Intended Use:
Protective mask is intended to protect against particles, protect against both solid and liquid aerosols
Warning:
1. Select the respirator type and determine the suitability of the respirator according to the intended use.
2. The respirator shall be abandoned when it is broken, becomes difficult to breathe, or polluted by body fluids.
3. This respirator is disposable, do not reuse or use for more than one shift.
4. Check the respirator before use. Make sure there are no broken ear loops, no holes or other damages on the respirator. If any damage detected, please do not use and replace with a new one immediately.
5. The effectiveness will be reduced if facial hair or beard passes under the face seal.
6. The respirator should only be used in intended environment. Failure to follow the instruction may cause harm to health or reduce the effectiveness.
7. The respirator cannot provide intended protection if the mask is not correctly donned or is damaged.
8. Never alter or modify this respirator in any way.
9. This respirator does not supply oxygen. Use only in adequately ventilated areas containing sufficient oxygen to support life.
10. The respirator has a shelf-life of 3 years. After the expiry date the respirator must no longer be used.
11.Do not use in explosive atmospheres.
Storage Condition and Methods:
1.The unused masks need to be stored in the environment that is well ventilated, dry and away from light. Keep away from fire, pollutants and potential source of pollution.
2. Storage temperature: -20℃ to 50 ℃ and less than 80% relative humidity.
3. Keep the masks in original packaging boxes during transportation to protect against mechanical damage and contamination.
Limitations of Use: None
The respirator complies with Regulation (EU) 2016/425, and conforms with EN 149:2001+ A1:2009.
Application Instruction:

1. Take out the respirator and pull on adjustable buckle with wire or other small gadget.
2. Place respirator over face with the nose clamp sitting upward as shown. Position both ear loops behind your ears, making sure your lower chin is covered.
3. Hold the nose clamp to the shape of your nose to ensure it fits your face comfortably.
4. Adjust the respirator with adjustable buckle to a proper position and ensure there is no leakage.
5. To test fit: cup both hands over the respirator and inhale vigorously,a negative pressure should be felt. If air leaks at the edge,reposition the mask.

Notified Body: UNIVERSAL CERTIFICATION and SURVEILLANCE SERVICES Trade Co. NB:2163
Address: Necip Fazıl Bulvarı Keyap Sitesi E2 Blok No:44/84 Yukarı Dudullu Ümraniye/İSTANBUL – TURKEY

Humidity limitation | Temperature limit | Keep Dry | Away from direct sunshine | Single Use | See information supplied by the manufacturer | CE 2163

Manufacturer:
Zhende Medical Co., Ltd.
Tel: +86-0575-88771820
Gaobu Town, 312035, Shaoxing City, Zhejiang Province, China
EU Declaration of Conformity Internet download address: www.zhende.com

The Manufactured date: Month/ Year
The Expiry time: 3 years
Version: 02
Rev. Date: 2020.08.11

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