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
400mm

200mm

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FFP2 NR
PARTICULATE RESPIRATORS

- en Instructions for Use
- fr Mode d’emploi
- nl Gebruiksaanwijzingen
- de Gebrauchsanweisung
- it Istruzioni per l’uso
- es Instrucciones de uso
- pt Instruções de utilização

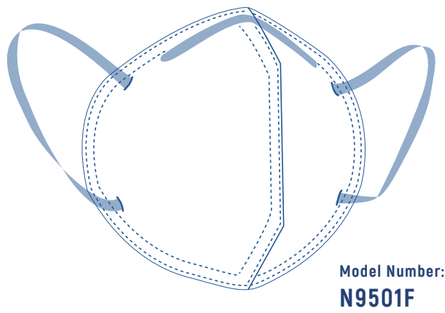Model Number:
N9501F

Notified Body: BSI Group The Netherlands B.V
Address: Say Building, John M.Keynesplein 9,1066 EP Amsterdam P.O. Box 74103, 1070 BC Amsterdam The Netherlands

Manufacturer:
Zhende Medical Co., Ltd.
Gaobu Town, 312035, Shaoxing City, Zhejiang Province, China
EU Declaration of Conformity Internet download address:
http://www.zhende.com/EN/News/1/94.html

European Respresentative:
MedPath GmbH
Mies-van-der-Rohe-Strasse 8,
80807 Munich, Germany

ZD logo marking on the product indicates Zhende Medical Co., Ltd.
This IFU will be translanted in all languages of the intended countries for sales or copies of the translations.

Version: 03
Rev. Date: 2021.09.25

Instructions for Use

Product Name: FFP2 NR Particulate Respirators
FFP2 means the performance level of the mask, and NR means non-reusable

Model: N9501F

Intended Use:
Filtering half masks, as RPE, use filters to remove contaminants from the air being breathed in, and to protect against particles except for escape purposes.

Warning:
1. Select the respirator type and determine the suitability of the respirator according to the intended use.
2. The respirator shall be abandoned when it is broken, becomes difficult to breathe, or polluted by body fluids.
3. This respirator is disposable, do not reuse or use for more than one shift.
4. Check the respirator before use. Make sure there are no broken ear loops, no holes or other damages on the respirator. If any damage detected, please do not use and replace with a new one immediately.
5. The effectiveness will be reduced if facial hair or beard passes under the face seal.
6. The respirator should only be used in intended environment. Failure to follow the instruction may cause harm to health or reduce the effectiveness.
7. The respirator cannot provide intended protection if the mask is not correctly donned or is damaged.
8. Never alter or modify this respirator in any way.
9. This respirator does not supply oxygen. Use only in adequately ventilated areas containing sufficient oxygen to support life.
10. The respirator has a shelf-life of 3 years. After the expiry date the respirator must no longer be used.
11.Do not use in explosive atmospheres.

Storage:
1.The unused masks need to be stored in the environment that is well ventilated, dry and away from light. Keep away from fire, pollutants and potential source of pollution.
2. Storage temperature: -20℃ to 50 ℃ and less than 80% relative humidity.
3. Keep the masks in original packaging boxes during transportation to protect against mechanical damage and contamination.

Limitations of Use: Do not use the product in explosive atmospheres. Discard and replace respirator if it becomes damaged or breathing resistance becomes excessive.
The respirator complies with Regulation (EU) 2016/425, and conforms with EN 149:2001+ A1:2009.

Application Instruction:

Notified Body: BSI Group The Netherlands B.V
Address: Say Building, John M.Keynesplein 9,1066 EP Amsterdam P.O. Box 74103, 1070 BC Amsterdam The Netherlands

Manufacturer:
Zhende Medical Co., Ltd.
Gaobu Town, 312035, Shaoxing City, Zhejiang Province, China
EU Declaration of Conformity Internet download address:
http://www.zhende.com/EN/News/1/94.html

European Respresentative:
MedPath GmbH
Mies-van-der-Rohe-Strasse 8,
80807 Munich, Germany

ZD logo marking on the product indicates Zhende Medical Co., Ltd.
This IFU will be translanted in all languages of the intended countries for sales or copies of the translations.

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Organismo notificado: BSI Group The Netherlands B.V
A marcação do logotipo ZD no produto designa a Zhende Medical Co., Ltd.