en MARATHON® XL

NO-STING CYANOACRYLATE SKIN PROTECTANT

Product Information Data Sheet

FOR EXTERNAL USE ONLY.

Product Description:
Marathon® XL No-Sting Cyanoacrylate Skin Protectant is a non-cytotoxic, fast drying, liquid barrier film for the protection of damaged or intact skin. It is applied as a liquid and dries, within a minute, adhering to the contours of the skin to form a transparent flexible film. Marathon® XL No-Sting Cyanoacrylate Skin Protectant will wear off naturally as the skin regenerates.

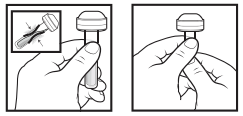
Intended Use:
- Marathon® XL No-Sting Cyanoacrylate Skin Protectant is intended to protect intact or damaged skin from the effects of moisture, friction (rubbing), shear (tearing).

Indications For Use:
- Marathon® XL No-Sting Cyanoacrylate Skin Protectant helps protect skin exposed to irritation from moisture such as urine, feces, digestive juices, perspiration and wound drainage.
- Marathon® XL No-Sting Cyanoacrylate Skin Protectant can also be used in areas that are exposed to friction and shear from bedding, clothing, shoes, or any material that would rub against the skin.
- Marathon® XL No-Sting Cyanoacrylate Skin Protectant helps protect the skin against irritation caused by adhesive products.

Contraindications:
- Do NOT apply directly to deep, open, or bleeding wounds.
- Do NOT apply to chronic wounds (persistent or unhealing wounds).
- Do NOT apply to second or third degree burns.
- Do NOT apply to infected areas.
- Do NOT use on individuals with a known sensitivity to cyanoacrylates or formaldehyde.

Instructions For Use:
- Clean and dry the skin before applying Marathon® XL No-Sting Cyanoacrylate Skin Protectant. Make sure that the skin is thoroughly dry and that all traces of moisturizers, lotions, lotion soaps or lanolin are removed from the application site.
- Protect any sensitive areas, such as stomas and mucous membranes (e.g. the gums and inside of mouth).
- Remove the applicator from the packaging. While holding the applicator upright (with the sponge tip to the top), firmly squeeze between the finger and thumb to crush the sealed inner tube.

Correct Activation:

One handed method

When two hands are needed

INCORRECT ACTIVATION:
DO NOT try to bend or snap the tube to activate.

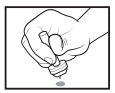
Application Methods:
Standard:
1) Turn ampule upside down (with tip to the floor) and gently squeeze to prime the tip.
2) Once the tip is primed with Marathon® XL, begin to gently apply in a thin, even coat that extends at least an inch (2.5 cm) beyond desired area.

Drip Technique:
1) Turn ampule upside down (with the tip to the floor) and gently squeeze to prime the tip.
2) Once tip is primed, position the ampule over the desired application area, just above the surface, and squeeze applicator to drip Marathon® XL onto the application site.
3) Use the tip to spread the Marathon® XL droplets evenly into a thin, even coat that extends at least an inch (2.5cm) beyond the desired area.

Standard method

Drip method

- If Marathon® XL No-Sting Cyanoacrylate Skin Protectant is used at room temperature it will dry within a minute. Some people may feel a temporary warming sensation as the cyanoacrylate dries.
- Do not over apply or rub the Marathon® XL No-Sting Cyanoacrylate Skin Protectant in like a cream as the coverage will not be as effective and the warming sensation will be increased. It is most effective when the coating is paper-thin.
- Ensure that Marathon® XL No-Sting Cyanoacrylate Skin Protectant has dried, completely, before covering with clothes, applying adhesive products such as tapes and dressings, or incontinence products.
- Marathon® XL No-Sting Cyanoacrylate Skin Protectant will last up to 3 days depending on skin type, activity level, and type of fluid that comes in to contact with the barrier. Inspect the film regularly (at least once a day) and if any cracking or flaking appears, apply another layer on top, until the area is completely covered once again. It is not necessary to remove the cracked and flaking layer, however, if it is applied too frequently excess build-up will lead to a thick inflexible film.
- One Marathon® XL No-Sting Cyanoacrylate Skin Protectant applicator will cover an area of approximately 8 in x 10 in (20 cm x 25 cm).
- Do not use cleansers that leave a protective barrier, such as dimethicone, petrolatum, mineral oil, etc.

Warnings:
- If irritation or pain develops (redness, inflammation, heat, or change in skin color) immediately discontinue product use. If irritation persists, seek medical attention.
- Only use for up to 30 days on damaged skin.
- For external use only.
- Do not re-use in whole or part as it may compromise the performance of the protectant.
- If Marathon® XL No-Sting Cyanoacrylate Skin Protectant is applied to a moist skin surface it will set more rapidly and a noticeable, but temporary warming sensation will be felt.
- Do not use under drug delivery patches.
- Avoid contact with mucous membranes such as the eyes and mouth.
- Keep out of reach of children.

Precautions:
- To avoid unintentional adherence, Marathon® XL No-Sting Cyanoacrylate Skin Protectant should dry for at least one minute before coming in to contact with skin and materials (gloves, instruments, bedding, clothing, etc.)
- This device is intended to be used by healthcare professionals or other users as directed.

Adverse Reactions:
- If Marathon® XL No-Sting Cyanoacrylate Skin Protectant enters the eye, immediately flush with warm water. Do not rub, or try to separate the eyelids. Seek medical attention.
- If Marathon® XL No-Sting Cyanoacrylate Skin Protectant accidentally contacts an open wound, flush with warm tap water. Removal is not necessary.
- If the skin or other materials become stuck to the area, gently roll apart to separate. If unable to separate easily, seek medical attention.
- If pain occurs, and persists, seek medical attention.

Removal:
- Excess Marathon® XL No-Sting Cyanoacrylate Skin Protectant can be removed before it dries by wiping with a clean tissue. Once it has dried, unwanted cyanoacrylate can be softened with petroleum jelly and removed. It will naturally wear off within 1-3 days as the skin renews itself.

Storage:
- Controlled room temperature (5-25 °C/41-77 °F).
- For shelf life refer to expiration date on product.

Contents:
- Each applicator contains 1.5g of a cyanoacrylate formulation.
- Each box contains 5 applicators of Marathon® XL No-Sting Cyanoacrylate Skin Protectant.

Disposal:
- Precautions should be taken when discarding this device and disposal of the device shall be made in accordance with applicable hospital or national regulations for biologically hazardous waste.

CAUTION: Federal (USA) law restricts this device to sale by or on the order of a Physician or Practitioner.

If a serious incident has occurred in relation to the device, it should be reported to Medline and the competent authority of the Member State in which the user and/or patient is established.

Reuse of single-use devices creates potential risk for the patient or user. It may cause contamination and/or impairment of function that may lead to injury, illness and/or death of a patient.

REF **MSC093001XL**

MD **Medical device**
Producto Sanitario
Dispositif médical
Medizinprodukt
Medisch hulpmiddel
Dispositivo Medico

www.medline.com/awc
Medline and Marathon® are registered trademarks of Medline Industries, LP.
Made in UK for: Medline Industries, LP,
Three Lakes Drive, Northfield,
IL 60093 USA. **1-800-MEDLINE** V5 RF23MLG

CE UKCA **UK representative and importer in UK:**
Medline Industries Ltd, WA4 6HL, UK

EC REP Medline International France SAS,
5 rue Charles Lindbergh, 44110
Châteaubriant, France

CH REP **CH Rep and importer:** Medline International Switzerland SÀRL, c/o MN & Associés SA, 1 Place de Longemalle, 1204 Genève, Switzerland

Medline International B.V., Nieuwe Stationsstraat 10, Arnhem, 6811 KS, Netherlands

Symbol Definition

Note: The following symbols and definitions come from ISO15223.-1: 2016

TEMPERATURE LIMITATION

NOT MADE WITH NATURAL RUBBER LATEX

MANDATORY ACTION TO READ THE INSTRUCTIONS FOR USE

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