
ASID BONZ
A MEDLINE COMPANY

en Impervious Cover Gown, Latex Free
fr Blouse de protection imperméable, Sans latex
de Flüssigkeitsundurchlässiger Schutzkittel, Latexfrei
it Camice da copertura impervio, Senza lattice
es Bata de aislamiento impermeable, Sin látex
pt Bata de isolamento impermeável, Não contém látex
nl Ondoordringbare beschermjas, LatexVrij
sv Skyddsrock, vätsketät, Latexfri
da Uigennemtrængelig overtrækskittel, Latexfri
fi Läpäisemätön suojatakki, Lateksiton
el Αδιάβροχη ρόμπα κάλυψης, Χωρίς λάτεξ
no Ugjennomtrengelig beskyttelsesfrakk, Latexfri
hu Vízálló fedoköpeny, Latex mentes
ro Halat impermabil pentru vizitator, Fara latex
pl Nieprzepuszczalna koszula ochronna, Nie zawiera lateksu
cs Nepropustný krycí plášť, Bez latexu
sl Neprepustna vrhnja halja, Brez lateksa
sk Nepriepustný krycí plášť, Bez latexu

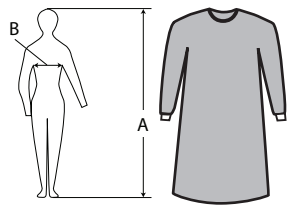Finess REF	ASID REF	SIZE	A	B
FN-IG0102060201	760510-N	130 cm x 140 cm	165-190cm	95-120cm
FN-IG0102060202	760520-N	130 cm x 160 cm	165-190cm	95-120cm

EU declaration of conformity available on:http://www.finess-healthcare.com/

CE Medical Device class I non-sterile under MDR 2017/745
2777 PPE category III under Regulation (EU) 2016/425
Certified by Satra Technology Europe Ltd, Bracetown Business Park, Clonee, Dublin 15 Dublin Ireland

Imported/Distributed by: ASID BONZ GmbH
Hertzstraße 3, DE-71083 Herrenberg, Germany
www.asid-bonz.de

Finess Healthcare Group Co., Ltd
NO.1503 A1 Block, NO.203 Lantian Road, 315000 Ningbo, PEOPLE'S REPUBLIC OF CHINA

EC REP ZOUSTECH S. L., Pso, Castellana
141-Planta 19,28046 Madrid, SPAIN

LAB190414 - 2022.02

en

Personal Protective Equipment of Complex Design Category III, in compliance with Regulation (EU) 2016/425 - Type [PB]-6

The shelf life is given on the product label. The product is not known to cause skin reactions. Type PB [6] partial body protection has not been tested to the whole suit test.

These gowns are made from materials with good breathability properties to reduce heat stresses and designed for comfortability.

Intended use:
Isolation gown is made of PP Non woven and PE film lamination garment,worn by healthcare providers or visitors to protect the wearer from the spread of infection or illness if the wearer comes in contact with potentially infectious liquid material.They may also be used to help prevent the gown wearer from transferring microorganisms that could harm vulnerable patients, such as those with weakened immune systems.This is a dispoable garment."

Storage conditions: 5 - 38°C.

Precautions and warnings:
For single use only, Do not reuse. Flammable material. Keep away from fire.

Disposal:
The used protective gowns have to be disposed as hazardous waste. Special marked containers have to be used.
Local regulations must be considered.
Gowns must remain in their primary packaging until required

Declaration of conformity: http://www.finess-healthcare.com/

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